



PRODUCTION EXECUTION

RELEASE 10.02.00

USER MANUAL

PUBLICATION PSPE-UM010C-EN-E-SEPTEMBER-2021

Supersedes publication PSPE-UM010B-EN-E



Contact Rockwell See contact information provided in your maintenance contract.

Copyright Notice © 2021 Rockwell Automation Technologies, Inc. All rights reserved.
This document and any accompanying Rockwell Software products are copyrighted by Rockwell Automation Technologies, Inc. Any reproduction and/or distribution without prior written consent from Rockwell Automation Technologies, Inc. is strictly prohibited. Please refer to the license agreement for details.

Trademark Notices FactoryTalk, PharmaSuite, Rockwell Automation, Rockwell Software, and the Rockwell Software logo are registered trademarks of Rockwell Automation, Inc.

The following logos and products are trademarks of Rockwell Automation, Inc.:

FactoryTalk Shop Operations Server, FactoryTalk ProductionCentre, FactoryTalk Administration Console, FactoryTalk Automation Platform, and FactoryTalk Security.

Operational Data Store, ODS, Plant Operations, Process Designer, Shop Operations, Rockwell Software CPGSuite, and Rockwell Software AutoSuite.

Other Trademarks ActiveX, Microsoft, Microsoft Access, SQL Server, Visual Basic, Visual C++, Visual SourceSafe, Windows, Windows 7 Professional, Windows 10, Windows Server 2008, Windows Server 2012, and Windows Server 2016 are either registered trademarks or trademarks of Microsoft Corporation in the United States and/or other countries.

Adobe, Acrobat, and Reader are either registered trademarks or trademarks of Adobe Systems Incorporated in the United States and/or other countries.

ControlNet is a registered trademark of ControlNet International.

DeviceNet is a trademark of the Open DeviceNet Vendor Association, Inc. (ODVA).

Ethernet is a registered trademark of Digital Equipment Corporation, Intel, and Xerox Corporation.

OLE for Process Control (OPC) is a registered trademark of the OPC Foundation.

Oracle, SQL*Net, and SQL*Plus are registered trademarks of Oracle Corporation.

All other trademarks are the property of their respective holders and are hereby acknowledged.

Warranty This product is warranted in accordance with the product license. The product's performance may be affected by system configuration, the application being performed, operator control, maintenance, and other related factors. Rockwell Automation is not responsible for these intervening factors. The instructions in this document do not cover all the details or variations in the equipment, procedure, or process described, nor do they provide directions for meeting every possible contingency during installation, operation, or maintenance. This product's implementation may vary among users.

This document is current as of the time of release of the product; however, the accompanying software may have changed since the release. Rockwell Automation, Inc. reserves the right to change any information contained in this document or the software at any time without prior notice. It is your responsibility to obtain the most current information available from Rockwell when installing or using this product.

Chapter 1	Production Execution.....	1
	Typographical Conventions	2
	Screen Areas.....	3
	Header Bar.....	3
	Work Area.....	5
	Notification Panel Area.....	6
	Control Bar	6
	Framework for Pre-defined Workflows	8
	Graphical Elements	10
	Framework.....	10
	Cockpit	13
	Notification Panel	21
	Common.....	22
	Specific.....	26
	Cockpit	29
	Navigator.....	38
	Selection Window	41
	Execution Window.....	42
	Exception Window.....	44
	Phase-specific Exception Window	44
	Operation-specific Exception Window	46
	Overview Exception Window	47
	Notification Panel	48
	Detail Information Window	50
	Input Devices	51

Keyboard.....	51
Scanner.....	56
Touch Screen.....	57
Basic Operations.....	57
Start, Login, Logout, and Password Change	57
Executing a Unit Procedure.....	62
Handling Exceptions.....	63
System-triggered Exception	64
User-triggered Exception.....	66
Aborting and Repairing Phases	76
Accessing Cockpit Actions.....	78
Pausing Unit Procedures.....	84
Appending Workflows to Orders	86
Handling Escalations	89
Detaching Unit Procedures	91
Aborting and Reactivating Unit Procedures.....	92
Detaching Operations.....	93
Aborting Operations.....	94
Signature Requests.....	95
Reference Selectors.....	99
Help Access.....	102
Chapter 2 Production Processes.....	105
Batch Processing.....	105
Material Processing.....	107
Workflow Processing	109
Chapter 3 Inventory Workflows.....	111
Material Receipt	111
What Is Material Receipt?	111
Material Receipt: Filter Materials	112
Material Receipt: Select a Material	113

Select or Create a Batch	114
Material Receipt: Filter Storage Locations.....	116
Material Receipt: Select a Storage Location.....	118
Material Receipt: Create Sublots	118
Material Receipt: Print Labels	119
Material Issue	119
What Is Material Issue?	119
Material Issue: Filter Sublots	120
Material Issue: Select Sublots	122
Material Issue: Issue Sublots.....	123
Relocation	123
What Is Relocation?	123
Relocation: Filter Sublots	123
Relocation: Select Sublots	126
Relocation: Relocate Sublots.....	127
Sublot Split	127
What Is Sublot Split?	127
Sublot Split: Filter Sublots	128
Sublot Split: Select a Sublot.....	130
Perform an Even or Uneven Split	130
Sublot Split: Print Labels.....	131
Inventory Correction	132
What Is Inventory Correction?.....	132
Inventory Correction: Filter Sublots	132
Inventory Correction: Select a Sublot	134
Index	135

Figure 1: Screen layout	3
Figure 2: Header bar	5
Figure 3: Header bar with treatment ID	5
Figure 4: Notification panel header bar with alarms	5
Figure 5: Work area control bar	7
Figure 6: Notification panel control bar	8
Figure 7: Screen layout	8
Figure 8: About dialog.....	30
Figure 9: About dialog - Details.....	30
Figure 10: Register at Station dialog	31
Figure 11: Change User dialog	32
Figure 12: Configure Screen Size dialog.....	33
Figure 13: Cockpit.....	34
Figure 14: Running processes	37
Figure 15: Navigator	40
Figure 16: Data panel of a non-order-related workflow	40
Figure 17: Selection Window.....	41
Figure 18: Phase representation during execution.....	43
Figure 19: Phase with required completion signature	43
Figure 20: Phase-specific Exception Window	45
Figure 21: Operation-specific Exception Window	46
Figure 22: Overview Exception Window.....	47
Figure 23: Detail Information Window	51
Figure 24: Keyboard mapping of navigation bar buttons.....	53
Figure 25: Scan in progress	56
Figure 26: PharmaSuite webstart page.....	57
Figure 27: Login form.....	58

Figure 28: PharmaSuite welcome page.....	59
Figure 29: Change password.....	61
Figure 30: Error and error details dialogs.....	64
Figure 31: Information, question, and exception dialogs.....	65
Figure 32: Exception recording panel with free-text entry.....	68
Figure 33: Exception recording panel with pre-defined text options.....	68
Figure 34: Phase repair mode with Option List editor.....	76
Figure 35: Actions dialogs of order unit procedures	79
Figure 36: Actions dialogs of workflow unit procedures	80
Figure 37: Actions dialog of operations	81
Figure 38: Actions dialog of an event-triggered operation	83
Figure 39: Appendable Unit Procedures.....	87
Figure 40: Appendable Workflows.....	89
Figure 41: Single electronic signature.....	95
Figure 42: Single electronic signature with pre-defined comment.....	96
Figure 43: Double electronic signature to support witness role	96
Figure 44: Single phase completion signature - active	97
Figure 45: Single phase completion signature - completed	97
Figure 46: Sequential phase completion signature - active	98
Figure 47: Sequential phase completion signature - completed.....	98
Figure 48: Phase action signature.....	98
Figure 49: Single-reference box with selector button	99
Figure 50: Search criteria.....	100
Figure 51: Sorting criteria	101
Figure 52: Single-reference selector.....	102
Figure 53: Help window navigation	103

Production Execution

PharmaSuite for Production Execution is a system designed to control the production on the shop floor by managing its individual processes.

This section contains important information about the basic principles of working with PharmaSuite for Production Execution. Please read this section carefully, because it provides a solid background for all operations you may wish to perform with your system.

Later sections will explain how to work with the system. We assume you are familiar with the conventions described in the following sections and the fundamentals of working with a touch-screen or keyboard-only application.

Typographical Conventions

This documentation uses typographical conventions to enhance the readability of the information it presents. The following kinds of formatting indicate specific information:

Bold typeface	Designates user interface texts, such as <ul style="list-style-type: none">■ window and dialog titles■ menu functions■ panel, tab, and button names■ box labels■ object properties and their values (e.g., status).
<i>Italic typeface</i>	Designates technical background information, such as <ul style="list-style-type: none">■ path, folder, and file names■ methods■ classes.
CAPITALS	Designate keyboard-related information, such as <ul style="list-style-type: none">■ key names■ keyboard shortcuts.
Monospaced typeface	Designates code examples.

Screen Areas

The basic screen layout of Production Execution with EBR consists of three areas:

- Header bar (page 3)
- Work area (page 5)
- Notification panel area (page 6)
- Control bar (page 6)

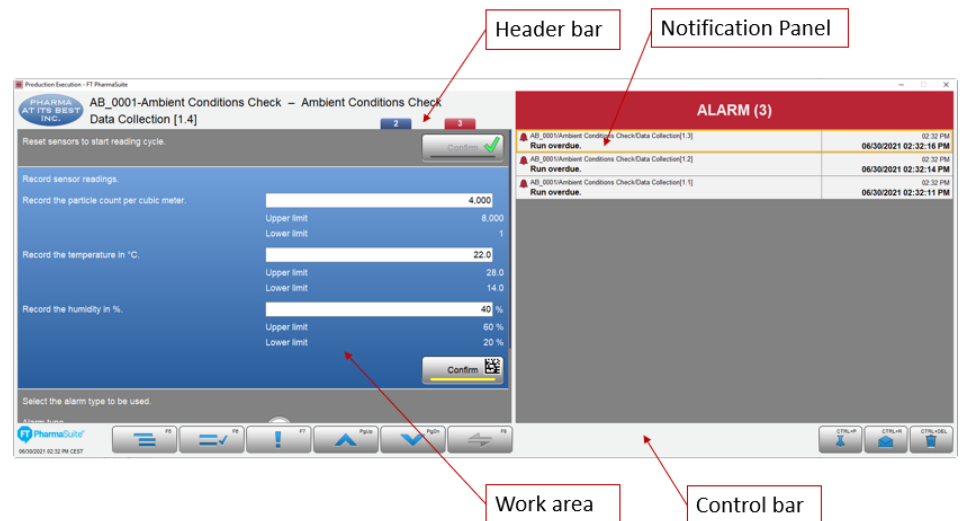


Figure 1: Screen layout

Header Bar

The **header bar** is located at the top of the screen and contains

- a company logo
- a screen title with up to three lines of text that refers to the content of the work area:
 - for the Cockpit (page 29) and the Navigator (page 38)
 - first line: surname and login name of the logged-in user.
 - second line: identifiers of the current work center and station.

- for the Execution Window (page 42)
 - first line for order steps:
identifier of the current order step, with its unit procedure count if it has been reactivated, followed by the product batch.
first line for workflows:
identifier of the current workflow with its processing name.
 - second line:
identifier of the operation with its count if it has been run more than once.
 - third line (only if a treatment ID is defined for the order or workflow)
treatment ID.
- for the Detail Information Window (page 50)
 - first line for order steps:
identifier of the current order step, with its unit procedure count if it has been reactivated, followed by the product batch.
first line for workflows:
identifier of the current workflow with its processing name.
 - second line:
identifier of the operation, with its count if it has been run more than once, followed by the identifier of the selected phase.
- for the phase-specific Exception Window (page 44)
 - first line:
identifier of the current order step or workflow, with its unit procedure count if it has been reactivated, for order steps followed by the product batch.
 - second line:
identifier of the operation with its count if it has been run more than once and phase identifier.
- for the overview Exception Window (page 47)
 - first line:
identifier of the current order step or workflow, for order steps followed by the product batch.
 - second line:
identifier of the unit procedure from which it has been accessed with its count if it has been reactivated.
- for a Selection Window (page 41)
 - selection-specific title and optional information, such as identifiers of affected objects or the affected work center and station

- a notification bar to display escalation stages (page 89) for the runs of event-triggered operations.
- an alarm indication bar (ALARM including the number of active alarms on a red background) to display that at least one alarming notification is present in the notification panel.

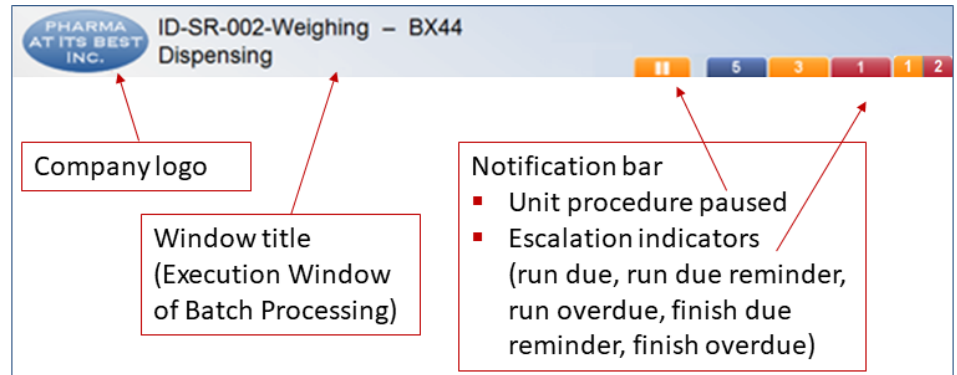


Figure 2: Header bar

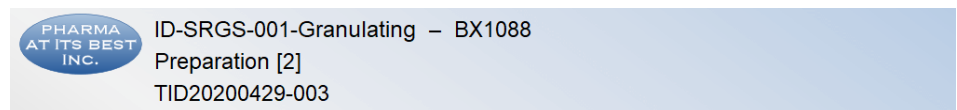


Figure 3: Header bar with treatment ID

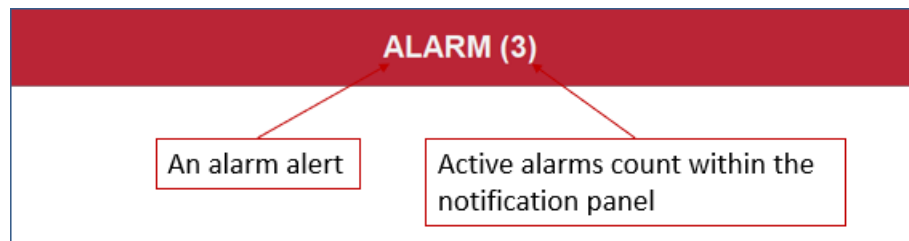


Figure 4: Notification panel header bar with alarms

Work Area

The **work area** is the main area in the middle of the screen and displays all dynamic data and controls.

The system indicates the currently focused control in yellow, as background color for input boxes and indicator bar for buttons.

- For the Cockpit (page 29), it contains the lists of **Startable** and **Running** processes and workflows.
- For the Navigator (page 38), it contains the history of executed steps with information and action buttons.

- For the Execution Window (page 42), it contains the phases of the active operation, with display text, input boxes, buttons, and all controls required for processing.
- For the Exception Window (page 44), it displays all exceptions and comments and provides an exception recording panel.
- For the Detail Information Window (page 50), it displays the relevant section of the batch report, along with existing exceptions, comments, and an exception recording panel.
- For a Selection Window (page 41), it displays the list of objects available for selection. Above the list, some Selection Windows provide filter functions to restrict the number of listed objects.

Notification Panel Area

The **Notification Panel** (page 48) is located on the right side of the screen and contains a notification table with two columns.

- for the Notification text column, it contains:
 - the context to which the notification belongs, described by a path,
 - the notification text defined by the sender of the notification,
 - indicators for pinned notifications and alarming notifications.
- for the Notification value column, it contains:
 - the timestamp when the notification was created or modified,
 - the timestamp defined by the sender, or
 - the timer that counts down in minute increments.

Control Bar

The **control bar** is located at the bottom of the screen and contains

- the PharmaSuite logo
- the current date and time with time zone information
- the control buttons for working with the system:
 - **Cockpit** button to access the Cockpit (page 29)
 - **Navigator** button to access the Navigator (page 38)
 - **Exception** button to register an exception (page 63) in the process

TIP

In the Cockpit, Navigator, or Exception Windows the respective buttons are replaced by a **Back** button to return to the Execution Window.

- **Page up/Page down** buttons to page through the content of the work area
- **Switch** button to switch between concurrently running operations

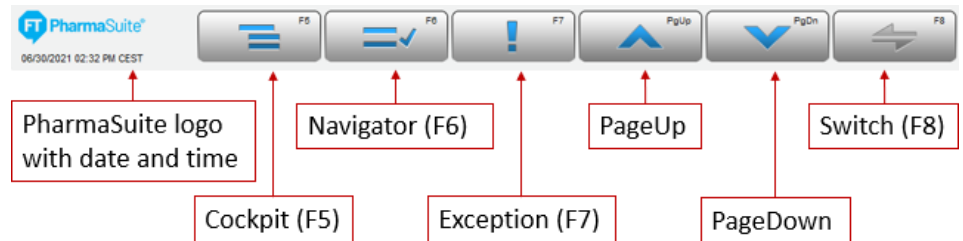


Figure 5: Work area control bar

- the notification panel buttons:
 - Pin and unpin notification toggle buttons:
 - The **Pin notification** button (page 22) moves the selected notification to the top of the list, on the current device. A pinned notification displays the Pinned icon.
 - The **Unpin notification** button (page 22) moves the selected notification down in the list, on the current device. The **Pinned** icon is removed when a notification is unpinned.
 - Mark notification as read/unread toggle buttons:
 - The **Mark notification as read** button (page 21) marks selected notification as read, on the current device by changing the background color of its table gray, as well as its text.
 - The **Mark notification as unread** button (page 21) restores the initial coloring of the table row of the selected notification to indicate that it is not read on the current device.
 - The **Remove notification** button (page 22) removes selected notification from the list, on the current device.

Enabled buttons show a blue indicator bar or graphic, enabled and focused buttons turn their indicators yellow, while disabled buttons are grayed out.

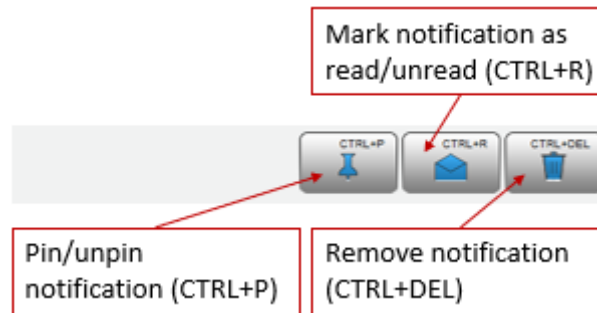


Figure 6: Notification panel control bar

Framework for Pre-defined Workflows

The basic layout of the framework for pre-defined workflows consists of:

- status bar with work center information
- work area with access to Cockpit (page 29) and pop-up panel for help
- navigation bar with action buttons

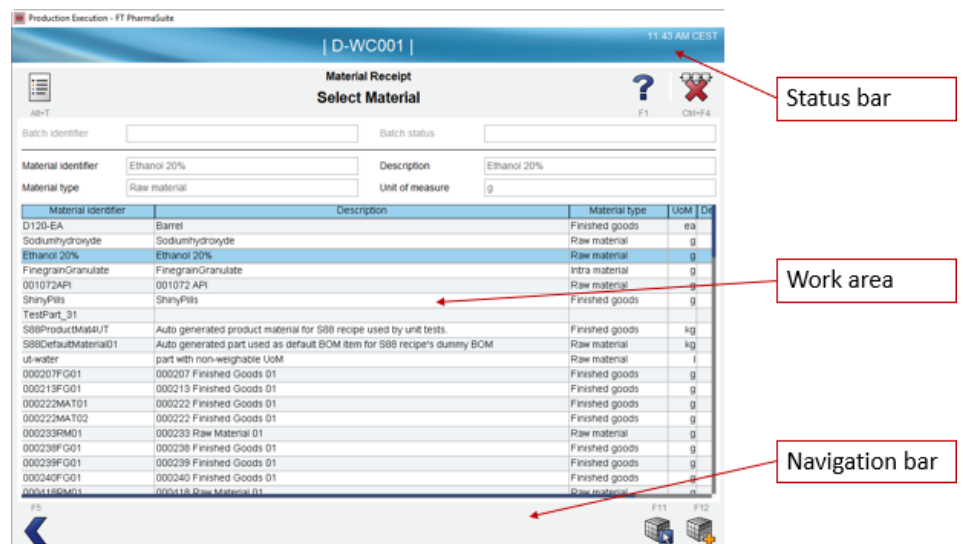


Figure 7: Screen layout

STATUS BAR

The status bar provides workflow-independent information such as the work center identifier and time and time zone information.

WORK AREA

The work area provides all information and controls a user needs for his current workflow step. It may display input boxes and grids for data display and selection, as well as graphics and graphical controls, such as a scale control, etc.

In the upper left corner of the work area, the system provides the **Open cockpit** button to open and close the Cockpit for starting workflows or switching between workflows. The upper right corner holds two buttons, the **Help** button to access the context-sensitive help system and the **Cancel workflow** button to stop processing the current workflow and close the form.

NAVIGATION BAR

The navigation bar provides the controls to move through a workflow. The **Back** button to the far left of the bar takes the user one step back in the workflow, provided this is allowed by the workflow structure. The **Next** button to the far right of the bar confirms the completion of the active step and takes the user one step forward in the workflow. If there are several options available for how to proceed from a workflow step, additional buttons indicate the different options, sometimes even replacing the **Next** button.

In addition to the workflow control buttons, the navigation bar can display indicator icons, which point to specific system statuses or events like "system ready to scan" or "scanning is mandatory".

Graphical Elements

The following lists contain the graphical buttons and icons used in PharmaSuite for Production Execution.


Framework

The framework provides the following buttons and graphical elements to control and navigate processes and workflows:



Back

In pre-defined workflows, returns to the previous step, if this is permitted by the workflow.

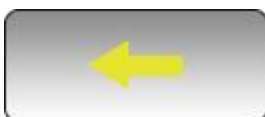

Back (F7)

In the control bar, returns to the window from which you have accessed the currently active window. The **Back** button can replace the **Cockpit**, **Navigator**, or **Exception** buttons, respectively.



Cancel workflow (CTRL+F4)

In pre-defined workflows, cancels the current workflow. This action depends on the workflow step and is not always available.



Cockpit (F5)

In the control bar, opens the Cockpit (page 29) to display the lists of **Startable** and **Running** processes and workflows.



Exception (F7)

In the control bar, opens the Exception Window (page 44) to record an exception to the regular execution process.

**Help (F1)**

In pre-defined workflows, opens the pop-up window with the context-sensitive help system.

**Navigator (F6)**

In the control bar, opens the Navigator (page 38) to list the current and completed phases of the operation.

**Next**

In pre-defined workflows, completes the current step and moves to the next step in the workflow.

5**Notification bar - due**

In the header bar, attached to the top right corner of the work area, indicates the number of runs of an event-triggered operation that are ready for processing.

3**Notification bar - due reminder**

In the header bar, attached to the top right margin of the work area, indicates the number of runs of an event-triggered operation that need to be processed as soon as possible as they have been due for a certain amount of time.

1**2****Notification bar - Finish due reminder and Finish overdue**

In the header bar, attached to the top right margin of the work area, it shows both the number of runs that are due for finishing and of those that are overdue.

- The left (orange) part of the bar indicates the number of runs of an event-triggered operation that needs to be executed and completed in due time.
- The right (red) part of the bar indicates the number of runs of an event-triggered operation that has not yet been completed but is overdue.

1**Notification bar - overdue**

In the header bar, attached to the top right margin of the work area, indicates the number of runs of an event-triggered operation that are overdue and need to be processed immediately.



Notification bar - paused unit procedure

In the header bar, attached to the top right margin of the work area, indicates that at least one unit procedure running at the current work center is currently paused.



OK

In pre-defined workflows, completes the last step of a workflow and finishes the workflow itself.



Open Cockpit (ALT+T)

In pre-defined workflows, switches to the Cockpit to display the lists of **Startable** and **Running** processes and workflows.



Page down (PGDN)

In the control bar, pages down in the content of the work area.



Page up (PGUP)

In the control bar, pages up in the content of the work area.



Running at current work center

In the **Appendable Workflows** Selection Window, controls if the list of appendable workflows (page 88) only shows the workflows that have been or are running at the current work center. By default, the button is selected, thus restricting the number of listed workflows.



Switch (F8)

In the control bar, switches between concurrently running operations.

Cockpit

The Cockpit provides the following buttons and icons:



Abort and reactivate unit procedure

On the Actions dialog (page 78) of a unit procedure, aborts and reactivates a unit procedure thus making it available for restarting.



Abort operation

On the Actions dialog (page 81) of an operation, aborts the operation.



Actions



Actions (active operation)

In the list of **Running** processes of the Cockpit (page 29), to the right of a unit procedure, an operation, or template, it opens the Actions dialog to display the additional actions that are available.



Append

On the Actions dialog (page 78) of a unit procedure, opens the **Appendable Workflows** Selection Window to filter and select a workflow available for appending to the unit procedure.



Cancel all runs

On the Actions dialog (page 82) of an event-triggered operation, cancels all runs of the operation that have not been started yet.



Cancel run

On the Actions dialog (page 82) of an event-triggered operation, cancels the run of the operation.

Only available if the run has not been started yet.



Cancel workflow

On the Actions dialog (page 78) of a workflow unit procedure, cancels the workflow. Only available if the workflow is not appended to an order.



Currently active operation

In the list of **Running** processes of the Cockpit (page 29), displays the identifier of the operation, the identifier of the device on which it is running, and the start time of its unit procedure.



Detach



Detach (active operation or unit procedure)

In the list of **Running** processes of the Cockpit (page 29), to the right of a detachable operation or unit procedure, detaches the operation (page 93) from the current station or the unit procedure (page 91) from the current work center, thus holding its processing and making it available for resuming at another station or work center, respectively.



Event-triggered operation run before

expiry

In the list of **Running** processes of the Cockpit (page 29), indented below its unit procedure, displays the identifier of the operation with the count of its run, the **X** marker indicating that the expiry of the run is the next stage in its escalation progress (page 89) with the time when this occurs, and the start time of its unit procedure.



Event-triggered operation run before

overdue


In the list of **Running** processes of the Cockpit (page 29), indented below its unit procedure, displays the identifier of the operation with the count of its run, the **!!!** marker indicating that an overdue marker is the next stage in its escalation progress (page 89) with the time when this occurs, and the start time of its unit procedure.



Event-triggered operation run before

reminder


In the list of **Running** processes of the Cockpit (page 29), indented below its unit procedure, displays the identifier of the operation with the count of its run, the ! marker indicating that a run due reminder is the next stage in its escalation progress (page 89) with the time when this occurs, and the start time of its unit procedure.

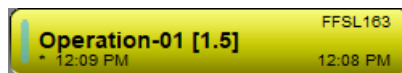

Event-triggered operation run before finish due

reminder

In the list of **Running** processes of the Cockpit (page 29), indented below its unit procedure, displays the identifier of the operation with the count of its run, the * marker indicating that a finish due reminder is the next stage in its escalation progress (page 89) with the time when this occurs, and the start time of its unit procedure.

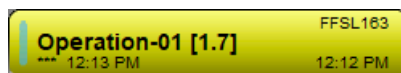

Event-triggered operation run before finish

overdue

In the list of **Running** processes of the Cockpit (page 29), indented below its unit procedure, displays the identifier of the operation with the count of its run, the *** marker indicating that a finish overdue is the next stage in its escalation progress (page 89) with the time when this occurs, and the start time of its unit procedure.



Exceptions

On the Action dialogs of a unit procedure (page 82) or an operation (page 81), opens the overview Exception Window to display all exceptions and comments that have so far been recorded for the operations and phases of the current run of the unit procedure or the phases of the current run of the operation, respectively.




New run

On the Actions dialog (page 82) of a template or an event-triggered operation, starts a new run of the operation.



Open in Viewer

On the Actions dialog of the unit procedure (page 78), operation (page 81) or event-triggered operation (page 82), it opens the selected component in Production Execution Viewer, provided the Production Execution Viewer application is already running in the background.

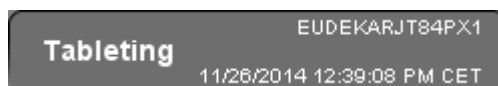

Operation running at different station

In the list of **Running** processes of the Cockpit (page 29), displays the identifier of the operation, the identifier of the device on which it is running, and the start time of its unit procedure.



Pause unit procedure

On the Actions dialog (page 78) of a unit procedure, pauses the execution of a unit procedure by informing all active phases of this fact. Trigger phases react by holding their triggering activities until the unit procedure is resumed.



Remove template

On the Actions dialog (page 82) of a template, removes the template from the list of running processes, thus ending the cycle of the event-triggered operation.



Resume unit procedure

On the Actions dialog (page 78) of a paused unit procedure, resumes the execution of the unit procedure by informing all active phases of this fact. Trigger phases react by resuming their triggering activities.



Run due reminder

In the list of **Running** processes of the Cockpit (page 29), to the left of an event-triggered operation's run button, indicates that the run needs to be processed as soon as possible as it have been due for a certain amount of time.



Run overdue

In the list of **Running** processes of the Cockpit (page 29), to the left of an event-triggered operation's run button, indicates that the run needs to be processed immediately.



Finish due reminder

In the list of **Running** processes of the Cockpit (page 29), to the left of an event-triggered operation's run button, indicates that the run needs to be finished as soon as possible as it has been due for completion for a defined amount of time.




Finish overdue

In the list of **Running** processes of the Cockpit (page 29), to the left of an event-triggered operation's run button, indicates that the run has not been completed within a defined amount of time and is declared overdue.



Running operation


In the list of **Running** processes of the Cockpit (page 29), indented below its unit procedure, displays the identifier of the operation, the identifier of the device on which it is running, and the start time of its unit procedure.

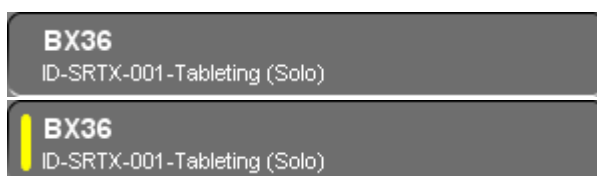

Running unit procedure of Batch Processing

In the list of **Running** processes of the Cockpit (page 29), displays the identifiers of the batch and of the unit procedure started from Batch Processing (page 105). Indented below the unit procedure button, the system lists those of its operations that are ready to be processed.

If the operator returns from processing, the button of the respective unit procedure is marked with a yellow indicator bar.

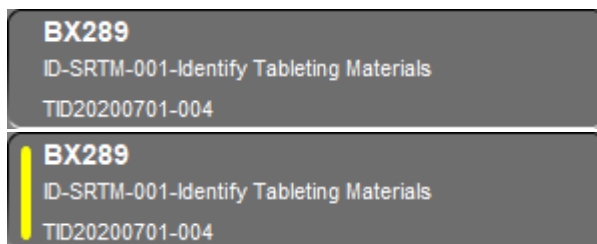


Running unit procedure of Batch Processing with treatment ID

In the list of **Running** processes of the Cockpit (page 29), displays the identifiers of the batch and of the unit procedure started from Batch Processing (page 105). As a third line the button shows the treatment ID. Indented below the unit procedure button, the system lists those of its operations that are ready to be processed.

If the operator returns from processing, the button of the respective unit procedure is marked with a yellow indicator bar.

TIP

Please note that the system can also be configured to display other information in the third line, such as a material description or a planning date.

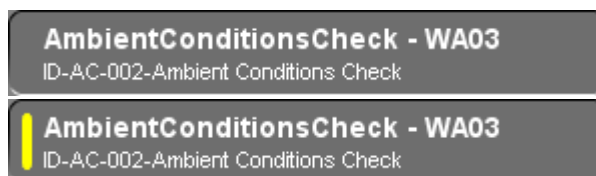

Running unit procedure of Workflow Processing

In the list of **Running** processes of the Cockpit (page 29), displays the name of the workflow, its identifier, and the identifier of the unit procedure started from Workflow Processing (page 109). Indented below the unit procedure button, the system lists those of its operations that are ready to be processed.


If the operator returns from processing, the button of the respective unit procedure is marked with a yellow indicator bar.



Running unit procedure of a one-click startable workflow

In the list of **Running** processes of the Cockpit (page 29), displays the name of the workflow, its identifier, and the identifier of the started unit procedure. Indented below the unit procedure button, the system lists those of its operations that are ready to be processed.

If the operator returns from processing, the button of the respective unit procedure is marked with a yellow indicator bar.

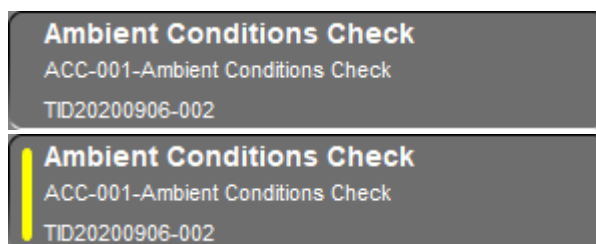


Running unit procedure of a one-click startable workflow with treatment ID

In the list of **Running** processes of the Cockpit (page 29), displays the name of the workflow, its identifier, and the identifier of the unit procedure started from Workflow Processing (page 109). As a third line the button shows the treatment ID. Indented below the unit procedure button, the system lists those of its operations that are ready to be processed.

If the operator returns from processing, the button of the respective unit procedure is marked with a yellow indicator bar.


TIP

Please note that the system can also be configured to display other information in the third line, such as a planning date.

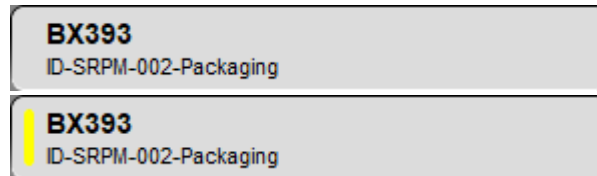**Running unit procedure of Batch Processing marked by Cockpit search**

In the list of **Running** processes of the Cockpit (page 29), displays the identifiers of the batch and of the unit procedure started from Batch Processing (page 105). If found as match of a Cockpit search, the display of the coloring is inverted and shows dark text on light background.

If the operator returns from processing, the button of the respective unit procedure is marked with a yellow indicator bar.

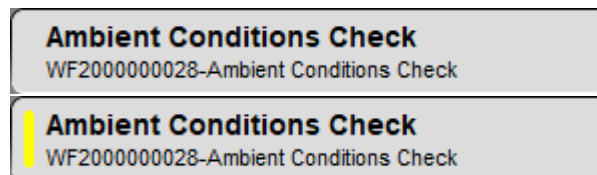**Running unit procedure of one-click startable workflow marked by Cockpit search**

In the list of **Running** processes of the Cockpit (page 29), displays the name of the workflow, its identifier, and the identifier of the unit procedure started from Workflow Processing (page 109). If found as match of a Cockpit search, the display of the coloring is inverted and shows dark text on light background.

**Sort alphabetically**

Above the list of **Running** processes of the Cockpit (page 29), toggles the sorting type of the running processes.

The icon displayed on the button indicates the current sorting type.

**Sort ascending**

Above the list of **Running** processes of the Cockpit (page 29), toggles the sort order of the running processes.

The icon displayed on the button indicates the current sort order.



Sort by timestamp

Above the list of **Running** processes of the Cockpit (page 29), toggles the sorting type of the running processes.

The icon displayed on the button indicates the current sorting type.



Sort by treatment ID

Above the list of **Running** processes of the Cockpit (page 29), toggles the sorting type of the running processes.

The icon displayed on the button indicates the current sorting type.

The sorting criterion is only available if your system is configured to support treatment ID-related features.



Sort descending

Above the list of **Running** processes of the Cockpit (page 29), toggles the sort order of the running processes.

The icon displayed on the button indicates the current sort order.



Startable batch orders

In the list of **Startable** processes of the Cockpit (page 29), initiates execution. It opens the list of startable batch order steps for selecting an order step to be processed.



Startable material for processing

In the list of **Startable** processes of the Cockpit (page 29), initiates execution. It opens the list of materials that are available for processing.



Startable operation

In the list of **Running** processes of the Cockpit (page 29), indented below its unit procedure, displays the identifier of the operation and the start time of its unit procedure.




Startable workflow

In the list of **Startable** processes of the Cockpit (page 29), initiates the execution of the selected one-click startable or pre-defined workflow.



Startable workflows

In the list of **Startable** processes of the Cockpit (page 29), initiates execution. It opens the list of startable workflow steps for selecting a workflow step to be processed.

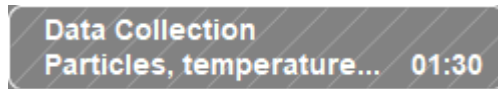

Template

In the list of **Running** processes of the Cockpit (page 29), indented below its unit procedure, displays information on the triggerable operations in two lines of text. The first line holds the identifier of the operation. The second line shows the description of the operation and in the bottom right the hours and minutes until the next operation will be triggered.

TIP

The countdown is only available for operations triggered by a time-based trigger. It shows hours and minutes while seconds are rounded down.

A template as such cannot be started for processing. Only when a run of the operation has been triggered, either by the system or by manual run creation, is the template button replaced by a tappable operation button.



Unit procedure paused

In the list of **Running** processes of the Cockpit (page 29), to the left of a paused unit procedure, indicates that the unit procedure is paused and some phases, such as trigger phases, may be on hold.

Notification Panel

The Notification panel provides the following buttons and icons:



Mark notification as read (CTRL+R)

In the control bar, marks selected notification as read.



Mark notification as unread (CTRL+R)

In the control bar, removes the read mark from the selected notification.




Pin notification (CTRL+P)

In the control bar, moves selected notification to the top of the list.



Remove notification (CTRL+DEL)

In the control bar, removes selected notification from the list.

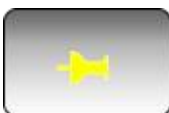

Unpin notification (CTRL+P)

In the control bar, moves selected notification from the top of the list to the bottom of the list.

Common

The following buttons and icons are common to several phases or workflows:



Alarm

For phases that visualize a status, indicates that an alarm status has occurred.



Bookmark

In the exception recording panel of the Exception Window, toggles the bookmark indicator for the exception selected in the list of exceptions.



Cancel order step

In pre-defined workflows, permanently terminates the selected order step.



Checkmark

For phases that visualize a status, indicates that the status is as intended.

**Enable**

For complex phases that require several actions to be completed before a follow-up action, such as confirming the phase, can be started, an **Enable** toggle button is available to explicitly enable the button that triggers the follow-up action.

**Error**

For phases that visualize a status, indicates that an error status has occurred.

**Exception at value**

In phases that display values as table, indicates that an exception has been recorded for the value given in the marked table cell.

**First page**

Phases that can display multi-page documents provide a document viewer with navigation buttons. The **First page** button moves to the beginning of the document.

**Identify batch manually**

In pre-defined workflows, allows the manual identification of a batch as fallback if the batch label cannot be scanned.

**Identify load carrier manually**

In pre-defined workflows, allows the manual identification of a load carrier as fallback if the load carrier label cannot be scanned.

**Identify subplot manually**

In pre-defined workflows, allows the manual identification of a subplot as fallback if the subplot label cannot be scanned.



Last page

Phases that can display multi-page documents provide a document viewer with navigation buttons. The **Last page** button moves to the end of the document.



Loading

For phases that load external data for display, such as a historical data chart, the animated icon indicates that loading is in progress.



Move all down

In pre-defined workflows, used for assignments or for detailed selection between two data tables. All data objects of the upper data table are moved to the lower data table.



Move all up

In pre-defined workflows, used for assignments or for detailed selection between two data tables. All data objects of the lower data table are moved to the upper data table.



Move selected down

In pre-defined workflows, used for assignments or for detailed selection between two data tables. The data objects selected in the upper data table are moved to the lower table.



Move selected up

In pre-defined workflows, used for assignments or for detailed selection between two data tables. The data objects selected in the lower data table are moved to the upper table.



Next page

Phases that can display multi-page documents provide a document viewer with navigation buttons. The **Next page** button moves to the beginning of the next page in the document.

**Phase activity ongoing**

For phases that are active for a period of time, such as monitoring phases, the animated icon indicates that the phase activity is ongoing.

**Previous page**

Phases that can display multi-page documents provide a document viewer with navigation buttons. The **Previous page** button moves to the beginning of the previous page in the document.

**Scanner ready - mandatory**

In pre-defined workflows, indicator icon that shows that the connected scanner is ready for use. The identification action is mandatory at this point.

**Scanner ready - optional**

In pre-defined workflows, indicator icon that shows that the connected scanner is ready for use. The identification action is optional at this point.

**Scroll down**

Phases that can display multi-page documents provide a document viewer with navigation buttons. The **Scroll down** button scrolls linewise down through the document.

**Scroll up**

Phases that can display multi-page documents provide a document viewer with navigation buttons. The **Scroll up** button scrolls linewise up through the document.

**Selected**

In pre-defined workflows, indicator icon used for multi-selection in data tables. The data objects are marked selected for further processing.



Skip step

In pre-defined workflows, skips the current workflow step.



Toggle unit of measure

For input values in phases that can deal with various units of measure, the selected unit is displayed in the center of the button. Clicking the button toggles to the next unit, which shows as preview in the top right corner of the button. The following units are available: **t, g, kg, lb, mg, oz.**



Unselected

In pre-defined workflows, indicator icon used for multi-selection in data tables. The data objects are no longer marked selected for further processing.

Specific

The following buttons and icons only occur during processing of specific workflows and phases:



Consume subplot

In the **Weigh** phase of a Dispense operation in Recipe Processing (page 105) or Material Processing (page 107), registers a subplot as consumed, thus removing it from all inventory lists.



Even split

In the **Sublot Split** workflow (page 127), provides a form to execute an even split on a subplot.



Manual scale indicator

In the **Select scale** phase of a Dispense or Weighing operation, on a scale button, indicates that the scale is configured as manual scale and that all values have to be entered manually.



New batch

In the **Material Receipt** workflow (page 111), provides a form to create a new batch.



Release scale

In the **Weigh** phase of a Dispense operation in Recipe Processing (page 105) or Material Processing (page 107), allows an operator to release the scale after an underweight condition during **Net removal** weighing.



Scale loaded indicator

In the **Select scale** phase of a Dispense or Weighing operation, on a scale button, indicates that the scale is loaded, typically with a container, and must not be zeroed.



Scan action indicator for container

In the **Identify container** phase of a Dispense or Weighing operation, in Recipe Processing (page 105) or Material Processing (page 107), indicates that the scan target for phase confirmation is a container.



Scan action indicator for printer

In the **Print report** phase of a Dispense or Weighing operation, in Recipe Processing (page 105) or Material Processing (page 107), indicates that the scan target for phase confirmation is a printer.



Scan action indicator for scale

In **Tare** and **Release scale** phases of a Dispense or Weighing operation, in Recipe Processing (page 105) or Material Processing (page 107), indicates that the scan target for phase confirmation is a scale.



Scan action indicator for scale with value

In the **Weigh** phase of a Dispense or Weighing operation, in Recipe Processing (page 105) or Material Processing (page 107), indicates that the scan target for phase confirmation is a scale.



Scan action indicator for subplot

In the **Identify material** phase of a Dispense or Weighing operation, in Recipe Processing (page 105) or Material Processing (page 107), indicates that the scan target for phase confirmation is a subplot.



Scan action indicator for undefined scale

In the **Select scale** phase of a Dispense or Weighing operation, in Recipe Processing (page 105) or Material Processing (page 107), indicates that the scan target for phase confirmation is a scale.



Select batch

In the **Material Receipt** workflow (page 111), provides a form to select a batch from the list of existing batches.



Uneven split

In the **Sublot Split** workflow (page 127), provides a form to execute an uneven split on a subplot.

Cockpit

The Cockpit represents an operator's entry point to the execution process. Specific to each station and user, the Cockpit provides a list of all startable, resumable, and running processes, from which the operator can select unit procedures and operations to start or switch to. The running processes also include those operations of the same unit procedure that are running at another station.

In its header bar, the Cockpit displays two lines of text. The first line shows the surname and login name of the logged-in user, the second line shows the identifiers of the current work center and station.

TIP

If you wish to process Dispense operations at your work center, make sure it is assigned to a room of the **Room (RS)** equipment type. Work centers and room equipment entities are maintained in Data Manager.

In the left part of its work area, the Cockpit displays buttons for the **Startable** processes and workflows, sorted by application area in the following sections:

TIPS

Please note that the system can be configured to display additional sections, such as **Quality** or **Maintenance** or to show only selected application sections, processes, and workflows at specific work centers.

Which processes are displayed also depends on the access rights of the logged-in user and the infrastructure of your system. If PharmaSuite is configured to communicate with Warehouse Management, the **Relocation** inventory workflow (page 111) is not available.

Processes that are based on confidential master recipes or master workflows to which the user does not have the required access rights remain invisible.

- Production (page 105)
- Inventory (page 111)
- System
 - **Help** system (page 102)
 - **About** function

It opens the **About PharmaSuite** dialog to display system-related information, such as the current system version and build, the logged-in user, work center, and database-related information.

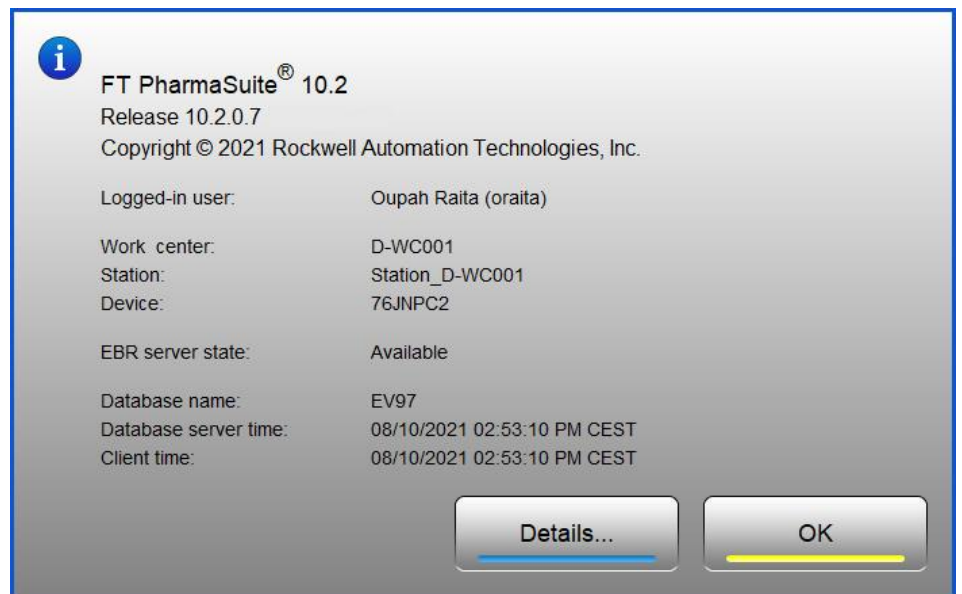


Figure 8: About dialog

TIP

Please note that the dialog also indicates the **EBR server state**. Only when the EBR server is available can you process recipes or EBR workflows with PharmaSuite for Production Execution.

Click the **Details** button to view more specific technical information on the system and its environment. From the **PharmaSuite Installation Details** dialog, you can copy the listed detail data to the clipboard.

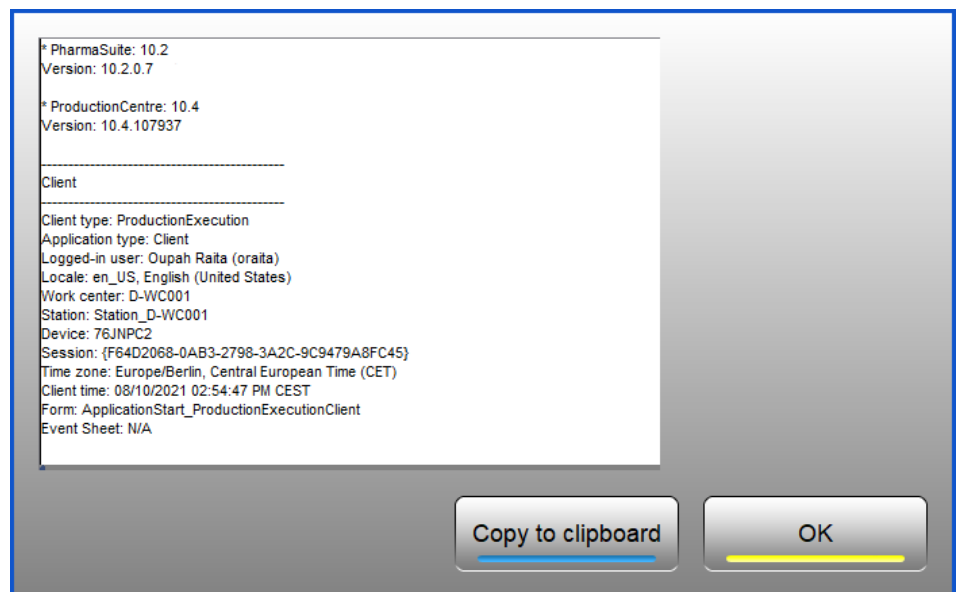


Figure 9: About dialog - Details

TIP

Please note that the path to the PharmaSuite log files is given in the last section of the listed detail data.

■ **Register at Station** function

It is primarily intended for mobile devices and allows you to register your device at another station. By doing so, the system automatically unregisters your device from where it was previously registered.

On the **Register at Station** dialog you can filter the number of station identifiers displayed in the **Station** option list. The filter becomes effective with the second character you type.

By default, the search is not case-sensitive and the search string may occur anywhere in the station identifiers. Click the magnifying glass icon to set the search criteria to be case-sensitive, allow the use of wildcards or regular expressions, or match from start or exactly.

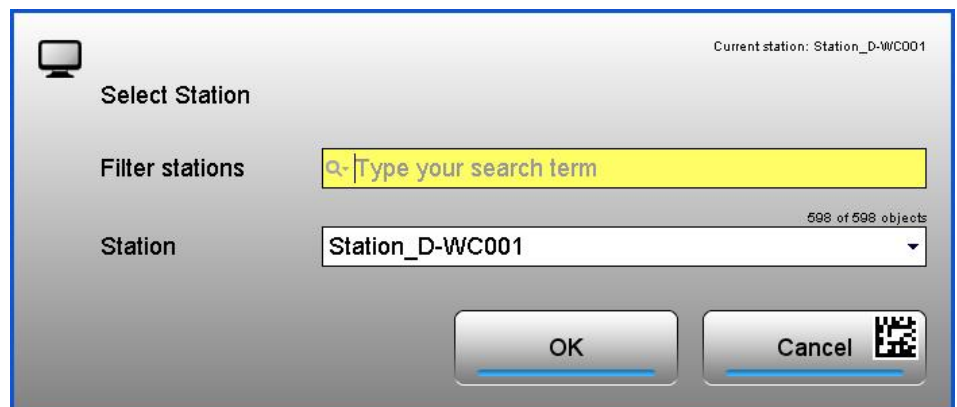


Figure 10: Register at Station dialog

TIPS


Please note that the system supports registration of a device at a station by scanning the station's barcode.

The system rejects a user's registration at the station if there are confidential processes running at the station to which the user does not have the required access rights.

Please note that a station change can also be triggered from an external system.

■ **Change User** function

It allows you to change to another user without having to exit the system, for example at shift change. The **Change User** dialog does not have a **Cancel** button and can only be completed with a valid login name and password. If the function was called unintentionally, it is possible to log in again by typing the login name and password of the currently logged-in user. Thus, you can also use the dialog to lock the screen and prevent unauthorized access to the system, while order or workflow processing continues in the background. For this reason, the **Change User** function also serves to provide a configurable auto-lock mechanism that detects if the system is not in active use by an operator and consequently locks the screen after a configurable period of inactivity, typically 15 minutes. Any kind of user interaction with the system by touch, keyboard, scanner, or mouse prevents the system from auto-locking.

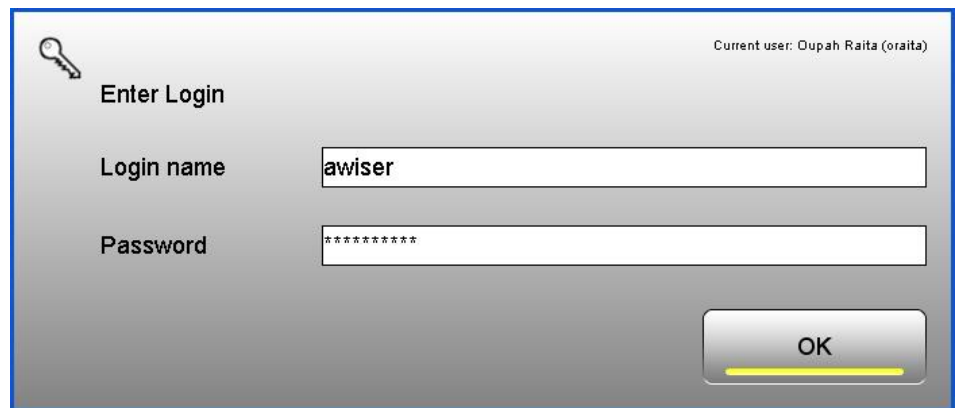

The screenshot shows a 'Change User' dialog box with a key icon in the top left. The title is 'Enter Login'. In the top right corner, it says 'Current user: Dupah Raita (oraita)'. There are two input fields: 'Login name' with the text 'awiser' and 'Password' with masked characters '*****'. An 'OK' button is located at the bottom right of the dialog.

Figure 11: Change User dialog

TIP

Please note that when a user change is initiated while an order or workflow is running, the new user needs to have the access rights required for processing the running order or workflow.

■ **Device Configuration** function

It allows you to extend the vertical screen size of PharmaSuite for Production Execution. The **Configure Scree Size** dialog provides a vertical slider that extends the height of the application window from its default size up to the full height of your current screen (minus the height of the taskbar, if visible). The size graphic to the right of the slider indicates the selected screen size as yellow inner frame, relative to the blue outer frame that represents the available screen size of your current device.

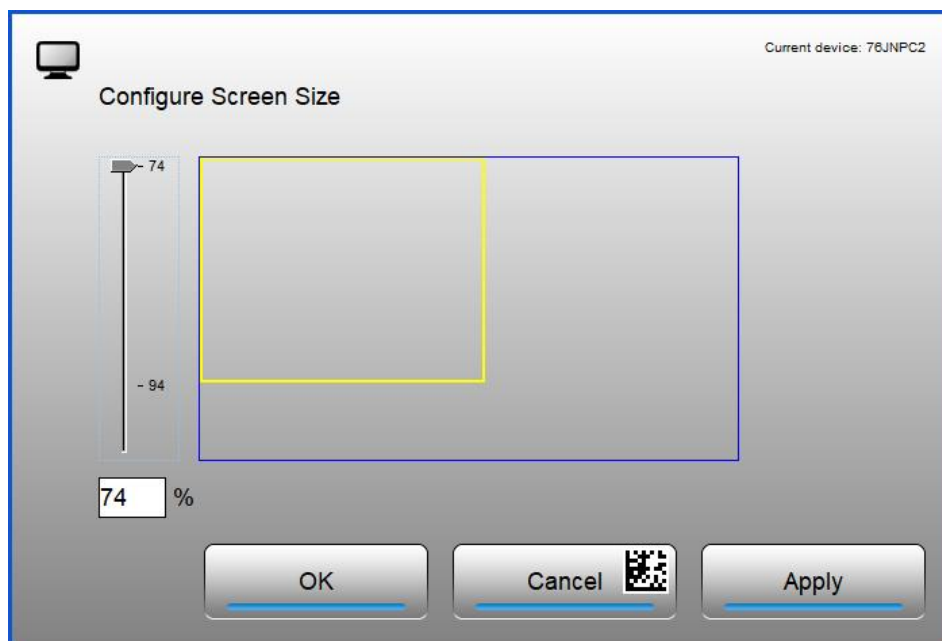


Figure 12: Configure Screen Size dialog

■ **Exit** function (page 60)

When you return to the Cockpit while one or more processes are running, the list of startable processes and workflows displays started workflows that cannot be started multiple times as grayed out. Only **Batch Processing**, **Material Processing**, and **Workflow Processing** can be started multiple times.

TIP

Please note that instead of selecting a batch operation from the Cockpit, you can also scan the barcode of its batch or its treatment ID. A batch's barcode is typically located on the cover page of its batch report.

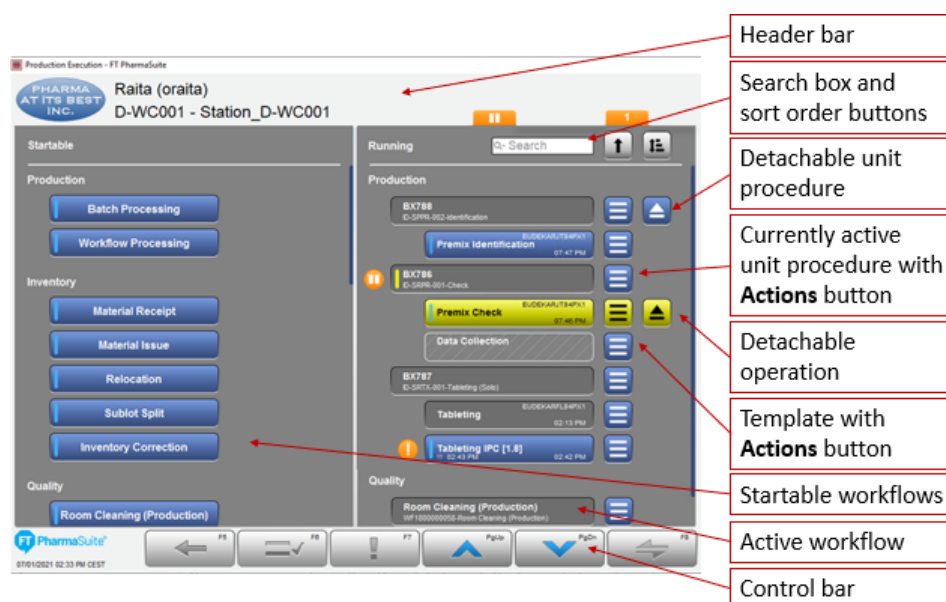


Figure 13: Cockpit

The list of **Running** processes and workflows is organized by application area, in the same way as the list of startable processes and workflows.

Above the list, the system provides a search box for marking all unit procedures that contain the given search term anywhere in the text displayed on the unit procedure buttons of the cockpit. The search becomes effective when you press the ENTER key. The system marks the matching unit procedures by inverting their coloring to dark text on light background and scrolls to show the first matching unit procedure at the top of the list of running unit procedures.

To the right of the search box, the system provides buttons to change the sort order of the listed unit procedures. The list can be sorted in ascending or descending order, alphabetically, by timestamp, or by treatment ID if the system is configured to support treatment ID-related features. The icon on the respective button indicates the current sort order.

TIP

Please note that the unit procedures are sorted within their respective application area sections.

- Alphabetically ascending displays the unit procedures sorted by their identifiers from A to Z from top to bottom.
- Alphabetically descending displays the unit procedures sorted by their identifiers from Z to A from top to bottom.
- Ascending by start timestamp displays the unit procedures sorted earliest to latest from top to bottom.
- Descending by start timestamp displays the unit procedures sorted latest to earliest from top to bottom.
- Ascending by treatment ID displays the unit procedures sorted alphanumerically from A to Z or 1 to the highest number, respectively, from top to bottom.
- Descending by treatment ID displays the unit procedures sorted alphanumerically from Z to A or the lowest number to 1, respectively, from top to bottom.

The system displays a running process with its unit procedure as gray header button and its startable or running operations as indented tappable buttons.

Unit procedure buttons in Batch Processing show the batch they are producing and their order step identifier with the unit procedure's count if it has been reactivated.

Unit procedure buttons of workflows show the unit procedure's identifier.

If a treatment ID is defined for an order or workflow, its unit procedure button also displays the treatment ID.

TIP

Please note that the system can also be configured to display other information in the third line, such as a material description or a planning date.

When you return to the Cockpit from processing, the button of the unit procedure from which you have come is marked with a yellow indicator bar. The list of running unit procedures scrolls to show the marked unit procedure automatically at the top.

Unit procedures of Batch Processing can provide a button to access further actions (page 78), such as pausing or resuming it (page 84), viewing and processing its exceptions (page 47), aborting and reactivating it (page 92), appending workflows (page 88) (only for orders), or canceling (page 78) (only for workflows).

Some unit procedures can be detached from their current work center (page 91) for being continued at another work center. A detachable unit procedure displays a **Detach** button (page 14) to the right of its unit procedure button.

Operations of the same unit procedure that are already running at another station and are thus unavailable at the current station, are disabled and gray. The currently active operation, from which you have opened the Cockpit, is represented as yellow button, while other startable or running operations have blue buttons.

Event-triggered operations (page 82) first show as template graphics (page 21) and later as tappable operation buttons (page 14) once a run has been created. Event-triggered operations can enforce their timely execution by implementing an escalation scenario (page 89), whose stages are indicated by the notification tabs in the header bar (page 3) and escalation markers to the left of the template or runs of the event-triggered operation.

TIPS

Please note that phases and operations of looped processes can be run and completed several times. In these cases, the system appends a count to their identifiers that indicates how often they have been run.

Event-triggered operations display a two-digit count where the first digit represents the number of times the operation has been run and completed, while the second digit represents the count of runs created as long as its template was active.

For some of the listed running operations, the system provides additional actions. They display an **Actions** button (page 13) to their right, which opens their Actions dialog (page 82) to access the available actions.

Once you have started an operation on a device, it is attached to the device as well as to the station at which the device is registered. A running operation shows its activity by pulsing its indicator bar. Additionally, it displays the identifier of the device on which it is being run in small font in the upper right corner of its operation button.

There are running operations that can be detached (page 93) from their current device and station for being continued at another station. A detachable operation displays a **Detach** button (page 14) to the right of its operation button.

TIP

The list of **Running** processes also contains all operations that have been started in a previous session, have been closed before completion, and can be resumed. This can be especially important if, for example, a power failure causes the system to shut down unexpectedly.

When you process an uninterruptible operation such as a Dispense operation and experience an unexpected shutdown, you must resume and process this operation first after restart. As indicator, the system displays the affected operation as currently active, with its operation button in yellow.

- To page through the listed processes and workflows, tap the **Page up/Page down** buttons in the control bar.
- To start any of the startable processes or workflows, tap it.
For processing of recipes, materials, and EBR workflows, the system displays the started unit procedure as running process with its list of startable operations. For pre-defined workflows, the system directly switches to the Execution Window (page 42) for processing the workflow.
- To switch to the Execution Window (page 42) and process the phases of an operation, tap the respective blue operation button.
- To return to the Execution Window (page 42) and proceed with processing the phases of an operation, tap the respective yellow operation button or the **Back** button in the control bar.
- To detach a running unit procedure in order to make it available for being resumed at another work center, tap its **Detach** button (page 14).
- To detach a running operation in order to make it available for being resumed at another device and station, tap its **Detach** button (page 14).
- To access additional actions that are available for an item listed as **Running**, tap its **Actions** button (page 13).
- To return to the form from which you have accessed the Cockpit, tap the **Back** button in the control bar.
- To access the Navigator (page 38) with the execution history, tap the **Navigator** button in the control bar.

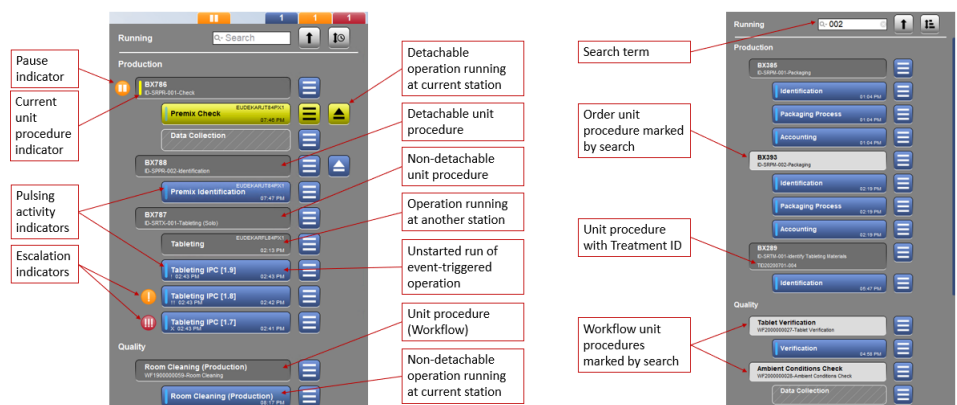


Figure 14: Running processes

Navigator

The Navigator represents an operator's information and service point during the execution process. Specific to each unit procedure, the Navigator provides a historical view of the executed process steps. For each step, the operator can access basic and detail information as well as additional actions.

In its header bar, the Navigator displays two lines of text. The first line shows the surname and login name of the logged-in user, the second line shows the identifiers of the current work center and station.

Event-triggered operations can enforce their timely execution by implementing an escalation scenario (page 89), whose stages are indicated by the notification tabs in the header bar (page 3).

The work area of the Navigator displays the data panel at the top and three columns below it:

- For order-related execution, the data panel displays two data frames. The left frame holds the identifiers of the order, the produced batch, and the unit procedure with its count if it has been reactivated. The right frame shows the identifier and short description of the produced material.
For workflow execution, only the left data frame is filled. It holds the identifier of the workflow, its processing name, and the identifier of the unit procedure.
- The left column (phase column) lists all phases of the unit procedure that have been executed so far, ending with the phase or phases of the currently active operation. A phase displays the identifier of the operation to which the phase belongs as well as its phase identifier. Phases that have been fully executed and confirmed are shown with a gray background and can no longer be directly accessed for processing. Phases that have been processed at another station have a dashed phase frame. The currently active phase or phases are represented as tappable, yellow buttons.

TIP

Please note that phases and operations of looped processes can be run and completed several times. In these cases, the system appends a count to their identifiers that indicates how often they have been run.

Event-triggered operations display a two-digit count where the first digit represents the number of times the operation has been run and completed, while the second digit represents the count of runs created as long as its template was active.

If a process includes phases that are skipped on the user interface, they are nevertheless included in the count, since the system executes them in the background as hidden phases.

- The center column (information column) provides buttons to access the Detail Information Window (page 50) for the individual phases that have been fully executed and confirmed. Detail information is provided by a phase-specific report. From here you can also see all exceptions and comments (page 63) that have been recorded for the phase so far and add new exceptions or comments, if necessary.

TIP

Phases for which an exception has been recorded are displayed with a red upper left corner as exception marker.

- The right column (action column) provides buttons to access up to three actions. For active phases the **Abort** and **Repair** buttons are available to resolve phase definition-related issues (page 76).
For phases that have been fully executed and confirmed, there may be exception situations that require a post completion correction of a value or reprinting of a label are available. Tapping the respective action button will thus take you to the Exception Window (page 44) where you can execute and confirm the action in the same way as when you record a user-triggered, phase-specific exception (page 70).
Phases, on which you have performed an action, display the red exception maker.
- To page through the listed phases, tap the **Page up/Page down** buttons in the control bar.
- To switch to the Execution Window (page 42) for processing a phase, tap the respective yellow phase button in the phase column.
- To switch to the Detail Information Window (page 50) of a phase, tap the respective button in the information column.

- To return to the form from which you have accessed the Navigator, tap the **Back** button in the control bar.
- To open the Cockpit (page 29) with access to other processes, workflows, or system functions, tap the **Cockpit** button.

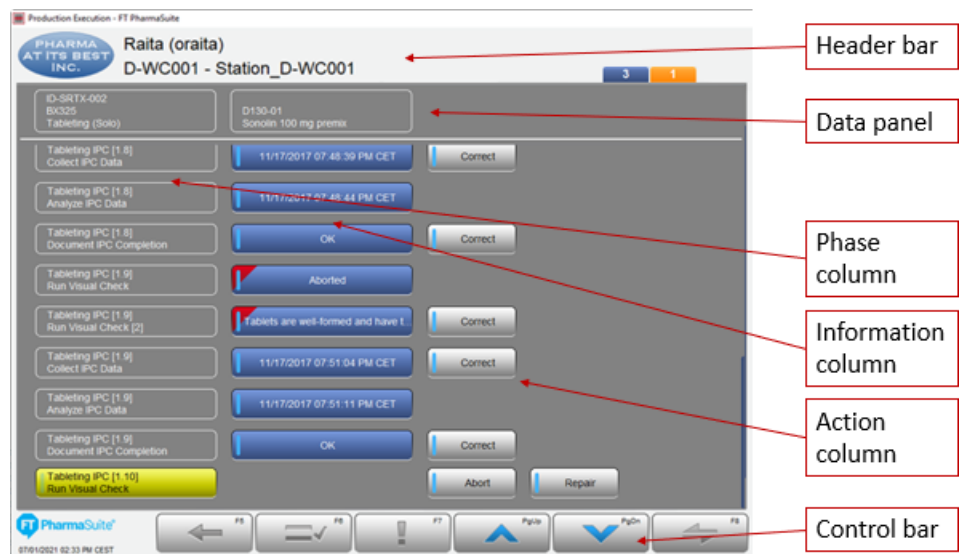


Figure 15: Navigator



Figure 16: Data panel of a non-order-related workflow

Selection Window

A Selection Window supports an operator whenever the process requires him to choose one or several objects, such as order steps for starting (page 107) or workflows for appending (page 86).

In its header bar, a Selection Window displays a selection-specific title and the identifiers of affected objects, if applicable.

Event-triggered operations can enforce their timely execution by implementing an escalation scenario (page 89), whose stages are indicated by the notification tabs in the header bar (page 3)

The work area contains the list of objects that are available for selection. For lists that can potentially have a large number of items the system provides filter tools to restrict the number of listed items. The filter becomes effective with the second character you type. By default, the search applies to all table columns, is not case-sensitive, and the search string may occur anywhere within the content of a cell. Click the magnifying glass icon to set the search criteria to search only in specific columns, be case-sensitive, allow the use of wildcards or regular expressions, or match from start or exactly.

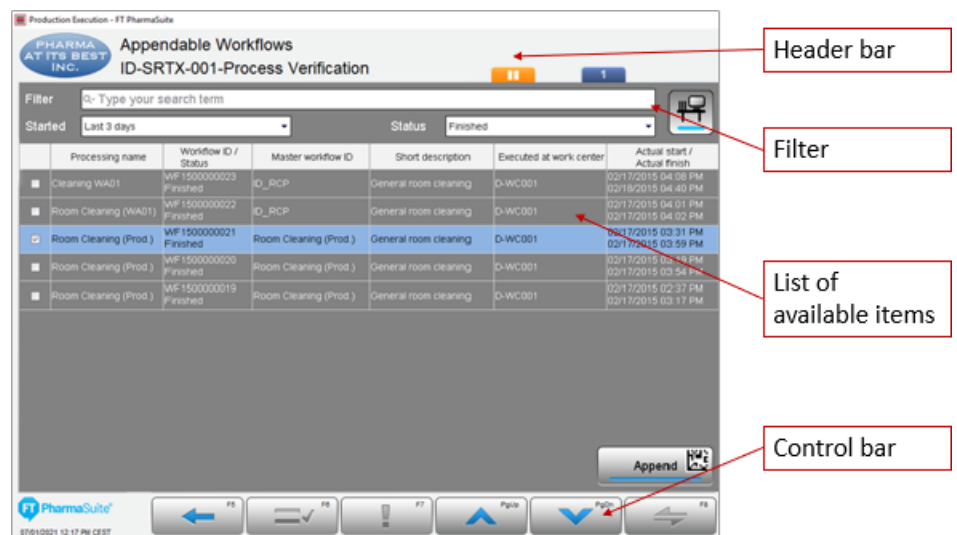


Figure 17: Selection Window

Execution Window

The Execution Window represents an operator's workspace during the execution process. Here, the operator performs the execution-related actions and registers process-relevant data.

In its header bar, the Execution Window displays in the first line for order steps, the identifier of the current order step with its unit procedure count if it has been reactivated, followed by the product batch. In the first line for workflows, it displays the identifier of the current workflow with its processing name. The second line holds the identifier of the operation with its count if it has been run more than once.

Event-triggered operations can enforce their timely execution by implementing an escalation scenario (page 89), whose stages are indicated by the notification tabs in the header bar (page 3).

The work area represents the dynamic part of the Execution Window and reflects the execution progress:

1. Typically, the first operation you perform in processing with EBR is to select what you wish to process. This can either be a specific order step, in **Batch Processing**, a workflow step in **Workflow Processing**, or one or several materials, in **Material Processing**.
2. After you have selected the order step or materials for processing the system returns you to the Cockpit (page 29) where all executable operations are listed as **Running**. Now you can tap to start the operation you wish to process.
3. The system displays the phases of the operation in their order of processing:
 - Previously executed phases are in read-only mode, with a checkmark on their **Confirm** buttons to indicate that they have been completed.
 - Currently active phases are accessible for processing.
Some phases may require an electronic signature before they can be completed.
 - Future phases that will only become active after one or all currently active phases have been executed are in read-only mode.
In a strictly sequential operation, the system displays all phases. If the operation includes one or more decision points, however, it will only show the future phases up to the next decision point.

Exceptions (page 44) that occur during processing can be recorded at all times.

TIP

Please note that some phases may rely on information or input that is not yet available when the preview is created. The preview of these phases only shows a corresponding system message.

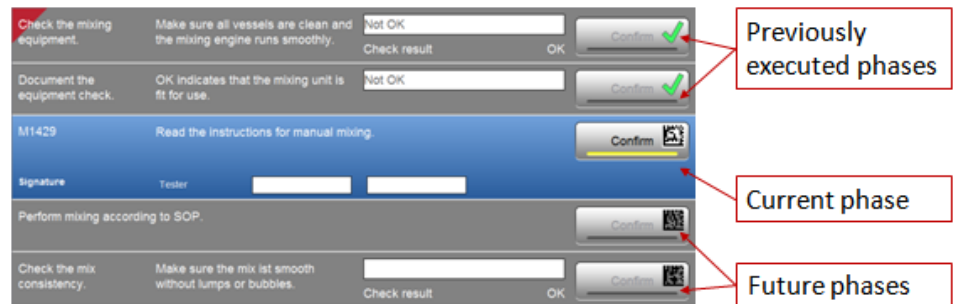


Figure 18: Phase representation during execution



Figure 19: Phase with required completion signature

4. Execute an accessible phase and tap its **Confirm** button.
 5. Move on to the next accessible phase until you have all phases processed and confirmed.
- When you have tapped the **Confirm** button of the last phase, the system returns you to the Cockpit for selecting a new operation or workflow.

- To page through the listed phases, tap the **Page up/Page down** buttons in the control bar.
- To access the Navigator (page 38) with the execution history, tap the **Navigator** button in the control bar.
- To open the Cockpit (page 29) with access to other processes, workflows, or system functions, tap the **Cockpit** button.

Exception Window

The Exception Window represents an exception handling (page 63) point where an operator can record, comment, and bookmark all exceptions that may occur during execution on the shop floor.

The system provides three variants of the Exception Window:

- the phase-specific Exception Window (page 44), which you can access during phase processing by tapping the **Exception** (F7) button (page 10) located in the control bar.
- the operation-specific Exception Window (page 46) that lists all exceptions and comments that have been recorded for the phases of the operation. You can access the operation-specific Exception Window from the Cockpit:
 1. Tap the **Actions** button (page 13) to the right of an operation.
The system opens the Actions dialog of the operation.
 2. Select the **Exceptions** button (page 15).
- the overview Exception Window (page 47) is only available for unit procedures started as order steps with recipe processing. It lists all exceptions that have been recorded for the phases and operations of the unit procedure. You can access the overview of all exceptions from the Cockpit:
 1. Tap the **Actions** button (page 13) to the right of a unit procedure.
The system opens the Actions dialog of the unit procedure.
 2. Select the **Exceptions** button (page 15).

Phase-specific Exception Window

In its header bar, the phase-specific Exception Window has a red triangle in the upper left corner as exception marker and displays in the first line, the identifier of the current order step or workflow with its unit procedure count if it has been reactivated, for order steps followed by the product batch, and in the second line, the identifier of the operation with its count if it has been run more than once and the phase identifier.

Event-triggered operations can enforce their timely execution by implementing an escalation scenario (page 89), whose stages are indicated by the notification tabs in the header bar (page 3).

The work area represents the dynamic part of the Exception Window and lists all exceptions and comments that have so far been recorded for the phase. Below the list, the system displays the exception recording panel that provides all controls to record and bookmark exceptions and add comments.

TIP

If the phase provides phase-specific exception (page 70) options, they are displayed above the list of recorded exceptions and the exception recording panel.

- To page through the listed exceptions and comments, tap the **Page up/Page down** buttons in the control bar.
- To return to the form from which you have accessed the Exception Window, tap the **Back** button in the control bar.

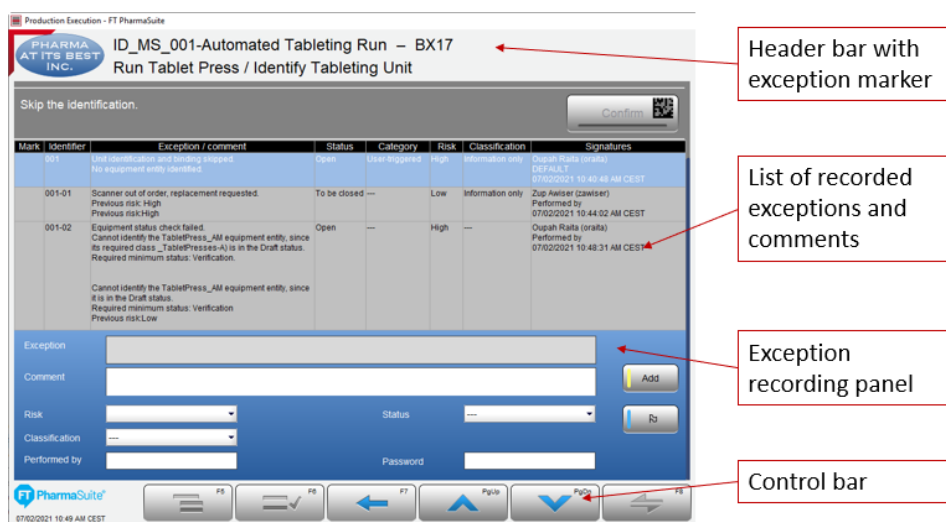


Figure 20: Phase-specific Exception Window

Operation-specific Exception Window

In its header bar, the operation-specific Exception Window has a red triangle in the upper left corner as exception marker and displays in the first line, the identifier of the current order step or workflow, for order steps followed by the product batch, and in the second line, the identifier of the unit procedure from which it has been accessed with its count if it has been reactivated.

Event-triggered operations can enforce their timely execution by implementing an escalation scenario (page 89), whose stages are indicated by the notification tabs in the header bar (page 3).

The work area represents the dynamic part of the Exception Window and lists all exceptions and comments that have so far been recorded for the phases of the current run of the operation. Below the list, the system displays the exception recording panel that provides all controls to record and bookmark exceptions and add comments.

- To page through the listed exceptions and comments, tap the **Page up/Page down** buttons in the control bar.
- To return to the Cockpit, tap the **Back** button in the control bar.

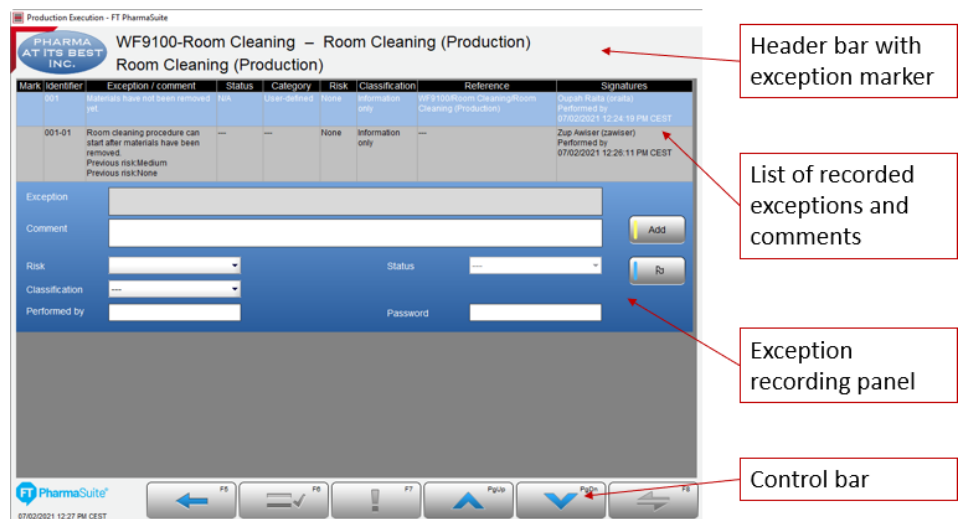


Figure 21: Operation-specific Exception Window

Overview Exception Window

In its header bar, the overview Exception Window has a red triangle in the upper left corner as exception marker and displays in the first line, the identifier of the current order step or workflow, for order steps followed by the product batch, and in the second line, the identifier of the unit procedure from which it has been accessed with its count if it has been reactivated.

Event-triggered operations can enforce their timely execution by implementing an escalation scenario (page 89), whose stages are indicated by the notification tabs in the header bar (page 3).

The work area represents the dynamic part of the Exception Window and lists all exceptions and comments that have so far been recorded for the operations and phases of the current run of the unit procedure. Below the list, the system displays the exception recording panel that provides all controls to record and bookmark exceptions and add comments.

TIP

Please note that the list of exceptions does not include the exceptions of workflows that are appended to the unit procedure.

- To page through the listed exceptions and comments, tap the **Page up/Page down** buttons in the control bar.
- To return to the Cockpit, tap the **Back** button in the control bar.

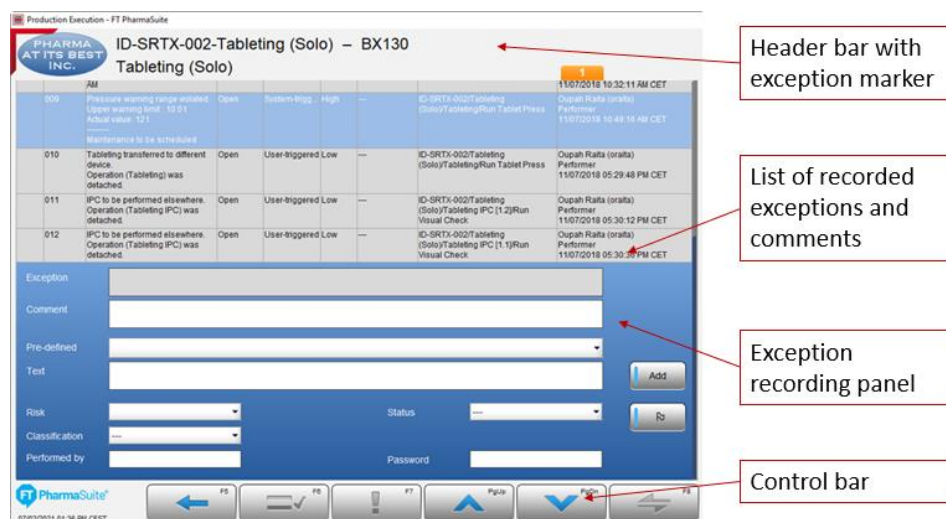


Figure 22: Overview Exception Window

Notification Panel

The Notification Panel lists station-specific notification messages and alarms available for the operator, on events, or actions that require the operator's attention.

TIP

Any changes and actions taken on the listed notifications for unit procedure or operation will only be displayed on the current station.

Notifications contain the following information:

- Order context (path to the relevant unit procedure or operation)
- A message text
- A timestamp indicating the creation time of the notification (determines the sort order of the notifications)
- A value as
 - a timestamp or
 - a duration as countdown timer, calculated by subtracting the current timestamp from the provided timestamp.
The timer counts down in minute increments.
The seconds are rounded up to next minute.
The timer stops when it reaches 0 minutes.
- An alarm (if triggered), that appears on the header bar (page 3)

The system allows to manage available notification via the toolbar buttons:

- Toggle button to:
 - Pin the selected notification (page 22)
A pinned notification moves on the top of the list.
 - Unpin the selected notification (page 22)
An unpinned notification moves down from the list.
- Toggle button to:
 - Mark the selected notification as read (page 21)
Selected notification is grayed out.
 - Mark the selected notification as unread (page 21)
Selected notification is back available.

- Remove the selected notification from the list (page 22)
Selected notification is removed from the list on the current station.

All notifications are scrollable. In order to perform an action, a notification must be selected. Only one notification can be selected at a time. A selected notification is visually distinguished by a yellow frame around its table row.

Notifications can be sorted by:

- Pin notification button,
Selected notification is moved further up in the list.
- Unpin notification button,
Selected notification is moved further down the list
- Creating new or updating existing notification,
Newer notification are further up in the list.

TIP

Notification text messages that have been updated or created for a preceding stage in the escalation scenario are overwritten only if an alarm/notification text has been configured for that stage in the escalation scenario.

Important notifications can be defined to be alarming notifications. As long as the notification panel displays unread notifications with an alarm indicator, the Notification Panel indicates this by displaying the Header Bar above it with a red background and an alarm text with the count of active alarms.

TIPS

Please note that, if an alarming notification is modified (text or value), it activates the alarm panel again, even if the original notification has already been marked as read.

Once the order or workflow to which a notification refers is finished, all corresponding notifications are removed from the Notification Panel.

Detail Information Window

The Detail Information Window represents an information point where an operator or reviewer can view phase-specific data of completed phases and comment all exceptions that have been recorded during execution. If necessary, he can also add further comments and exceptions.

In its header bar, the Detail Information Window displays in the first line for order steps, the identifier of the current order step with its unit procedure count if it has been reactivated, followed by the product batch. In the first line for workflows, it displays the identifier of the current workflow with its processing name. The second line holds the identifier of the operation, with its count if it has been run more than once, followed by the identifier of the selected phase.

Event-triggered operations can enforce their timely execution by implementing an escalation scenario (page 89), whose stages are indicated by the notification tabs in the header bar (page 3).

The work area represents the dynamic part of the Detail Information Window. In its upper section, it displays the phase-specific part of the batch report with detailed phase information, such as

- the batch and order identifier to which the phase belongs
- the material that will be produced with the order
- the start and finish timestamps of the phase.

In its lower section, the work area holds the list of exceptions and comments that have been recorded for the phase and the exception recording panel.

- To page through the listed information, tap the **Page up/Page down** buttons in the control bar.
- To return to the form from which you have accessed the Detail Information Window, tap the **Back** button in the control bar.

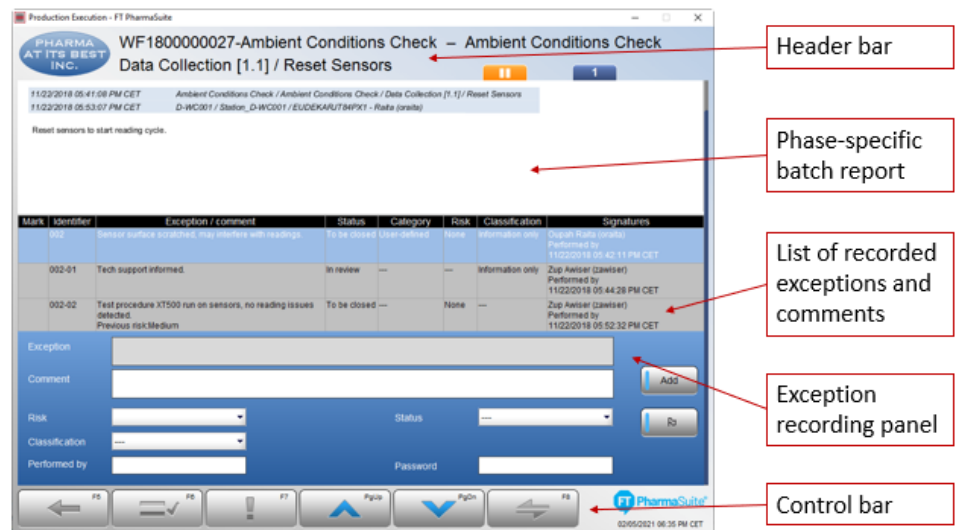


Figure 23: Detail Information Window

Input Devices

Production Execution with EBR is designed for use with various input devices to support the most suitable equipment for all expected application environments.

Keyboard

A physical keyboard may be available at your work station for typing data and values in input fields. The system also supports keyboard navigation and control with keyboard shortcuts.

Depending on the configuration of your system and environment, a virtual keyboard may be available instead of a physical one.

FRAMEWORK SHORTCUTS

Use the following keyboard shortcuts to trigger the functions of the control bar.

- ALT+F4
Closes a modal window or dialog, such as the help window, that may be open on top of the application.
If there is no other window or dialog open, it logs the current user out and closes the application.
- ALT+T
In the non-EBR framework, toggles the display of the Cockpit.
- CTRL+F4
Cancels the currently active workflow.
- F1
Opens a web browser to display the context-sensitive help (page [102](#)) of PharmaSuite for Production Execution.
- F5
In the control bar, toggles the display of the Cockpit (page [29](#)).
- F6
In the control bar, toggles the display of the Navigator (page [38](#)).
- F7
In the control bar, toggles the display of the phase-specific Exception Window (page [44](#)).
- F8
In the control bar, switches between concurrent operations of a running unit procedure.
- PAGE DOWN
In the control bar, pages the content of the work area down.
- PAGE UP
In the control bar, pages the content of the work area up.

NAVIGATION BAR KEYBOARD SHORTCUTS

In the navigation bar of pre-defined workflows, there are no specific assignments of keyboard shortcuts to the individual buttons located on the bar. Instead, the eight available button spaces are mapped to the function key blocks from F5 through F12. Press F5 through F8 to access the four button spaces to left of the center element in the navigation bar or press F9 through F12 to access the four spaces to the right of the center element.

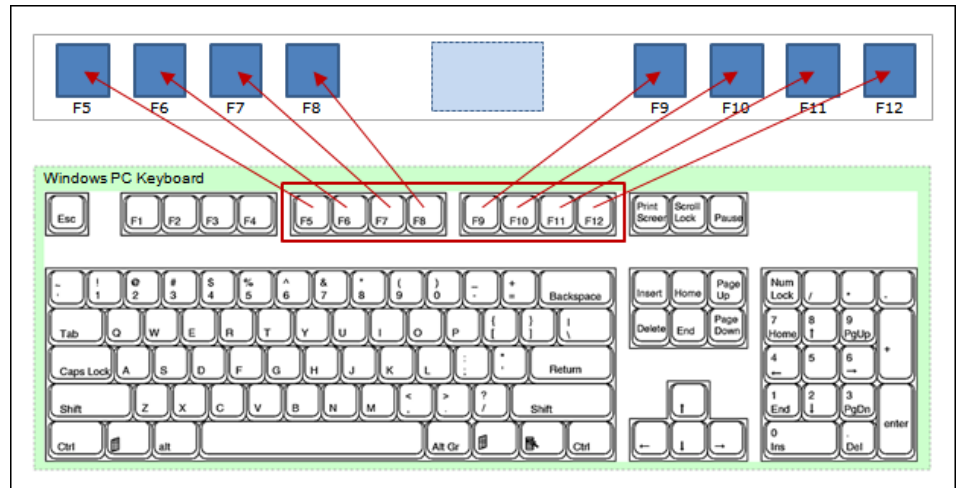


Figure 24: Keyboard mapping of navigation bar buttons

COCKPIT NAVIGATION

Use the following keys to navigate the Cockpit.

- TAB
Navigates along the defined navigation path, from top to bottom of the **Startable Workflows** list and then of the **Running Workflows** list.
- SHIFT+TAB
Navigates back along the defined navigation path from bottom to top of the lists of **Startable** or **Running** processes or workflows.

EXECUTION WINDOW NAVIGATION

Use the following keys to navigate the Execution Window.

- CTRL+TAB
For multi-line text boxes, navigates to the next control in the work area, along the defined navigation path.
- ENTER
For tables, moves to the next row and marks it for multi-selection, if applicable.
For option lists, selects the highlighted option.
- F9
Starts processing in the **Material Processing** workflow.
- SHIFT+TAB
Navigates back along the defined navigation path from control to control in the work area.
- SPACEBAR
For tables, marks the highlighted row for multi-selection.
For option lists, selects the highlighted option.
- TAB
Navigates along the defined navigation path, Typically, from the top-left corner to the bottom-right corner of the form, from control to control in the work area.
- UP/DOWN ARROW
Navigates between the rows of a table or the options of an option list.

DIALOG BOXES

The system displays information messages, warnings, error messages, and some types of signature requests as modal dialog boxes on top of the application. The following keys are available for navigating and operating dialog boxes:

- ENTER
For buttons, triggers a single-tap action on the default button.
- ESC
Cancels the dialog and closes the dialog box.
- SPACEBAR
For buttons, triggers a single-tap action on the focused button.
- TAB
Navigates along the defined navigation path, moving the focus between input boxes and buttons.

Scanner

The system can process scanner input from both 1D and 2D barcode scanners to

- identify materials or equipment during processing by scanning the barcode of a subplot, a container or other piece of equipment,
- start order or workflow processing by scanning the barcode of a batch or treatment ID.

Additionally, you can use a 2D barcode scanner to operate the system itself during processing by scanning the 2D barcodes provided on the buttons in the work area of the Execution Window.

When you operate the system with barcode scans, which is a typical scenario in Dispense operations, the system indicates that it has received an input from a scanner by changing the blue background gradient of the active phase to orange.

Pos.	Material	Allocated or identified batch / subplot / container	Planned (Original)	Recorded	Remaining	# Target	Cleaning	Status	ID only
10	D080-30 / Sonolin		100.0 g (100.0 g)				None	Not started	<input type="checkbox"/>
20	D001-03 / Aqua purificata		20.0 g (20.0 g)				None	Not started	<input type="checkbox"/>
30	D001-03 / Aqua purificata		70.0 g (70.0 g)				None	Not started	<input type="checkbox"/>
40	D005-02 / Magnesium stearate		800.0 g (800.0 g)				None	Not started	<input type="checkbox"/>
50	D003-04 / Ethanol 96 %		10.0 g (10.0 g)				None	Not started	<input checked="" type="checkbox"/>
60	D003-04 / Ethanol 96 %		10.0 g (10.0 g)				None	Not started	<input checked="" type="checkbox"/>
70	D003-04 / Ethanol 96 %		10.0 g (10.0 g)				None	Not started	<input checked="" type="checkbox"/>
80	D005-01 / Lactose		80.0 g (80.0 g)				None	Not started	<input type="checkbox"/>
90	D002-06 / Glucose-Monohydrate		20.0 g (20.0 g)				None	Not started	<input type="checkbox"/>
100	D005-04 / Talc		20.0 g (20.0 g)				None	Not started	<input type="checkbox"/>

Figure 25: Scan in progress

Touch Screen

The system is optimized for running on a touch-screen monitor. The buttons in the work area and the control bar are sized to be large enough for accurate tapping even when you wear protective gear such as gloves.

In addition to the **Page Up** and **Page Down** buttons for paging through the content of the work area, you can also pan the work area manually.

Basic Operations

The following sections describe basic and recurring operations and functions in production execution with PharmaSuite.

Start, Login, Logout, and Password Change

Before you can start working with PharmaSuite your system administrator must have created a user account for you. The PharmaSuite administrator will inform you of your login name and initial password.

Depending on your company policy you may be forced to change your password when you log in for the first time. In this case, the system will display a message that indicates that your password has expired. Then the system will prompt you to change it (page 60).

START PHARMASUITE

To start PharmaSuite double-tap the respective icon on the user interface or select it from the start menu. The system displays the webstart page in a browser window, from which you can select to start an application or view either the help system or the documentation.

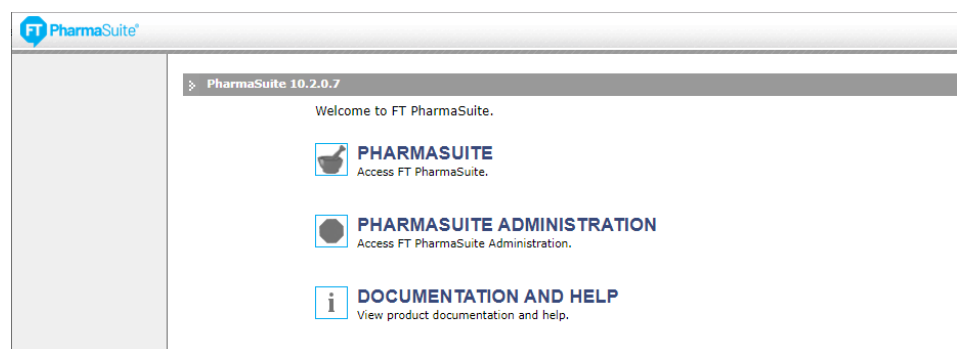


Figure 26: PharmaSuite webstart page

LOGIN

When you select to start PharmaSuite it runs through an initialization phase in the course of which you will also see the splash screen of Shop Operations, which is the internal platform of PharmaSuite. As soon as the initialization phase has been completed, the login form for user login appears.

The login form contains two mandatory fields, one for the login name and one for the password. Your login name and your password are unique for all PharmaSuite applications and are linked to your role and user privileges.

Type your login name and password in the respective boxes. Please note that your password is masked by asterisks (*). Tap the **OK** button to complete the login procedure. If your login attempt is not successful, a message appears, and you have to repeat the procedure.



Figure 27: Login form

After you have successfully logged in, the system displays the PharmaSuite welcome page. From here you can start the Production Execution, Production Execution Viewer, Production Management, Data Manager, Recipe and Workflow Designer, and Production Responses applications, change your password and work station, or access the system documentation and help.



Figure 28: PharmaSuite welcome page

TIP

Please note that logins can be linked to access rights, which means that you can only start an application if your system administrator has assigned the suitable access privileges to you.

Some logins, especially in the production execution environment, can be connected directly to a work station, an application and, within the application, to one or more processes or workflows. If configured in this way, the system will skip the welcome page and start the application automatically with the specified process or workflow after you have successfully logged in.

LOGOUT

In PharmaSuite for Production Execution, open the Cockpit, navigate to the **System** section, and select the **Exit** function to quit PharmaSuite and return to the webstart page.

On the PharmaSuite welcome page, the **Logout** button is located in the top right corner. Tap it to return to the webstart page.

If you decide to log out from a running application, the system will request you to confirm the decision.

PASSWORD CHANGE

You can access the function for changing your password from the PharmaSuite welcome page.

1. Tap the **Change Password** link to open the **Change Password** form.
When your password expires the system will open the form automatically. This can also happen when you log in for the first time to force you to change the initial password, which your system administrator defined for you.
2. On the **Change Password** form, the **User Name** box is output-only and contains your login name.
3. Type your current password in the **Old Password** box.

4. Type your new password first in the **New Password** box and then in the **Confirm New Password** box.
For security reasons, passwords are masked by asterisks (*).
5. Tap the **OK** button to close the form.
From now on, use the new password to log in.

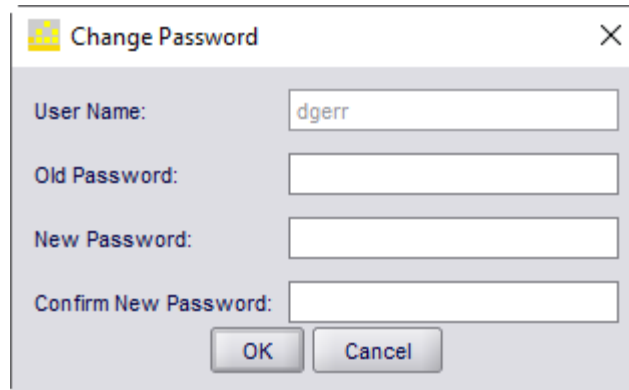


Figure 29: Change password

REQUIRED SERVERS

For providing its full functional scope, PharmaSuite relies on the following servers that are responsible for communication to external systems or between its applications.

- **Electronic Batch Recording (EBR) server**
It controls the execution of EBR recipes and workflows and can process incoming messages from a Distributed Control System.
- **Triggered Operation Management (TOM) server**
It manages event-triggered operations.
- **Operation Execution (OE) server**
It controls the execution of server-run operations.
- **Automation Integration (AI) server**
It controls the communication with automation-related systems.
- **Transition server**
It performs automatic, system-triggered status changes on objects, such as master recipes, master workflows, batches, orders, workflows, or equipment entities and can process incoming messages from external systems, such as a Quality Management System or Warehouse Management.

PharmaSuite runs a heartbeat check on the servers to monitor their availability. To see if there are any issues, switch to the Cockpit and open the **About PharmaSuite** dialog (page 29), which shows the status of the EBR server. For information on the other servers, open the **Details** dialog and refer to the section that indicates servers with heartbeat issues.

If one of the servers does not respond to the heartbeat check while execution is active, the system displays a message to indicate this potential issue.

Executing a Unit Procedure

To execute the unit procedure of recipes or EBR workflows, proceed as follows:

1. In the Cockpit (page 29), select the order step of the recipe or the EBR workflow you wish to execute from the list of **Startable** processes and workflows.
From the list of **Running** unit procedures, select the operation you wish to process.
The system starts the first phase of the operation and switches to the Execution Window.
2. In the Execution Window (page 42), complete and confirm the phases as indicated and available on the screen.
3. To monitor your progress and view the history or executed steps, open the Navigator (page 38).
4. Use the Exception Window (page 44) to record and view exceptions that may occur or have occurred during processing.
5. Continue to execute the phases until you have confirmed the last phase available.
The system returns to the Cockpit from where you can
 - start another operation in the current unit procedure, if available.
 - start another unit procedure.
 - resume an operation that was terminated in a previous session.
 - switch to another unit procedure that is running in the background.
6. When you have finished your work, select the **Exit** (page 60) function from the **System** section of the Cockpit.

TIP

Please note that you can only exit the system and thus log off from your current station when you have either completed your running processes or detached all running operations (page 93).

Handling Exceptions

Exceptions need to be recorded to document irregular circumstances that have occurred during processing. They are an important element of the processing report (batch, workflow). An exception can, on the one hand, be triggered by the phase that is being processed, for example when an operator enters values that exceed the limits defined for the phase. On the other hand, it can be the operator who decides that an exception needs to be recorded, either during processing, through the Exception (page 6) button on the control bar, or after a phase has been completed and confirmed, from the information column of the Navigator (page 38).

Once you have recorded an in-process exception for a phase and return to the Execution Window (page 42), the phase is displayed with a red upper left corner as exception marker. The same red exception marker also shows in the Navigator (page 38) for each phase for which an exception has been registered.

Exceptions can no longer be edited once they have been recorded. If you need to add further information or wish to indicate that the risk or status of the exception should be changed, for example when you have performed a review, you can add comments to any exception.

TIP

Please note that some exceptions include a mandatory comment, which you need to supply before you add the exception.

Exception identifiers are numbered consecutively per batch/order. Comment identifiers consist of the identifier of the exception to which they refer and a two-digit consecutive number.

Additionally, you can bookmark exceptions you wish to find more easily for a later review, for example if you have started the review and could not complete it, or for any other purpose.

System-triggered Exception

When a phase has restrictions or conditions defined for the values it records, deviations from the permitted values cause the system to respond with an information, error, question, or exception dialog. The dialogs indicate their level of severity by a colored frame (blue for information or question, yellow for exception, red for error) and a corresponding markup icon in the upper left corner. The severity itself depends on the inherent capabilities and requirements of each phase as well as the values and limits defined during phase configuration. To indicate the exact context of a message, the dialog displays context information in its upper right corner, which consists of the timestamp of the message's occurrence and the identifiers of the operation and phase to which it refers.

Error dialogs can only be confirmed with the **OK** button but can also display a **Details** button to access further information on the error, such as internal error messages sent by the database or from an external system. If you wish to retain the information listed in the **Details** dialog, click the **Copy to clipboard** button to copy the information and then paste it into a system-external text editor.

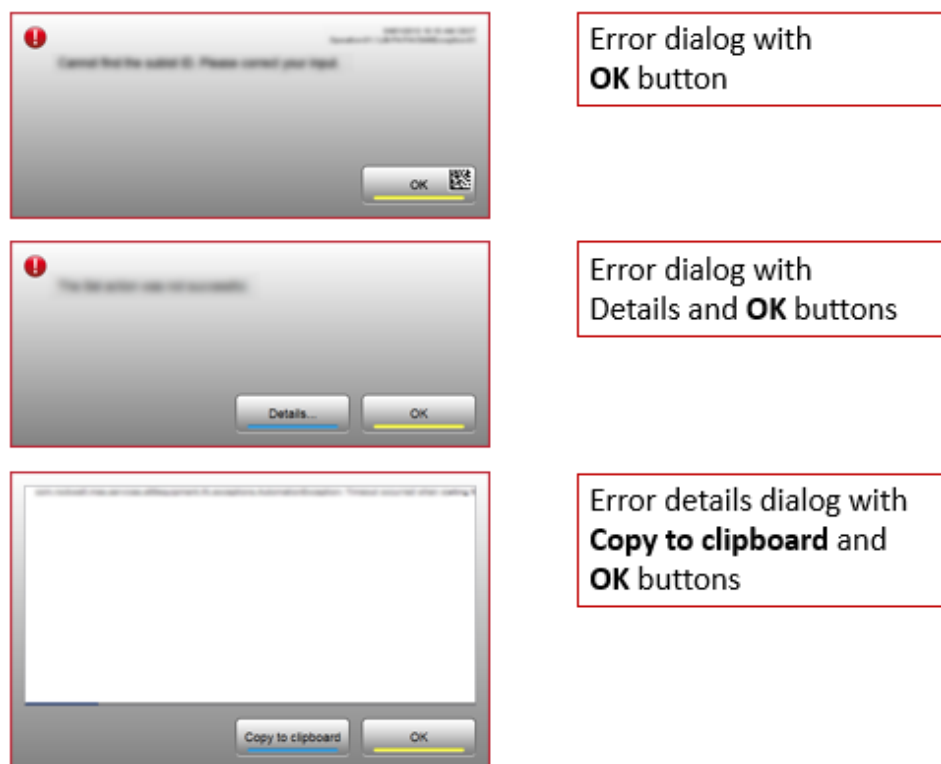


Figure 30: Error and error details dialogs

Information dialogs can only be confirmed with the **OK** button, questions can be answered with **Yes** or **No** and may initiate the recording of an exception.

In exception dialogs you always trigger the recording of an exception. If you accidentally triggered the exception, for example due to an input error, you can cancel this specific exception by tapping the **Back** button in the Exception Window. In this case, however, the system still requires you to sign an **Exception canceled** exception and thus ensures the integrity of all data entered during execution.

Some phases collect several check results to feed them into one dialog. For these cases, the dialogs provide a scrollbar and you can use panning to scroll through the dialog text.

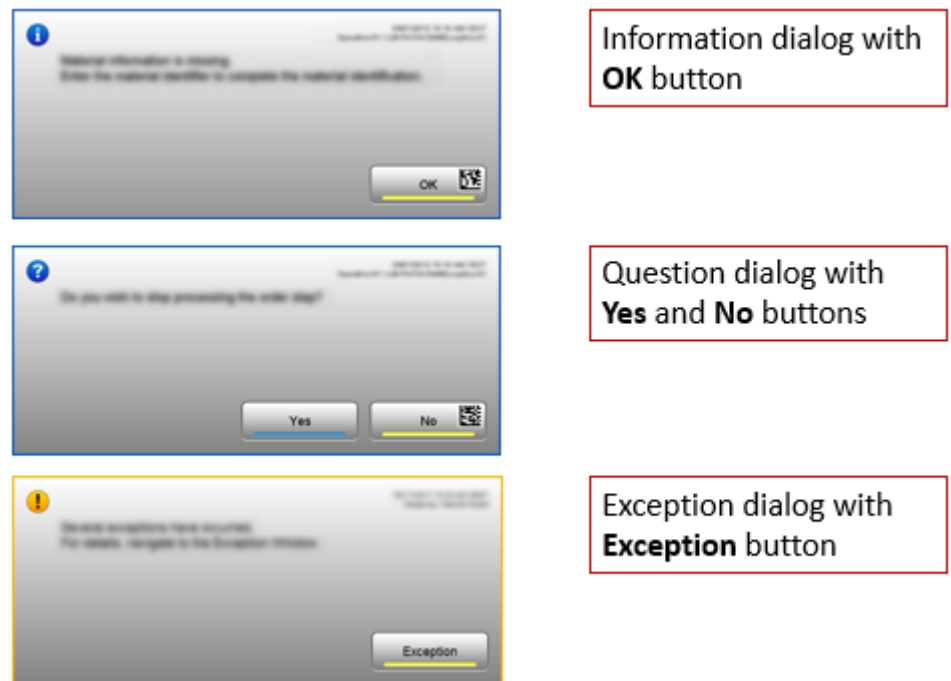


Figure 31: Information, question, and exception dialogs

A system-triggered exception consists of

- a read-only text box that is pre-filled by the system with a description why the current situation represents an exception,
- only if the exception includes a mandatory comment, a text box for entering the comment,
- only if a list of pre-defined exception texts is configured, a **Pre-defined** option list for selecting a text and a read-only **Text** box to display the selected text,
- a read-only option list that is pre-filled by the system with the risk that is attributed to the exception (**None**, **Low**, **Medium**, **High**),

- a read-only option list for the status of the exception (**Open**, **In review**, **To be closed**, **Closed**) that is initially preset with **Open**, and
- a signature panel with **Performed by** and **Password** text boxes to document who recorded the exception.

TIPS

Which of the options for **Risk** and **Status** are available depends on user rights and system configuration.

A **Classification** is only available during comment recording (page 73).

User-triggered Exception

Whenever you encounter an irregular situation during processing you can record it as exception. The system usually provides a free-text description box where you can specify the circumstances of the exception. In addition to free-text exceptions, some situations provide specific exception options to cover recurring irregular events. A typical example of a situation-specific exception could be that a faulty scanner forces you to perform a manual identification instead of the expected barcode scanning.

User-triggered exceptions are recorded from a variety of situations and contexts:

- When you record an exception during phase processing, it is linked to the active phase from which you have triggered the Exception Window. You can access it later for adding comments or bookmarks either from the phase-specific Exception Window (page 44) or after phase completion from the Detail Information Window (page 50).
- When you record an exception for an operation, for instance to document that the operation was detached from its processing device, it affects all active phases of the operation and is thus inserted for each phase. After processing has resumed on another device, you can access the exception later for adding comments or bookmarks either on the phase-specific Exception Window (page 44) of each phase or, after phase completion by way of the Navigator (page 38), on their Detail Information Windows (page 50).
- When you record an exception for a unit procedure, for instance to document that the unit procedure was paused, it is linked to the unit procedure. You can access it later by way of the unit procedure's Actions dialog (page 78), for adding comments or bookmarks on the overview Exception Window (page 47).

A user-triggered free-text exception consists of

- a text box where you can describe the exception circumstances,
- only if a list of pre-defined exception texts is configured, a **Pre-defined** option list for selecting a text and a read-only **Text** box to display the selected text,
- an option list from which you can select the risk to the process you attribute to the exception (**High, Medium, Low, None**),
- a read-only option list for the status of the exception (**Open, In review, To be closed, Closed**) that is initially preset with **Open**, and
- a signature panel with **Performed by** and **Password** text boxes to document who recorded the exception.

TIPS

Which of the options for **Risk** and **Status** are available depends on user rights and system configuration.

A **Classification** is only available during comment recording (page 73).

A user-triggered situation-specific exception consists of

- a read-only text box that is pre-filled by the system with a description why the current situation represents an exception,
- only if the exception includes a mandatory comment, a text box for entering the comment,
- only if a list of pre-defined exception texts is configured, a **Pre-defined** option list for selecting a text and a read-only **Text** box to display the selected text,
- a read-only option list that is pre-filled by the system with the risk that is attributed to the exception (**None, Low, Medium, High**),
- a read-only option list for the status of the exception (**Open, In review, To be closed, Closed**) that is initially preset with **Open**, and
- a signature panel with **Performed by** and **Password** text boxes to document who recorded the exception

TIPS

Which of the options for **Risk** and **Status** are available depends on user rights and system configuration.

A **Classification** is only available during comment recording (page 73).

RECORDING AN EXCEPTION

Exception recording differs depending on the circumstances that cause the exception, whether it is operator- or system-triggered, free-text or phase-specific. Some exceptions only support free-text entries for exception texts or comments, while others have a pre-defined option list of exception or comment texts configured.

Figure 32: Exception recording panel with free-text entry

Figure 33: Exception recording panel with pre-defined text options

To record a user-triggered free-text exception, proceed as follows:

1. Tap the **Exception** button in the control bar or the respective button in the information column of the Navigator to access the exception recording panel. At the top, the form displays a list of all exceptions and comments that have so far been recorded for the phase. If you are about to enter the first exception, the list will be empty.
At the bottom it shows the exception recording panel where you can enter the exception.
2. In the **Exception** box, type a description of the circumstances that cause you to record the exception.

TIP

For very long exception descriptions, only the first 1984 characters are displayed as exception text. Any description text characters beyond the limit are displayed as one or more comments to the exception.

3. If a list of pre-defined exception texts is configured, the system displays the **Pre-defined** option list, which shows the headers of the available text options. Select an option.
The system displays the full text of the selected option in the read-only **Text** box. Later on, the pre-defined comment text is added as exception comment to the related exception.
4. Select the suitable risk level you attribute to the exception.

TIP

Please note that the **Status** and **Classification** boxes are read-only. They are only available for editing when adding a comment (page 73) to the exception.

5. Type your login name and password in the respective boxes. Please note that your password is masked by asterisks (*).
For situations requiring a witness, the system will ask not only for a single but for a double signature. In these cases, two different users, typically with different qualifications, have to complete the signature.

TIP

Please note that for security reasons the system deletes the typed password if you tap any active box or button other than the **Add exception** button so you need to retype it.

6. Tap the **Add** button to record the exception.
The system will now display it in the list of recorded exceptions.
Once you have recorded the exception, any situation-specific exception options the phase may provide are disabled (grayed out). You can, however, add a comment to the exception or bookmark it.

TIP

Please note that you can only record one exception at a time. To record another exception, you have to start the exception recording process from scratch by accessing the Exception Window again.

7. Tap the **Back** button to return to the form from which you have accessed the Exception Window.

You can record user-triggered, situation-specific exceptions for phases, operations, or unit procedures. Proceed as follows:

1. Access the Exception Window to record the situation-specific exception.

■ For a phase:

Tap the **Exception** button in the control bar or the respective button in the information column of the Navigator to access the exception recording panel. At the top, you will see one or more pre-defined exception options that cover expected exception cases of the specific phase. Below the pre-defined exception options, the form displays a list of all exceptions and comments that have so far been recorded for the phase. If you are about to enter the first exception, the list will be empty.

Select the suitable exception option, fill in any information that may be required, and tap the **Confirm** button.

The system displays a pre-defined description of the exception. The system displays a pre-defined description of the exception situation and its risk level in read-only mode.

TIP

In case none of the pre-defined exception options fit your exception situation, the system also provides the exception recording panel where you can enter a free-text exception (page 68).

■ For an operation:

Performing an action on an operation that in the current context represents an exception situation, such as detaching the operation, automatically opens the Exception Window, showing only the exception recording panel. The system displays a pre-defined description of the exception situation and its risk level in read-only mode.

■ For a unit procedure:

Performing an action on a unit procedure that in the current context represents an exception situation, such as aborting and reactivating the unit procedure, automatically opens the Exception Window, showing only the exception recording panel. The system displays a pre-defined description of the exception situation and its risk level in read-only mode.

2. If a list of pre-defined exception texts is configured, the system displays the **Pre-defined** option list, which shows the headers of the available text options. Select an option.
The system displays the full text of the selected option in the read-only **Text** box. Later on, the pre-defined comment text is added as exception comment to the related exception.

TIPS

Please note that some exceptions include a mandatory comment, which you need to supply before you add the exception.

The **Status** and **Classification** boxes are read-only. They are only available for editing when adding a comment (page 73) to the exception.

3. Type your login name and password in the respective boxes. Please note that your password is masked by asterisks (*).
For situations requiring a witness, the system will ask not only for a single but for a double signature. In these cases, two different users, typically with different qualifications, have to complete the signature.

TIP

Please note that for security reasons the system deletes the typed password if you tap any active box or button other than the **Add exception** button so you need to retype it.

4. Tap the **Add** button to record the exception.
The system will now display it in the list of recorded exceptions.
Once you have recorded the exception, the exception options are disabled (grayed out). You can, however, add a comment to the exception or bookmark it.

TIP

Please note that you can only record one exception at a time. To record another exception, you have to start the exception recording process from scratch by accessing the Exception Window again.

5. Tap the **Back** button to return to the form from which you have accessed the Exception Window.

To record a system-triggered exception, proceed as follows:

1. Tap the **Exception** button on a warning or exception dialog to open the Exception Window. At the top, the form displays a list of all exceptions and comments that have so far been recorded for the phase. If you are about to enter the first exception, the list will be empty.

At the bottom it shows the exception recording panel where you can enter the exception.

The system displays a pre-defined description of the exception situation and its risk level in read-only mode.

TIP

Please note that some exceptions include a mandatory comment, which you need to supply before you add the exception.

The Exception Window only displays the system-triggered exception you are about to record. If the phase provides phase-specific, user-triggered exceptions (page 70), you cannot access them at this point. However, the list of recorded exceptions shows all exceptions that have been registered so far, regardless of their origin.

2. If a list of pre-defined exception texts is configured, the system displays the **Pre-defined** option list, which shows the headers of the available text options. Select an option.
The system displays the full text of the selected option in the read-only **Text** box. Later on, the pre-defined comment text is added as exception comment to the related exception.

TIP

Please note that the **Status** and **Classification** boxes are read-only. They are only available for editing when adding a comment (page 73) to the exception.

3. Type your login name and password in the respective boxes. Please note that your password is masked by asterisks (*).
For situations requiring a witness, the system will ask not only for a single but for a double signature. In these cases, two different users, typically with different qualifications, have to complete the signature.

TIP

Please note that for security reasons the system deletes the typed password if you tap any active box or button other than the **Add exception** button so you need to retype it.

4. Tap the **Add** button to record the exception.
The system will now display it in the list of recorded exceptions.
5. Tap the **Back** button to return to the form from which you have accessed the Exception Window.

RECORDING A COMMENT

TIP

Please note that an exception may be under review in another session of PharmaSuite for Production Execution or in the Exception Dashboard of PharmaSuite for Production Responses.

To record a comment, proceed as follows:

1. Tap the **Exception** button in the control bar or the respective button in the information column of the Navigator to access the exception recording panel. At the bottom, it shows the panel where you can enter the comment.
2. In the list of recorded exceptions, tap the exception to which you wish to add a comment.
3. In the **Comment** box, type a description of the circumstances that cause you to add the comment.
4. If a list of pre-defined exception texts is configured, the system displays the **Pre-defined** option list, which shows the headers of the available text options. Select an option.
The system displays the full text of the selected option in the read-only **Text** box. Later on, the pre-defined comment text is added as exception comment to the related exception.
5. To change the risk level (**None**, **Low**, **Medium**, **High**) or the status (**Open**, **In review**, **To be closed**) of the exception to which the comment refers, select the new risk level or status. The new value will be propagated to the exception once you have added the comment.

TIP

Please note that setting the risk level of an exception to **None** automatically sets its status to **N/A** only. The **N/A** status will not be available when the risk level is other than **None**.

6. If required, select a **Classification**, such as **Information only** or **Machine error** for the exception. The selected classification will be propagated to the exception once you have added the comment.

TIP

An exception can only be closed from the Exception Dashboard of PharmaSuite for Production Responses after it has been reviewed.

7. Type your login name and password in the respective boxes. Please note that your password is masked by asterisks (*).
8. Tap the **Add** button to record the comment.
The system will now display it in the list of recorded exceptions.
9. Tap the **Back** button to return to the form from which you have accessed the Exception Window.

BOOKMARKING AN EXCEPTION

TIP

Please note that an exception may be under review in another session of PharmaSuite for Production Execution or in the Exception Dashboard of PharmaSuite for Production Responses.

To bookmark an exception, proceed as follows:

1. Tap the **Exception** button in the control bar or the respective button in the information column of the Navigator to access the exception recording panel. At the bottom it shows the panel where you can bookmark the exception.
2. In the list of recorded exceptions, tap the exception you wish to bookmark.
3. Tap the **Bookmark** button to bookmark the exception.
The system will now display the marker it in the list of recorded exceptions.

TIP

Tapping the **Bookmark** button for a bookmarked exception will remove the marker.

4. Tap the **Back** button to return to the form from which you have accessed the Exception Window.

Aborting and Repairing Phases

To resolve processing issues that can result from a faulty or inconsistent process parameter definition during recipe or workflow creation, the system provides functions to abort or abort and repair a phase.

When a user aborts a phase the system terminates its execution and moves on to the next phase in the process. An aborted phase does not record data or provide phase outputs, which means that all subsequent phases that reference an output of the aborted phase will be impacted.

When a user aborts and repairs a phase, the system enters a specific repair mode in which it displays the phase's process parameters and allows the user to change them with suitable editors. Expressions, however, cannot be changed and there is no Expression editor available. This means that it is not possible to change an equipment object reference specified with an equipment requirement parameter. Apart from this restriction, however, the system provides the same editors that are used for defining recipes or workflows in the Process Parameter Panel of Recipe and Workflow Designer. For further information, please refer to "Editors" in Vol. 1 "Introduction and Basics" of "User Manual Recipe and Workflow Designer".

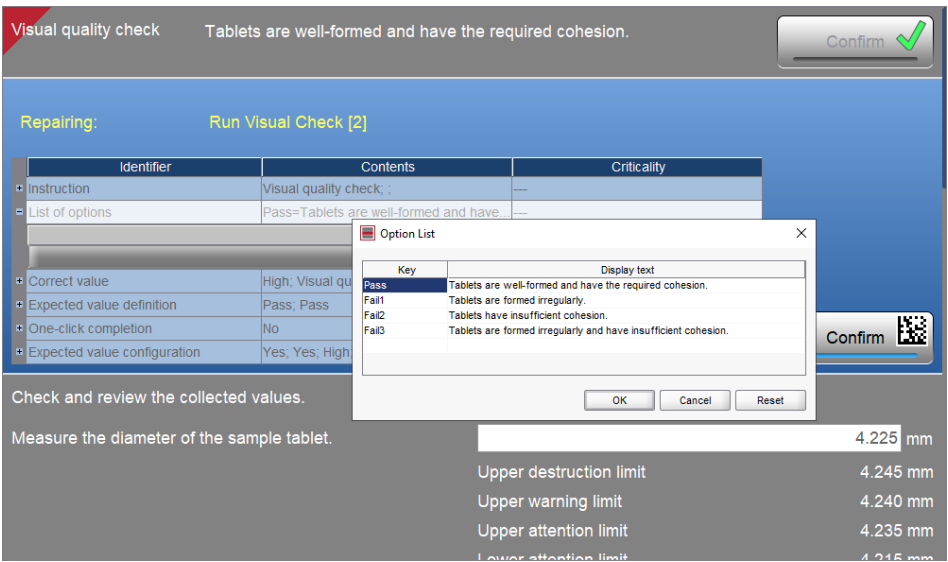


Figure 34: Phase repair mode with Option List editor

TIPS

Please note that the repair mode may not be available for all phases, the system can be configured to exclude specific phases from being repairable.

The repair mode is only available for phases for which no exceptions other than **Exception canceled** exceptions (page 64) and exceptions related to aborting and repairing have been recorded yet.

To abort a phase, proceed as follows:

1. In the right column of the Navigator (page 38), tap the **Abort** button.
The system displays a confirmation dialog.
2. Tap **Yes** to confirm the dialog.
The system displays the Exception Window (page 44).
3. Sign and add the **Phase execution aborted** exception.
4. In the Control Bar (page 6), tap the **Back** button to return to the Execution Window, where the aborted phase is shown as completed with the exception marker in the top left corner.

To abort and repair a phase, proceed as follows:

1. In the right column of the Navigator (page 38), tap the **Repair** button.
The system displays a confirmation dialog.
2. Tap **Yes** to confirm the dialog.
The system displays the Exception Window (page 44).
3. Sign and add the **Phase repair mode started** exception.
4. In the Control Bar (page 6), tap the **Back** button to return to the Execution Window, where the phase is shown as aborted with the exception marker in the top left corner.
The system displays the phase in the repair mode with yellow phase text. It shows the identifier of the phase in the top row and below that, the list of the phase's process parameters. In the list, each process parameter is displayed as a nested table, with the common parameter attributes as main table and each parameter's specific attributes in a collapsible sub-table.
To edit the values of the parameter attributes, proceed as follows:

1. Tap the table cell you wish to edit. The system displays a text cursor or a button to access a suitable editor for the required value type.
2. Enter the value of the parameter attribute.
3. Tap in another cell or the ENTER key to conclude the edit action.

When you have made all required changes, tap the **Confirm** button.
The system displays the Exception Window (page 44).

5. Sign and add the **Phase parameters repaired** exception.
6. In the Control Bar (page 6), tap the **Back** button to return to the Execution Window where the phase is shown as active with the exception marker in the top left corner. It is ready for processing with the changed values.

Accessing Cockpit Actions

From the Cockpit, you can access additional actions that may be available to items listed as **Running**, such as the actions required for managing event-triggered operations and their runs.

To access the additional actions of a running item, proceed as follows:

1. Open the Cockpit where the template or operation is displayed in the list of **Running** processes.
2. Tap its **Actions** button (page 13).
The system displays the respective **Actions** dialog that holds a tappable button for each action available.
Depending on the context from which you call the **Actions** dialog, some of the actions may be disabled.
3. Tap the action you wish to perform.
The system closes the dialog and performs the action.
To close the dialog and return to the Cockpit without performing an action, tap the **Close** button.

ACTIONS OF UNIT PROCEDURES

The Actions dialog of a unit procedure holds the following actions:

- **Pause**
Only for running unit procedures, pauses the unit procedure (page 84), thus preventing trigger phases from generating triggers.
- **Resume**
Only for paused unit procedures, resumes a unit procedure, which returns to operating and informs all affected phases of this fact.
- **Exceptions**
Opens the overview Exception Window (page 47) to display all exceptions and comments that have so far been recorded for the unit procedure.
- **Abort and reactivate**
Only available for unit procedures of orders.
Terminates processing of the unit procedure and makes it available for restart (page 92).
- **Open in Viewer**
Opens Production Execution Viewer with a selected unit procedure.

TIP

Please note that in order for this option to work, Production Execution Viewer must be running in the background.

■ Append

Opens the Selection Window for **Appendable Workflows** to filter and select a workflow (page 88) that needs to be appended to the unit procedure.

TIPS

Please note that you can append workflows only to unit procedures of batch processing, but not to unit procedures of workflows. On the **Actions** dialog of a workflow's unit procedure, the **Append** button is hence disabled. Workflows that are based on a confidential master workflow to which you do not have the required access rights are not listed and are thus not available for selecting and appending.

■ Cancel

Terminates a running workflow that cannot be completed by regular means.

TIP

Please note that you cannot cancel a workflow that is appended to a unit procedure.

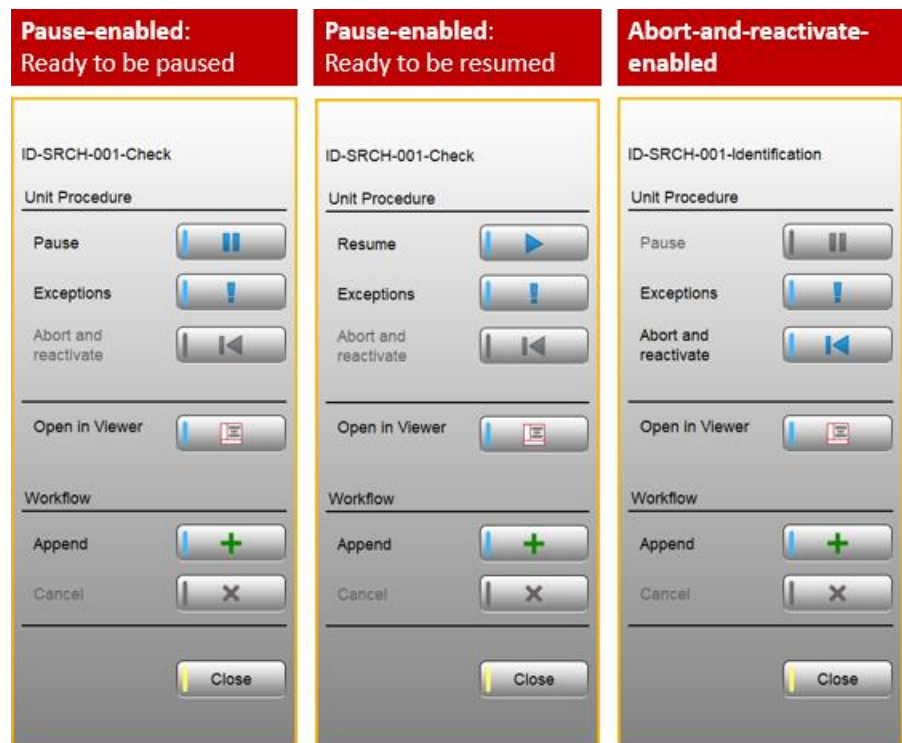


Figure 35: Actions dialogs of order unit procedures



Figure 36: Actions dialogs of workflow unit procedures

ACTIONS OF OPERATIONS

The Actions dialog of an operation holds the following actions:

- **Exceptions**
Opens the operation-specific Exception Window (page 46) to display all exceptions and comments that have so far been recorded for the operation.
- **Abort**
Terminates processing of the operation.
- **Open in Viewer**
Opens Production Execution Viewer with a selected operation

Please note that in order for this action to work, Production Execution Viewer must be running in the background.

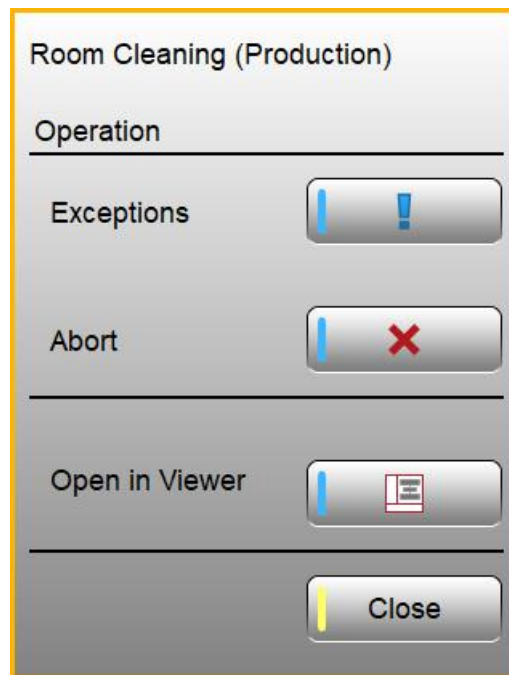


Figure 37: Actions dialog of operations

ACTIONS OF EVENT-TRIGGERED OPERATIONS

During processing, an event-triggered operation only exists as non-executable template until an external trigger or an operator action creates an executable run of the template. Once a run has been created, an operator can start it like any other operation to process its phases. When the operator confirms the last phase of the run, the template remains available for creating further runs. It is also possible to create additional runs of the template while there are already startable runs or runs that are being processed. As long as a template is displayed during processing and thus available for creating new runs, an operator cannot start succeeding operations for processing. For this reason, the system provides a function to remove the template, so that the operator can continue with processing.

Event-triggered operations are especially suited for processes that need to be performed repeatedly at definable intervals, such as in-process controls (IPC). They show in the Cockpit either as templates or as runs.

Their Actions dialog holds the following actions:

- **New run**
Available for both templates and runs.
Starts a new run of the operation.
- **Cancel run**
Only available if the run has not been started yet.
Cancels the run of the operation.
- **Cancel all runs**
Only available if there is more than one run available.
Cancels all runs of the operation that have not been started yet.

- Open in Viewer

Opens Production Execution Viewer with a selected event-triggered operation.

Please note that in order for this action to work, Production Execution Viewer must be running in the background.

- Remove template

Only available for templates.

Removes the template from the list of running processes, thus ending the cycle of the event-triggered operation.

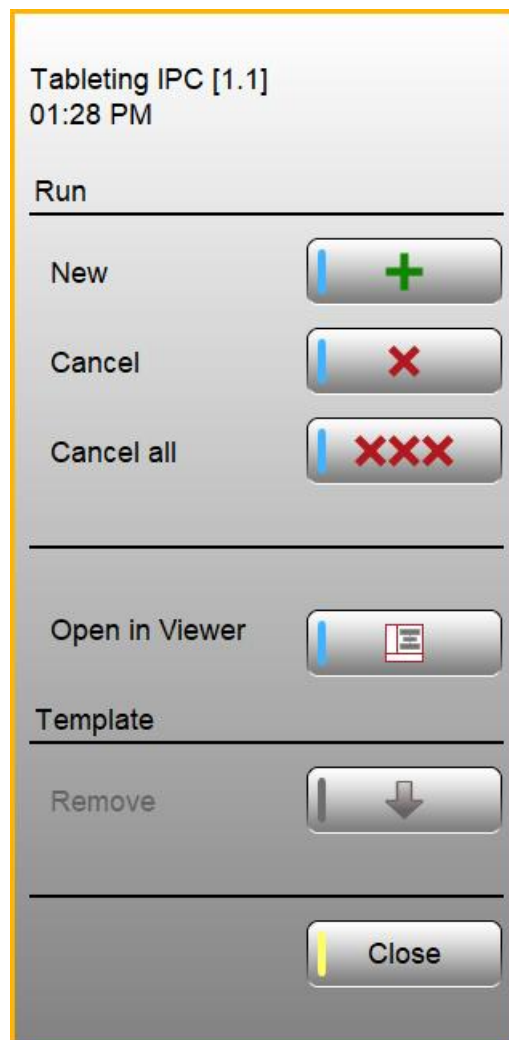


Figure 38: Actions dialog of an event-triggered operation

Pausing Unit Procedures

Some unit procedures allow an operator to pause their execution. When the operator sets a unit procedure to paused, the system informs all currently active phases of this fact. Phases that run without user interaction, such as trigger phases, react to this information by holding their triggering activities. This prevents the system from stacking runs of event-triggered operations while the execution itself is paused.

To pause a unit procedure, proceed as follows:

1. Open the Cockpit where the unit procedure is displayed in the list of **Running** processes.
2. Tap the **Actions** button of the unit procedure.
The system displays the Actions dialog of the unit procedure (page 78).
3. Tap the **Pause unit procedure** button (page 16).

If defined so in the recipe, the system opens the Exception Window to record the pausing as exceptional situation.

If not, the system immediately marks the unit procedure as paused by displaying the **Paused unit procedure** icon (page 21) to the left of its unit procedure button and by showing a corresponding marker tab (page 12) in the notification bar.

- If pausing the unit procedure is not critical to the process and is thus permitted without restrictions, the system immediately marks the unit procedure as paused by displaying the **Paused unit procedure** icon (page 21) to the left of its unit procedure button and by showing a corresponding marker tab (page 12) in the notification bar.
- If pausing the unit procedure represents an exception situation, the Exception Window opens and displays a corresponding exception, which you need to sign with an electronic signature.
The exception description is pre-defined and read-only, but you can add a comment if you need to provide additional information. You can also access it later by way of the unit procedure's Actions dialog (page 78), for adding comments or bookmarks on the overview Exception Window (page 47).
Tap the **Back** button to return to the Cockpit.

To resume a paused unit procedure, proceed as follows:

1. Tap the **Actions** button of the paused unit procedure you wish to resume.
The system displays the Actions dialog of the unit procedure (page 78).
2. Tap the **Resume unit procedure** button.
The system removes the marker icon and the marker tab. The unit procedure resumes its execution, including the generation of triggers.
3. Tap the **Back** button to return to the Cockpit.

TIP

Please note that pausing a unit procedure does not interfere with the escalation scenario defined for existing runs of an event-triggered operation. An existing run still passes through all stages of its escalation scenario, even if its unit procedure is paused.

Appending Workflows to Orders

Some workflows are executed within the context of an order, such as specific cleaning or preparation procedures. For this reason, they need to be included in the order's batch report. By appending a workflow to those unit procedures of an order to which their execution applies, you establish the connection and thus include the processing details of the workflow in the order's batch report.

TIPS

Please note that only those workflows whose master workflow is marked as **Production-relevant** and **Appendable during processing** are available for being appended to unit procedures of batch processing.

The system can be configured to request an electronic signature when an operator appends a workflow to a unit procedure.

Workflows that are based on a confidential master workflow to which you do not have the required access rights are not listed and are thus not available for selecting and appending.

When appending workflows to unit procedures, PharmaSuite for Production Execution supports several scenarios:

- At workflow start, the workflow can be appended to the only unit procedure that is currently running at the work center.
The system displays a corresponding question dialog. After you have selected whether or not you wish to append the workflow to the unit procedure, the appended or unappended workflow is available for processing in the list of **Running** processes of the Cockpit.
- At workflow start, the workflow can be appended to one of several unit procedures that are currently running at the work center.
 1. The system displays a corresponding question dialog.
Tap **Yes**, to open the Selection Window for **Appendable Unit Procedures** that lists all unit procedures that are currently running at your work center.
 2. To reduce the number of listed unit procedures, use the **Filter** box to type the characters you are looking for. The filter becomes effective with the second character you type.
By default, the search applies to all table columns, is not case-sensitive, and the search string may occur anywhere within the content of a cell. Click the magnifying glass icon to set the search criteria to search only in specific columns, be case-sensitive, allow the use of wildcards or regular expressions, or match from start or exactly.

3. Tap to select the unit procedure to which you wish to append the workflow.
4. Tap the **Append** button to append the workflow to the selected unit procedure.

The system returns to the Cockpit where the appended workflow is available for processing.

Tap the **Back** button of the Control Bar to return to the Cockpit and process the workflow without having appended it to the unit procedure.

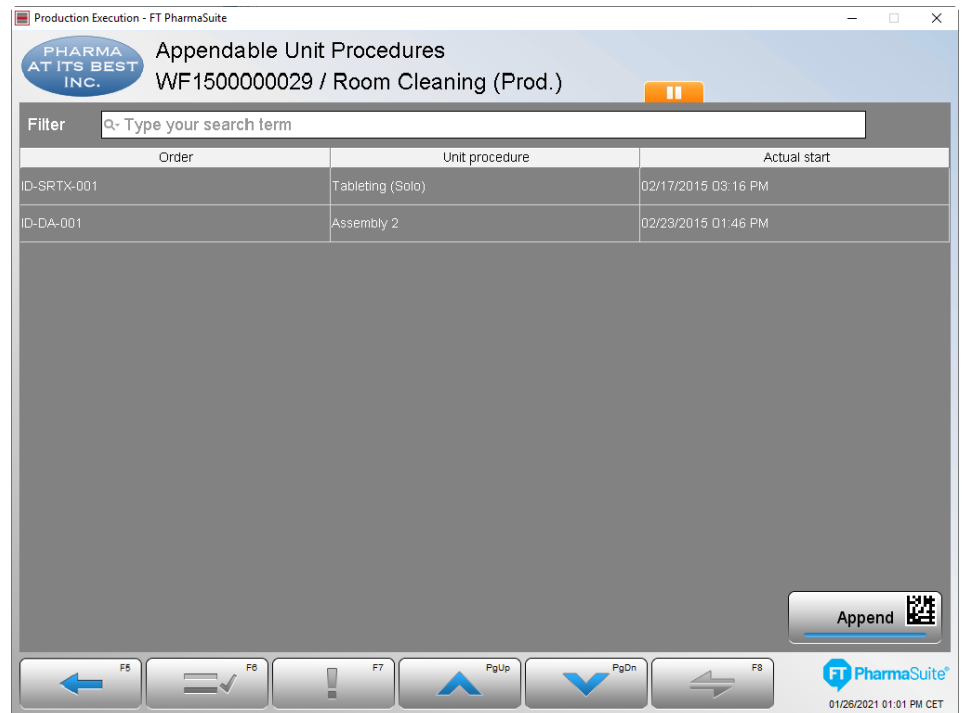


Figure 39: Appendable Unit Procedures

- During unit procedure processing, a workflow for appending can be selected from the list of workflows that are running or have been completed at any work center.

1. Tap the **Actions** button of the unit procedure to which you wish to append a workflow.

The system displays the Actions dialog of the unit procedure (page 78).

TIP

Please note that you can append workflows only to unit procedures of batch processing, but not to unit procedures of workflows. On the **Actions** dialog of a workflow's unit procedure, the **Append** button is hence disabled.

2. Tap the **Append** button (page 13).

The system opens the Selection Window for **Appendable Workflows** that displays a filtered list of workflows available for appending.

In the leftmost column of the list, a checkbox indicates for each workflow if it is already appended to the unit procedure.

TIP

Please note that with this action you can only append workflows, but not unappend workflows that have been appended before. Modifications of this kind can only be done in PharmaSuite for Production Management.

3. You can use the filter tools to restrict the number of listed workflows:

- In the **Filter** box, type the characters you are looking for. The filter becomes effective with the second character you type.
By default, the search applies to all table columns, is not case-sensitive, and the search string may occur anywhere within the content of a cell. Click the magnifying glass icon to set the search criteria to search only in specific columns, be case-sensitive, allow the use of wildcards or regular expressions, or match from start or exactly.
- From the **Started** option list, select how long ago the workflow may have been started.
The default option is **Last 8 hours**.
- From the **Status** option list, select the status the workflow you wish to append needs to have (**Canceled**, **Finished**, **In process**, **Production-reviewed**).
The default option is **Finished**.
- Use the **Running at current work center** button (page 12) to toggle between listing only the workflows that were or are running at your work center or at all work centers.
By default, the button is toggled on and thus the list contains only the workflows of your current work center.

4. Tap to select the workflow you wish to append to the unit procedure.
5. Tap the **Append** button to append the selected workflow to the unit procedure.

The system returns to the Cockpit.

Tap the **Back** button of the Control Bar to return to the Cockpit without having appended a workflow to the unit procedure.

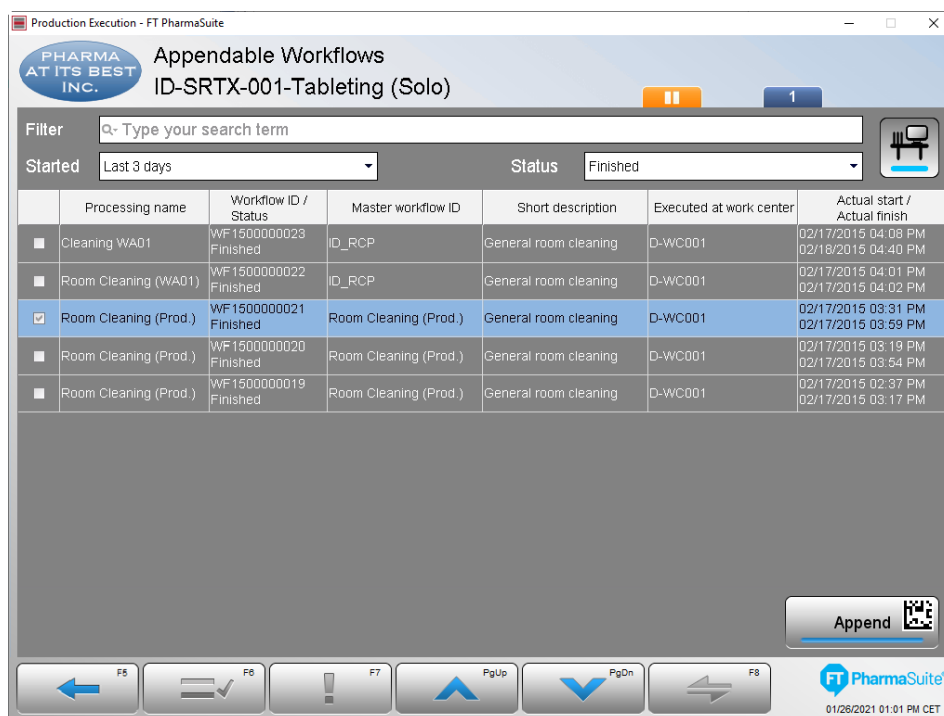


Figure 40: Appendable Workflows

Handling Escalations

Event-triggered operations (page 82) can enforce their timely execution by implementing a four-level escalation scenario, with the following stages:

1. Run due
When a run becomes due, the system displays its tappable operation button (page 15), which also indicates the next stage of the escalation progress and when it will occur. Additionally, the center-left (blue) tab of the notification bar (page 3) reflects this stage.
2. Run due reminder
After a run has been due for a certain amount of time, the system alerts its operator to this fact by displaying the **Run due reminder** icon (page 16) to the left of its tappable operation button (page 14), which also indicates the next stage of the escalation progress and when it will occur. Additionally, the center-right (orange) tab of the notification bar (page 3) reflects this stage.

3. Run overdue

When a run has been due for too long and needs to be processed immediately, the system displays the **Run overdue** icon (page 16) to the left of its tappable operation button (page 14), which also indicates the next stage of the escalation progress and when it will occur. Additionally, the center-right (red) tab of the notification bar (page 3) reflects this stage.

4. Run expired

After a run has been in the overdue stage for some time, it eventually expires and disappears from the list of **Startable** operations.

5. Finish due reminder

When a run is being executed and its completion has been due for a defined amount of time, the system alerts its operator to this fact by displaying the **Finish due reminder** icon (page 16) to the left of its tappable operation button (page 15). Additionally, the number in the left (orange) section of the right-most tab of the notification bar (page 3) reflects this stage.

6. Finish overdue

When a run has not been completed within a defined amount of time and declared overdue, the system alerts its operator to this fact by displaying the **Finish overdue** icon (page 17) to the left of its tappable operation button (page 15). Additionally, the number in the right (red) section of the right-most tab of the notification bar (page 3) reflects this stage.

Detaching Unit Procedures

The system allows an operator to detach a running unit procedure from its work center if

- the unit procedure is not paused and
- none of its operations are running.

Detaching a unit procedure puts it on hold until an operator selects it again for processing at another work center to which it is consequently attached. The operator resumes processing where it was left off, by starting the unit procedure, which continues at the next available operation. As putting a running unit procedure on hold can represent a critical situation that may compromise the quality of the materials that are being processed, a detach action may need to be recorded as exception or may even not be available at all.

To detach a unit procedure, proceed as follows:

1. Open the Cockpit where the unit procedure is displayed in the list of **Running** processes.
2. Tap its **Detach** button (page 14).
 - If detaching the unit procedure is not critical to the process and is thus permitted without restrictions, the button disappears, thus indicating that the system has detached the unit procedure from both the device and the station where the device is registered.
 - If detaching the unit procedure represents an exception situation, the Exception Window opens and displays a corresponding exception, which you need to sign with an electronic signature.
The exception description is pre-defined and read-only, but you can add a comment if you need to provide additional information. You can also access it later by way of the unit procedure's Actions dialog (page 78), for adding comments or bookmarks on the overview Exception Window (page 47). Tap the **Back** button to return to the Cockpit.
3. Now the unit procedure is available again for selection and processing from another work center.

When a unit procedure with detach exception is resumed again, the system automatically adds a comment to that effect to the exception.

Aborting and Reactivating Unit Procedures

The system allows an operator to abort and reactivate a running unit procedure. A new run of the unit procedure is then available for restart from the Selection Window for processing all of its operations and phases.

This action is intended for rework situations when an entire unit procedure needs to be repeated. Any equipment identified and bound to the unit procedure during its aborted run is unbound and available for identification again. As aborting and reactivating a running unit procedure can represent a critical situation that may compromise the quality of the materials that are being processed, the action may need to be recorded as exception or may even not be available at all.

To abort and reactivate a unit procedure, proceed as follows:

1. Open the Cockpit where the unit procedure is displayed in the list of **Running** processes.
2. Tap the **Actions** button of the unit procedure.
The system displays the Actions dialog of the unit procedure (page 78).
3. Tap the **Abort and reactivate** button (page 13).
 - If aborting and reactivating the unit procedure is not critical to the process and is thus permitted without restrictions, the button disappears, thus indicating that the system has detached the unit procedure from the work center.
 - If aborting and reactivating the unit procedure represents an exception situation, the Exception Window opens and displays a corresponding exception, which you need to sign with an electronic signature.
The exception description is pre-defined and read-only, but you can add a comment if you need to provide additional information.
Tap the **Back** button to return to the Cockpit.

Now a new run of the unit procedure is available again for selection and processing.

When a reactivated unit procedure is selected for processing, the system displays its instance count along with its identifier.

Detaching Operations

The system allows an operator to detach a running operation from its device and station, thus putting it on hold until it is selected for processing from another device and thereby attached to it and resumed. As putting a running operation on hold can represent a critical situation that may compromise the quality of the materials that are being processed, a detach action can for some phases lead to a user-triggered exception (page 66) or may even not be available from PharmaSuite for Production Execution at all.

To detach an operation, proceed as follows:

1. Open the Cockpit where the operation is displayed in the list of **Running** processes.
2. Tap its **Detach** button (page 14).
 - If detaching the operation is not critical to the process and is thus permitted without restrictions, the button disappears, thus indicating that the system has detached the operation from both the device and the station where the device is registered.
 - If detaching the operation represents an exception situation, the Exception Window opens and displays a corresponding exception, which you need to sign with an electronic signature.
The exception description is pre-defined and read-only, but you can add a comment if you need to provide additional information. The exception is added to each phase of the operation that is active when the operation is detached. You can access it later for adding comments or bookmarks either on the phase-specific Exception Window (page 44) of each phase or, after phase completion by way of the Navigator (page 38), on their Detail Information Windows (page 50).
Tap the **Back** button to return to the Cockpit.
3. Now the operation is available again for selection and processing from another device or station.

TIP

Please note that the system provides an emergency detach function for non-detachable operations from the **Unlock Objects** task in PharmaSuite for Production Management.

When an operation with detach exception is resumed again, the system automatically adds a comment to that effect to the respective exceptions of each affected phase.

Aborting Operations

The system allows an operator to abort any operation, regardless of whether it has not been started yet or if it is already running.

This action is primarily intended for testing purposes, to provide increased flexibility during the processing of recipes that are still under development.

When an operation is aborted the system records an exception at its unit procedure.

To abort an operation, proceed as follows:

1. Open the Cockpit where the operation is displayed in the list of **Running** processes.
2. Tap the **Actions** button of the operation.
The system displays the Actions dialog of the operation (page 81).
3. Tap the **Abort** button (page 13).
The system opens the Exception Window and displays a corresponding exception, which you need to sign with an electronic signature.
The exception description is pre-defined and read-only, but you can add a comment if you need to provide additional information.
Tap the **Back** button to return to the Cockpit.

TIP

Please note that aborting an operation does not undo the identification of equipment or materials that may have been performed already.

Signature Requests

When performing safety-sensitive or GxP-relevant functions the system may request you to enter an electronic signature, for example during material identification. Signatures are linked to user groups and access privileges, which means that the system will only accept the signature of a user who is qualified to perform the task in question. Unless the required signature data has been entered correctly, subsequent functions cannot be executed.

For situations requiring a witness, the system will ask not only for a single but for a double signature. In these cases, two different users, typically with different qualifications, have to complete the signature form before task processing can continue.

TIP

Please note that the system registers each signature with the timestamp when it has been verified successfully. Thus, the signature timestamps of a double signature will usually be different.

To perform an electronic signature, type your login name and password and click the **OK** button. Free-text comments can be optional or mandatory and may consist of up to 255 characters.

If a list of pre-defined comment texts is configured for a signature, the system displays an option list, which shows the headers of the available text options. Once an option has been selected, the system shows its full text in the read-only text box below the option list. If a pre-defined comment text is configured, it is mandatory. Later on, the pre-defined comment text is added as exception comment to the related exception.

Electronic Signature

Description: Reprint label

Reason: Confirmation

User Name: oraita


Password: [Yellowed out]

Comment: [Yellowed out]

OK Cancel

Identifier:WIP_ES_LABEL_REPRINT

Figure 41: Single electronic signature

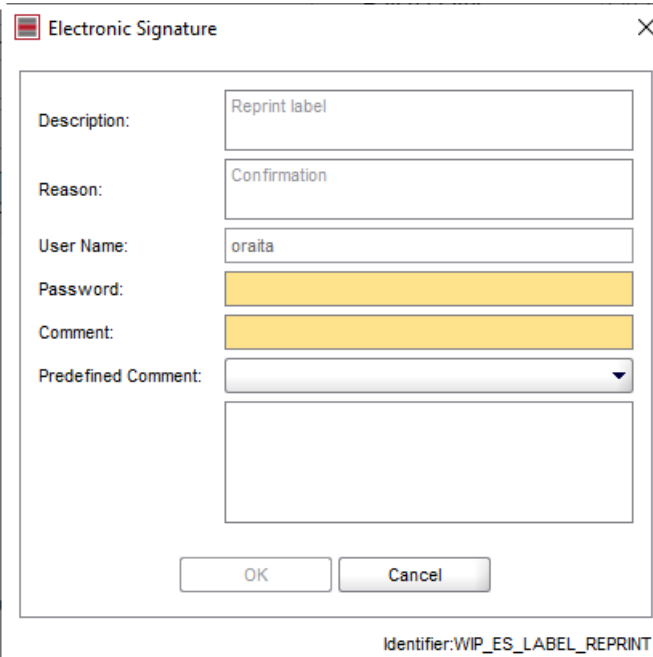


The dialog box is titled "Electronic Signature" and contains the following fields:

- Description: Reprint label
- Reason: Confirmation
- User Name: oraita
- Password: [Redacted]
- Comment: [Redacted]
- Predefined Comment: [Dropdown menu]

Buttons: OK, Cancel

Identifier: WIP_ES_LABEL_REPRINT

Figure 42: Single electronic signature with pre-defined comment

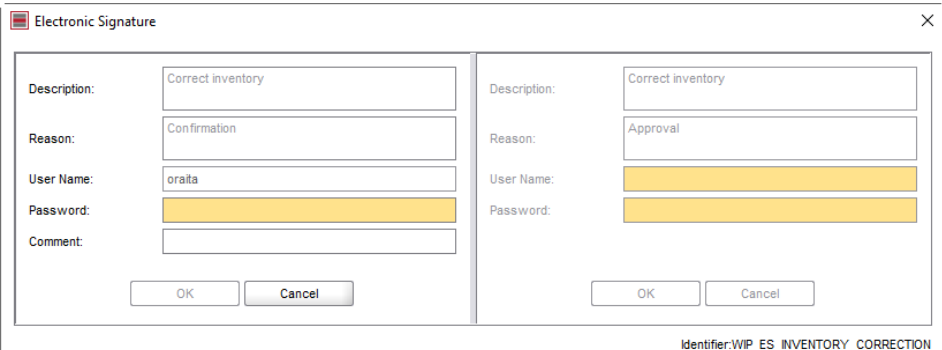

The dialog box is titled "Electronic Signature" and is split into two panes. The left pane contains the following fields:

- Description: Correct inventory
- Reason: Confirmation
- User Name: oraita
- Password: [Redacted]
- Comment: [Redacted]

Buttons: OK, Cancel

The right pane contains the following fields:

- Description: Correct inventory
- Reason: Approval
- User Name: [Redacted]
- Password: [Redacted]

Buttons: OK, Cancel

Identifier: WIP_ES_INVENTORY_CORRECTION

Figure 43: Double electronic signature to support witness role

INPUT BOX COLORS IN SIGNATURE DIALOGS

Input boxes can have different background colors that assist you with filling in the data of a form.

The different background colors indicate the following:



White

An unfilled input box. A blinking cursor indicates that it has the focus, and you can start to type your entry.



Yellow

A mandatory box that must be filled before the form can be saved or closed with the **OK** button. A blinking cursor indicates that it has the focus, and you can start to type your entry.

EBR processes can include mandatory, phase-specific phase completion signatures, which are displayed below the respective phases. You can only confirm such a phase when you have performed the signature.

After phase completion, the system displays the full name and login name of the person who performed the signature, along with the timestamp when it was performed. Phase completion signatures can be defined as single or as double signatures.

TIP

Please note that for security reasons the system deletes an already entered password unless you tap the confirm button within a reasonable amount of time after having entered the password.

A screenshot of a software interface for selecting a load type. It features three radio buttons: 'Get gross weight' (selected), 'Get net weight (Net)', and 'Get net weight (Gross)'. Below these are labels for 'Signature' and 'Operator' with a text input field containing 'oraita' and a yellow mandatory box. A 'Confirm' button with a QR code is on the right.

Figure 44: Single phase completion signature - active

A screenshot of the same software interface after completion. The 'Get gross weight' radio button is still selected. The 'Operator' field now shows 'Oupah Raita (oraita)' and a timestamp '05/09/2016 03:38:29 PM CEST'. The 'Confirm' button now has a green checkmark.

Figure 45: Single phase completion signature - completed

In addition to the phase completion signature for a specific phase, the system also provides a sequential phase completion signature to be performed as confirmation after one or more phases with completion signatures have been completed. The system requires a sequential phase completion signature to be performed by another user than its preceding plain phase completion signatures. Sequential signatures can only be defined as single signatures.

Figure 46: Sequential phase completion signature - active

Figure 47: Sequential phase completion signature - completed

Some phases provide specific phase action signatures that are triggered when an operator performs the related action during processing.

Figure 48: Phase action signature

Reference Selectors

When building complex data structures, it is often necessary to reference one or more objects from the object you are working on. For creating these references, the system supports you with single-reference selectors (page 99), for assigning one object to a base object.

If the number of objects that are available for selection is very large, the reference selectors include search and sort options (page 99) to restrict the list of objects available for referencing. If you do not specify any search criteria, the system will list all available objects of this type, which can take a while.

SINGLE-REFERENCE SELECTORS

The single-reference selector supports you when you need to assign another object to the base object you are currently working on. Whenever there is a reference to another object required and you know the object's identifier, you can simply type it in the respective box and the assignment is done. To undo the assignment, simply clear the box.

If you do not know the identifier, tap the ... button to open the reference selector for the assignment.



Figure 49: Single-reference box with selector button

The **Select Reference** dialog opens on top of the work area.

1. Specify search and sort criteria to reduce the number of objects from which you have to select your assignment object later. The system provides the following operators to refine your search criteria. Please note that all search criteria are case-sensitive.
 - ▣ **Starts with**
One or more characters or numbers the results will start with.
 - ▣ **Containing**
A sequence of one or more characters or numbers the results will contain.
 - ▣ **Ends with**
One or more characters or numbers the results will end with.
 - ▣ **Equal to**
Only if there is an exact match for the given sequence of characters or numbers will a search result be displayed. If you select this operator and do not type a criterion, no results will be found.
 - ▣ **Not equal to**
All search results other than the exact match for the given sequence of characters or numbers will be displayed.

▣ **Greater than**

For numbers: the search results will be higher than the given number.
For letters: the search results will consider all letters in the alphabet after the given letter.
For dates: the search results will be later than the given date.

▣ **Greater than or equal**

For numbers: the search results will be equal or higher than the given number.
For letters: the search results will consider all letters in the alphabet with and after the given letter.
For dates: the search results will be identical to or later than the given date.

▣ **Less than**

For numbers: the search results will be lower than the given number.
For letters: the search results will consider all letters in the alphabet before the given letter.
For dates: the search results will be earlier than the given date.

▣ **Less than or equal**

For numbers: the search results will be equal or lower than the given number.
For letters: the search results will consider all letters in the alphabet before and with the given letter.
For dates: the search results will be identical to or earlier than the given date.

▣ *** (asterisk)**

Select the asterisk operator without a criterion to see all objects.

▣ **Distinct only**

Check this option for the rare occurrence that your search criteria find variants of the same object several times. The respective objects are only displayed once then.

Basic Data	
Material Identifier	Starts with <input type="text"/>
Material description	Starts with <input type="text"/>
Material type	* <input type="text"/>
Number of rows	100

Figure 50: Search criteria

2. Choose the criterion you want to sort the list by. The system provides the following sorting operators to influence the order in which your search results are listed:

- ▣ **Ascending**

For numbers: smallest number first and rising.

For letters: from A to Z.

- ▣ **Descending**

For numbers: highest number first and going down.

For letters: from Z to A.

- ▣ **Disabled**

No Sorting.



The image shows a 'Sort by' label above a dropdown menu. The dropdown menu is open, displaying 'Material identifier' as the selected option. The text 'Material identifier' is in a blue font. To the right of the text is a small downward-pointing arrow icon.

Figure 51: Sorting criteria

3. Tap the **Search** button to start the search. The system lists all objects that match your search criteria in the **Results** list.
Tap the **Clear** button to clear all search criteria except the restriction to 100 rows.
4. In the **Results** list, tap the row with the object you wish to assign to your base object.

- Tap the **Cancel** button to close the reference selector and return to the work area without having assigned an object to your base object.

Select Reference - Topology

Filter by

Identifier

Starts with

Description

Starts with

Number of rows

100

☐ Distinct only

More

Storage area

Starts with

D010

Storage area description

Starts with

Sort by

Ascending

Clear

Search

Results: (16 / 16)

Identifier	Descr.	Storage area	Stor. area description
D010-10-10		D010-10	
D010-10-20		D010-10	
D010-10-30		D010-10	
D010-10-40		D010-10	
D010-20-10		D010-20	
D010-20-20		D010-20	
D010-20-30		D010-20	
D010-20-40		D010-20	
D010-30-10		D010-30	
D010-30-20		D010-30	
D010-30-30		D010-30	
D010-30-40		D010-30	
D010-40-10		D010-40	
D010-40-20		D010-40	

Select

Cancel

Figure 52: Single-reference selector

Help Access

The help system is accessible from the Cockpit via the **Help** button located in the **System** section of the **Startable** processes and workflows.

For keyboard use, press the F1 key to open the context-sensitive help in a web browser window.

In pre-defined workflows the **Help** button is also located in the upper right corner of the work area. The help system is context-sensitive on step level, which means that the active workflow step triggers the content of the help window.

The help window is non-modal and resizable.

In order to access other topics than the one directly related to the current context, use the navigation arrows located at the top of the help page, or related topic links located at the bottom of the page, if available.

The following additional features support your use of the help system:

- To access an overview of all available topics, open the **Contents** tab in the **Contents and Index** frame. The system additionally provides **Expand all** and **Collapse all** buttons to facilitate easier navigation in the contents tree.
- To access the index, open the **Index** tab in the **Contents and Index** frame.

- To use the **Search** function
 1. type the term you are looking for in the **Search** box and
 2. tap the **Search** button or press the ENTER key.

The system will display all occurrences of the search term in a third tab in the **Content and Index** frame.
- To print the page that is currently displayed in the help window, tap the **Print this page** button. The system displays a print preview of the page along with the default Windows **Print** dialog.

TIP

Please note that printing is only available from the stand-alone format of the help system and not from within the application.

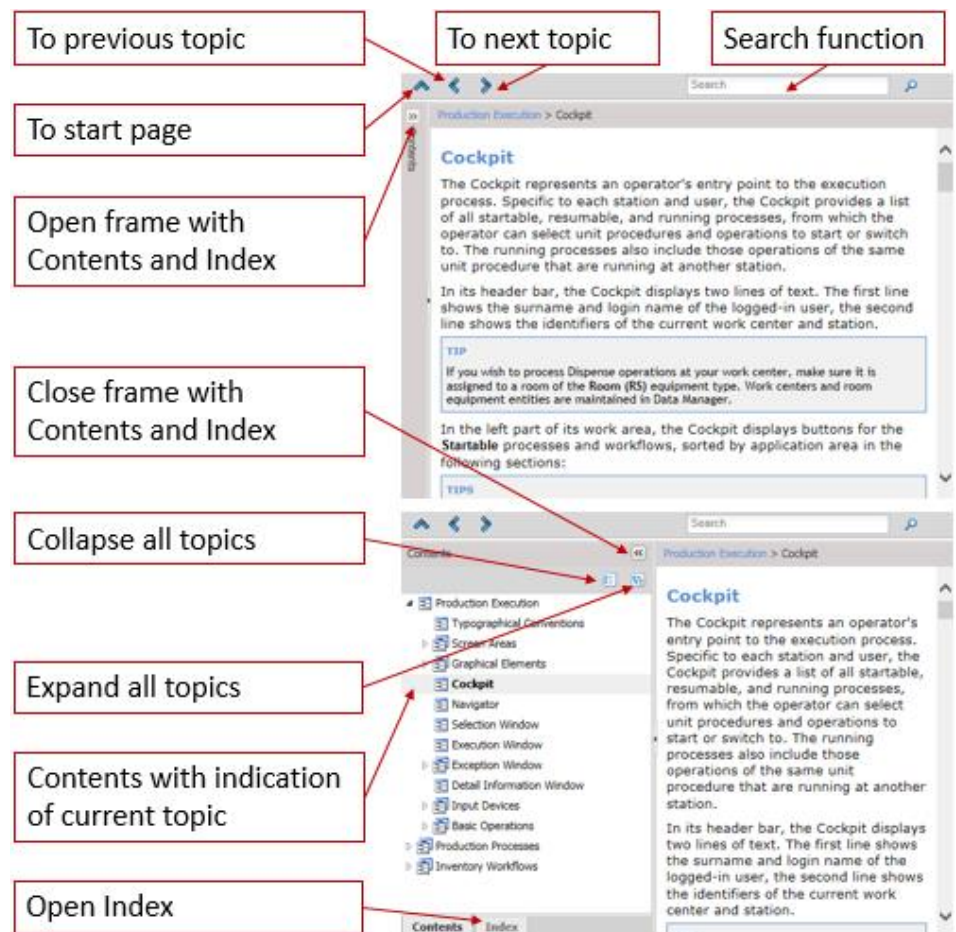


Figure 53: Help window navigation

-
-
- FT PharmaSuite® - Production Execution
-
-

Production Processes

The **Production** processes of PharmaSuite for Production Execution are accessible from the Cockpit (page 29).

TIP

Please note that the following list holds standard processes and workflows that are available unless configured to be hidden for this work center. Additional workflows for this application area type that can be created in the Workflow Designer application of PharmaSuite are added at the bottom of the list.

- Batch Processing (page 105)
- Material Processing (page 107)
- Workflow Processing (page 109)

Batch Processing

Batch Processing provides all functions and actions you need for executing the processes of a batch order.

It begins with the **Start Batch Processing** Selection Window where the system lists all order steps that are startable at your current work center.

TIP

Which order steps are displayed depends on the access rights of the logged-in user. Order steps of orders that are based on confidential master recipes to which the user does not have the required access rights remain invisible.

To restrict the number of listed order steps, you can filter for any string of characters in any column of the table of order steps. The filter becomes effective with the second character you type.

By default, the search applies to all table columns, is not case-sensitive, and the search string may occur anywhere within the content of a cell. Click the magnifying glass icon to set the search criteria to search only in specific columns, be case-sensitive, allow the use of wildcards or regular expressions, or match from start or exactly.

By default, the order steps are sorted alphabetically by their identifiers. It is, however, possible to reverse the sort order or sort by another column.

- To adjust the sorting, proceed as follows:
 1. Tap any column header to sort the table by this column in ascending order. The system indicates the sort order with a triangle pointing up.
 2. Re-tap the same column header to switch the sort order to descending. The system indicates the sort order with a triangle pointing down.
 3. CTRL-tap a yet unmarked column header to add this column as further sort level in ascending sort order, indicated by the triangle pointing up and a count number indicating the sort level.
 4. Re-CTRL-tap the same column header to switch the sort order to descending without changing its sort level.

When you have located the order step you wish to process, proceed as follows:

1. To start the order step, double-tap it.
The system returns you to the Cockpit (page 29) where it lists all executable operations as **Running**.

TIP

Please note that batch processes allow you to start processing by scanning the barcode of the batch or treatment ID for which production is to take place. A batch's barcode is typically located on the cover page of its order's batch report.
Scanning a batch or treatment ID barcode takes you directly to the Execution Window for phase processing if the order step's first operation is auto-startable.

2. Tap the operation you wish to execute.
The system displays the phases of the operation in their order of processing:
 - Previously executed phases are in read-only mode, with a checkmark on their **Confirm** buttons to indicate that they have been completed.
 - Currently active phases are accessible for processing.
Some phases may require an electronic signature before they can be completed.
 - Future phases that will only become active after one or all currently active phases have been executed are in read-only mode.
In a strictly sequential operation, the system displays all phases. If the operation includes one or more decision points, however, it will only show the future phases up to the next decision point.

Exceptions (page 44) that occur during processing can be recorded at all times.

3. Execute an accessible phase and tap its **Confirm** button.
4. Move on to the next accessible phase until you have all phases processed and confirmed.
When you have tapped the **Confirm** button of the last phase, the system returns you to the Cockpit for selecting a new operation or workflow.

Material Processing

Material Processing provides all functions and actions you need for processing the materials of several order steps together, thus building a weighing campaign. It is only available for Dispense work centers.

TIP

The materials of which order steps are displayed depends on the access rights of the logged-in user. Order steps of orders that are based on confidential master recipes to which the user does not have the required access rights remain invisible.

To restrict the number of listed materials, you can filter for any string of characters in any column of the table of order steps. The filter becomes effective with the second character you type.

By default, the search applies to all table columns, is not case-sensitive, and the search string may occur anywhere within the content of a cell. Click the magnifying glass icon to set the search criteria to search only in specific columns, be case-sensitive, allow the use of wildcards or regular expressions, or match from start or exactly.

Material processing begins with the **Select Materials** Selection Window where the system lists all materials of all order steps that are available for processing at your current work center.

By default, the materials are sorted alphabetically by their identifiers. It is, however, possible to reverse the sort order or sort by another column.

- To adjust the sorting, proceed as follows:
 1. Tap any column header to sort the table by this column in ascending order. The system indicates the sort order with a triangle pointing up.
 2. Re-tab the same column header to switch the sort order to descending. The system indicates the sort order with a triangle pointing down.
 3. CTRL-tap a yet unmarked column header to add this column as further sort level in ascending sort order, indicated by the triangle pointing up and a count number indicating the sort level.

4. Re-CTRL-tap the same column header to switch the sort order to descending without changing its sort level.

To support you with optimizing the sequence of materials you will process at your work center, the **Cleaning** column in the list displays the room cleaning demand, which the system determines for each material input item on the basis of the predecessor material processed at the work center. When you have located the material input items you wish to process together, proceed as follows:

1. Tap the checkbox in the leftmost column of a material to select it.

TIP

When order steps are available for processing at several work centers, some of the listed materials may already have been selected for processing at other work centers and are thus unavailable to you. These materials display their **Status** column with the identifier of the work center where they are being processed instead of their status.

2. When you have selected all materials you wish to process, tap the **Process** button. The system returns you to the Cockpit where it lists all operations that contain the materials you have selected from the **Running** processes.

TIP

Please note that you can use the **Back** button from the Cockpit to return to the **Select Materials** window and change your selection. The operations listed as **Running** will change accordingly.

3. Tap the operation you wish to execute.
The system displays the phases of the operation in their order of processing:
 - Previously executed phases are in read-only mode, with a checkmark on their **Confirm** buttons to indicate that they have been completed.
 - Currently active phases are accessible for processing.
Some phases may require an electronic signature before they can be completed.
 - Future phases that will only become active after one or all currently active phases have been executed are in read-only mode.
In a strictly sequential operation, the system displays all phases. If the operation includes one or more decision points, however, it will only show the future phases up to the next decision point.

Exceptions (page 44) that occur during processing can be recorded at all times.

4. Execute an accessible phase and tap its **Confirm** button.
5. Move on to the next accessible phase until you have all phases processed and confirmed.
When you have tapped the **Confirm** button of the last phase, the system returns you to the Cockpit for selecting a new operation or workflow.

Workflow Processing

Workflow Processing provides all functions and actions you need for executing the processes of a workflow.

It begins with the **Start Workflow Processing** Selection Window where the system lists all workflow steps that are startable at your current work center.

TIP

Which workflow steps are displayed depends on the access rights of the logged-in user. Workflow steps of workflows that are based on confidential master workflows to which the user does not have the required access rights remain invisible.

To restrict the number of listed order steps, you can filter for any string of characters in any column of the table of order steps. The filter becomes effective with the second character you type.

By default, the search applies to all table columns, is not case-sensitive, and the search string may occur anywhere within the content of a cell. Click the magnifying glass icon to set the search criteria to search only in specific columns, be case-sensitive, allow the use of wildcards or regular expressions, or match from start or exactly.

By default, the workflow steps are sorted alphabetically by their identifiers. It is, however, possible to reverse the sort order or sort by another column.

- To adjust the sorting, proceed as follows:
 1. Tap any column header to sort the table by this column in ascending order. The system indicates the sort order with a triangle pointing up.
 2. Re-tap the same column header to switch the sort order to descending. The system indicates the sort order with a triangle pointing down.
 3. CTRL-tap a yet unmarked column header to add this column as further sort level in ascending sort order, indicated by the triangle pointing up and a count number indicating the sort level.
 4. Re-CTRL-tap the same column header to switch the sort order to descending without changing its sort level.

When you have located the workflow step you wish to process, proceed as follows:

1. To start the workflow step, double-tap it.
The system returns you to the Cockpit (page 29) where it lists all executable operations as **Running**.

TIP

Please note that you can start processing of a workflow by scanning the barcode of its treatment ID.

Scanning a treatment ID barcode takes you directly to the Execution Window for phase processing if the workflow step's first operation is auto-startable.

2. Tap the operation you wish to execute.
The system displays the phases of the operation in their order of processing:
 - Previously executed phases are in read-only mode, with a checkmark on their **Confirm** buttons to indicate that they have been completed.
 - Currently active phases are accessible for processing.
Some phases may require an electronic signature before they can be completed.
 - Future phases that will only become active after one or all currently active phases have been executed are in read-only mode.
In a strictly sequential operation, the system displays all phases. If the operation includes one or more decision points, however, it will only show the future phases up to the next decision point.

Exceptions (page 44) that occur during processing can be recorded at all times.

3. Execute an accessible phase and tap its **Confirm** button.
4. Move on to the next accessible phase until you have all phases processed and confirmed.
When you have tapped the **Confirm** button of the last phase, the system returns you to the Cockpit for selecting a new operation or workflow.

Inventory Workflows

The pre-defined **Inventory** workflows of PharmaSuite for Production Execution are accessible from the Cockpit (page 29).

TIP

Please note that the following list holds pre-defined workflows that are available unless configured to be hidden for this work center. Additional one-click startable workflows for this type of application area, which you can create in the Workflow Designer application of PharmaSuite, are added at the bottom of the list.

- Material Receipt (page 111)
- Material Issue (page 119)
- Relocation (page 123)

TIP

If PharmaSuite is configured to communicate with Warehouse Management, the **Relocation** workflow is not available.

- Sublot Split (page 127)
- Inventory Correction (page 132)

Material Receipt

The **Material Receipt** workflow provides all functions and actions you need for introducing materials into an inventory.

What Is Material Receipt?

Material receipt is the pre-defined workflow for introducing materials into small storage areas in a shop-floor environment, such as buffer or line warehouses. It can be used for adding materials to existing batches or to generate new batches. At the end of a successful material receipt the fully identified and labeled material sublots are available for further processing on the shop floor. The material receives a batch identifier (page 114), the target storage location (page 116) is specified, its sublots (page 118) are registered individually, and one or more labels (page 119) are printed for sublot identification.

Material Receipt: Filter Materials

In order to receive your material into the system you have to pick it from the list of available materials. Since the number of materials in the system can be very large, it is often useful to apply a filter and thus reduce the number of materials listed.

1. The system provides the following operators to refine your search criteria. Please note that all search criteria are case-sensitive.

- ▣ **Starts with**

One or more characters or numbers the results will start with.

- ▣ **Containing**

A sequence of one or more characters or numbers the results will contain.

- ▣ **Ends with**

One or more characters or numbers the results will end with.

- ▣ **Equal to**

Only if there is an exact match for the given sequence of characters or numbers will a search result be displayed. If you select this operator and do not type a criterion, no results will be found.

- ▣ **Not equal to**

All search results other than the exact match for the given sequence of characters or numbers will be displayed.

- ▣ **Greater than**

For numbers: the search results will be higher than the given number.

For letters: the search results will consider all letters in the alphabet after the given letter.

For dates: the search results will be later than the given date.

- ▣ **Greater than or equal**

For numbers: the search results will be equal or higher than the given number.

For letters: the search results will consider all letters in the alphabet with and after the given letter.

For dates: the search results will be identical to or later than the given date.

- ▣ **Less than**

For numbers: the search results will be lower than the given number.

For letters: the search results will consider all letters in the alphabet before the given letter.

For dates: the search results will be earlier than the given date.

▣ **Less than or equal**

For numbers: the search results will be equal or lower than the given number.

For letters: the search results will consider all letters in the alphabet before and with the given letter.

For dates: the search results will be identical to or earlier than the given date.

▣ *** (asterisk)**

Select the asterisk operator without a criterion to see all objects.

2. Choose the criterion by which you wish to sort the list. The system provides the following sorting operators to influence the order in which your search results are listed:

▣ **Ascending**

For numbers: smallest number first and rising.

For letters: from A to Z.

▣ **Descending**

For numbers: highest number first and going down.

For letters: from Z to A.

▣ **Disabled**

No Sorting.



3. Tap the **Next** button to complete the workflow step and move to the next.

Material Receipt: Select a Material

The system lists the materials according to your search and sorting criteria.

1. Select the material you are about to receive.
2. Decide whether you need to create a new batch for your material or add the material to an existing batch.



Tap the **New batch** button to create a new batch.



Tap the **Select batch** button to select a batch from the list of available batches.

Select or Create a Batch

When materials are received into the system they are assigned to a batch. Typically, a new batch is created to receive a new material. In some cases, however, you will receive material that belongs to an already existing batch. Then you can select the batch and add the received material to it.

MATERIAL RECEIPT: CREATE A NEW BATCH

You have decided to create a new batch. To define the batch data, proceed as follows:

1. Type the potency if it is relevant for the material.
2. Select retest and expiry dates, if necessary.
3. Tap the **Create batch** button to have the system generate an identifier for the new batch.



4. Tap the **Next** button to complete the workflow step and move to the next.

TIP

If you need to use a pre-defined batch identifier instead of a system-generated one, type it in the **Batch identifier** box.
If your system is configured to support non-unique batch identifiers, the box will hold not only the batch identifier but also the material identifier.

MATERIAL RECEIPT: FILTER BATCHES

You have decided to add your material to an existing batch. Since the number of batches in the system can be very large, it is often useful to apply a filter and thus reduce the number of batches listed.

1. The system provides the following operators to refine your search criteria. Please note that all search criteria are case-sensitive.
 - **Starts with**
One or more characters or numbers the results will start with.
 - **Containing**
A sequence of one or more characters or numbers the results will contain.
 - **Ends with**
One or more characters or numbers the results will end with.

▣ **Equal to**

Only if there is an exact match for the given sequence of characters or numbers will a search result be displayed. If you select this operator and do not type a criterion, no results will be found.

▣ **Not equal to**

All search results other than the exact match for the given sequence of characters or numbers will be displayed.

▣ **Greater than**

For numbers: the search results will be higher than the given number.

For letters: the search results will consider all letters in the alphabet after the given letter.

For dates: the search results will be later than the given date.

▣ **Greater than or equal**

For numbers: the search results will be equal or higher than the given number.

For letters: the search results will consider all letters in the alphabet with and after the given letter.

For dates: the search results will be identical to or later than the given date.

▣ **Less than**

For numbers: the search results will be lower than the given number.

For letters: the search results will consider all letters in the alphabet before the given letter.

For dates: the search results will be earlier than the given date.

▣ **Less than or equal**

For numbers: the search results will be equal or lower than the given number.

For letters: the search results will consider all letters in the alphabet before and with the given letter.

For dates: the search results will be identical to or earlier than the given date.

▣ *** (asterisk)**

Select the asterisk operator without a criterion to see all objects.

2. Choose the criterion by which you wish to sort the list. The system provides the following sorting operators to influence the order in which your search results are listed:

▣ **Ascending**

For numbers: smallest number first and rising.
For letters: from A to Z.

▣ **Descending**

For numbers: highest number first and going down.
For letters: from Z to A.

▣ **Disabled**

No Sorting.



3. Tap the **Next** button to complete the workflow step and move to the next.

MATERIAL RECEIPT: SELECT A BATCH

The system lists the batches according to your search and sorting criteria.

1. Select the batch to which you want to add the material.



2. Tap the **Next** button to complete the workflow step and move to the next.

Material Receipt: Filter Storage Locations

In order to define the target storage location for your material you have to pick it from the list of available storage locations. Since the number of storage locations in the system can be very large, it is often useful to apply a filter and thus reduce the number of storage locations listed.

1. The system provides the following operators to refine your search criteria. Please note that all search criteria are case-sensitive.
 - ▣ **Starts with**
One or more characters or numbers the results will start with.
 - ▣ **Containing**
A sequence of one or more characters or numbers the results will contain.
 - ▣ **Ends with**
One or more characters or numbers the results will end with.

▣ **Equal to**

Only if there is an exact match for the given sequence of characters or numbers will a search result be displayed. If you select this operator and do not type a criterion, no results will be found.

▣ **Not equal to**

All search results other than the exact match for the given sequence of characters or numbers will be displayed.

▣ **Greater than**

For numbers: the search results will be higher than the given number.

For letters: the search results will consider all letters in the alphabet after the given letter.

For dates: the search results will be later than the given date.

▣ **Greater than or equal**

For numbers: the search results will be equal or higher than the given number.

For letters: the search results will consider all letters in the alphabet with and after the given letter.

For dates: the search results will be identical to or later than the given date.

▣ **Less than**

For numbers: the search results will be lower than the given number.

For letters: the search results will consider all letters in the alphabet before the given letter.

For dates: the search results will be earlier than the given date.

▣ **Less than or equal**

For numbers: the search results will be equal or lower than the given number.

For letters: the search results will consider all letters in the alphabet before and with the given letter.

For dates: the search results will be identical to or earlier than the given date.

▣ *** (asterisk)**

Select the asterisk operator without a criterion to see all objects.

2. Choose the criterion by which you wish to sort the list. The system provides the following sorting operators to influence the order in which your search results are listed:

▣ **Ascending**

For numbers: smallest number first and rising.
For letters: from A to Z.

▣ **Descending**

For numbers: highest number first and going down.
For letters: from Z to A.

▣ **Disabled**

No Sorting.



3. Tap the **Next** button to complete the workflow step and move to the next.

Material Receipt: Select a Storage Location

The system lists the storage locations according to your search and sorting criteria.

1. Select the storage location where the material will be stored.



2. Tap the **Next** button to complete the workflow step and move to the next.

Material Receipt: Create Sublots

You have selected a target storage location for your material. Now you need to specify the how many sublots you will store there and how much material they contain.

1. Type the number of sublots and the quantity per subplot with its unit of measure in the respective boxes. Type a comment, if necessary.
2. Tap the **Create sublots** button.
3. If you need to add further sublots, repeat the first two steps.



4. Tap the **Next** button to complete the workflow step and move to the next.

Material Receipt: Print Labels

You have specified all sublots you are about to store. The last step in the Material Receipt process is to print the labels for your sublots.

1. Type the number of labels you need to print per subplot in the respective box.
2. If you print more than one label per subplot, select whether you want to collate them or not.

TIP

Please note that the default system only allows one label per subplot. For this reason, the **Number of labels per subplot box** is preset with 1 and read-only. The **Printing Sequence** options are consequently disabled.

3. Tap the **Print labels** button.
4. If there was a problem during label printing that renders the printed labels unusable, tap the **Reprint labels** button to print again.



5. Tap the **OK** button to complete the step and finish the workflow.

Material Issue

The **Material Issue** workflow provides all functions and actions you need for removing materials from an inventory.

What Is Material Issue?

Material issue is the pre-defined workflow for removing materials, in the form of individual sublots, from the inventory. This can happen either because the material is required for processing and thus issued for a specific cost center or because the material is not fit for utilization and needs to be destroyed. The sublots are selected (page 122) individually and are issued (page 123) after the reason for their issue has been indicated.

Material Issue: Filter Sublots

In order to issue material, you have to pick the sublots concerned from the list of available sublots. Since the number of sublots in the system can be very large, it is often useful to apply a filter and thus reduce the number of sublots listed.



The **Scanner ready** icon in the navigation bar indicates that your scanner is ready for use. After you have successfully scanned a label, the system will move on to the **Select Sublots** step.

1. The system provides the following operators to refine your search criteria. Please note that all search criteria are case-sensitive.

- ▣ **Starts with**

One or more characters or numbers the results will start with.

- ▣ **Containing**

A sequence of one or more characters or numbers the results will contain.

- ▣ **Ends with**

One or more characters or numbers the results will end with.

- ▣ **Equal to**

Only if there is an exact match for the given sequence of characters or numbers will a search result be displayed. If you select this operator and do not type a criterion, no results will be found.

- ▣ **Not equal to**

All search results other than the exact match for the given sequence of characters or numbers will be displayed.

- ▣ **Greater than**

For numbers: the search results will be higher than the given number.

For letters: the search results will consider all letters in the alphabet after the given letter.

For dates: the search results will be later than the given date.

- ▣ **Greater than or equal**

For numbers: the search results will be equal or higher than the given number.

For letters: the search results will consider all letters in the alphabet with and after the given letter.

For dates: the search results will be identical to or later than the given date.

▣ **Less than**

For numbers: the search results will be lower than the given number.

For letters: the search results will consider all letters in the alphabet before the given letter.

For dates: the search results will be earlier than the given date.

▣ **Less than or equal**

For numbers: the search results will be equal or lower than the given number.

For letters: the search results will consider all letters in the alphabet before and with the given letter.

For dates: the search results will be identical to or earlier than the given date.

▣ *** (asterisk)**

Select the asterisk operator without a criterion to see all objects.

2. Choose the criterion by which you wish to sort the list. The system provides the following sorting operators to influence the order in which your search results are listed:

▣ **Ascending**

For numbers: smallest number first and rising.

For letters: from A to Z.

▣ **Descending**

For numbers: highest number first and going down.

For letters: from Z to A.

▣ **Disabled**

No Sorting.



3. Tap the **Next** button to complete the workflow step and move to the next.

Material Issue: Select Sublots

The system lists the sublots according to your search and sorting criteria.



The **Scanner ready** icon in the navigation bar indicates that your scanner is ready for use. After you have successfully scanned a label, the system will display the subplot in the list of selected sublots.

1. Select one or more sublots from the filtered list.



2. Tap the **Move selected down** button to move the sublots to the list of selected sublots.

If you have set your filter criteria to only list the sublots you want to issue:



Tap the **Move all down** button to move all sublots to the list of selected sublots.

To undo your selections and move sublots back up to filtered list of sublots:



Select one or more sublots in the list of selected sublots and tap the **Move selected up** button, or



Tap the **Move all up** button to move the sublots all back up.



3. Tap the **Next** button to complete the workflow step and move to the next.

Material Issue: Issue Sublots

You have selected the sublots you wish to issue.

1. Select the reason for the issue, **Destruction** or **Cost center**. If the issue is cost center-related, type the **Cost center identifier**.
2. Tap the **Issue sublots** button.



3. Tap the **OK** button to complete the step and finish the workflow.

Relocation

The **Relocation** workflow provides all functions and actions you need for moving materials within an inventory.

What Is Relocation?

Relocation is the pre-defined workflow for moving material, in the form of sublots, from one storage location to another. The sublots are selected (page 126) individually and then relocated (page 127) together.

Relocation: Filter Sublots

In order to relocate material, you have to pick the sublots from the list of available sublots. Since the number of sublots in the system can be very large, it is often useful to apply a filter and thus reduce the number of sublots listed.



The **Scanner ready** icon in the navigation bar indicates that your scanner is ready for use. After you have successfully scanned a label, the system will move on to the **Select Sublots** step.

1. The system provides the following operators to refine your search criteria. Please note that all search criteria are case-sensitive.
 - ▣ **Starts with**
One or more characters or numbers the results will start with.
 - ▣ **Containing**
A sequence of one or more characters or numbers the results will contain.
 - ▣ **Ends with**
One or more characters or numbers the results will end with.

▣ **Equal to**

Only if there is an exact match for the given sequence of characters or numbers will a search result be displayed. If you select this operator and do not type a criterion, no results will be found.

▣ **Not equal to**

All search results other than the exact match for the given sequence of characters or numbers will be displayed.

▣ **Greater than**

For numbers: the search results will be higher than the given number.
For letters: the search results will consider all letters in the alphabet after the given letter.
For dates: the search results will be later than the given date.

▣ **Greater than or equal**

For numbers: the search results will be equal or higher than the given number.
For letters: the search results will consider all letters in the alphabet with and after the given letter.
For dates: the search results will be identical to or later than the given date.

▣ **Less than**

For numbers: the search results will be lower than the given number.
For letters: the search results will consider all letters in the alphabet before the given letter.
For dates: the search results will be earlier than the given date.

▣ **Less than or equal**

For numbers: the search results will be equal or lower than the given number.
For letters: the search results will consider all letters in the alphabet before and with the given letter.
For dates: the search results will be identical to or earlier than the given date.

▣ *** (asterisk)**

Select the asterisk operator without a criterion to see all objects.

2. Choose the criterion by which you wish to sort the list. The system provides the following sorting operators to influence the order in which your search results are listed:

- ▣ **Ascending**

For numbers: smallest number first and rising.

For letters: from A to Z.

- ▣ **Descending**

For numbers: highest number first and going down.

For letters: from Z to A.

- ▣ **Disabled**

No Sorting.



3. Tap the **Next** button to complete the workflow step and move to the next.

Relocation: Select Sublots

The system lists the sublots according to your search and sorting criteria.



The **Scanner ready** icon in the navigation bar indicates that your scanner is ready for use. After you have successfully scanned a label, the system will display the subplot in the list of selected sublots.

1. Select one or more sublots from the filtered list.



2. Tap the **Move selected down** button to move the sublots to the list of selected sublots.

If you have set your filter criteria to only list the sublots you want to issue:



Tap the **Move all down** button to move all sublots to the list of selected sublots.

To undo your selections and move sublots back up to filtered list of sublots:



Select one or more sublots in the list of selected sublots and tap the **Move selected up** button, or



Tap the **Move all up** button to move the sublots all back up.



3. Tap the **Next** button to complete the workflow step and move to the next.

Relocation: Relocate Sublots

You have selected the sublots you wish to relocate.

1. In the **Target location** box, type the identifier of the storage location to which you wish to move the sublots. If you do not know the identifier, tap the ... button to open the reference selector (page 99) for storage locations that will support you with this workflow.

TIP

Please note that you can only relocate sublots that have the same source storage location.

If you need to relocate sublots from different source storage locations, please run the workflow for each of the source storage locations.

2. Tap the **Relocate sublots** button to book the sublots to their new storage location.



3. Tap the **OK** button to complete the step and finish the workflow.

Sublot Split

The **Sublot Split** workflow provides all functions and actions you need for splitting sublots into even or uneven parts.

What Is Sublot Split?

Sublots occasionally need to be split, for example to make material available simultaneously at two different locations. Splits can be even or uneven. An even split means that the base sublot is split in two or more new sublots with identical quantities, thus using up the base sublot. In an uneven split one or more sublots are split off from the base sublot, which itself also continues to hold material. The sublot to be split is selected (page 130), the number and quantities of the new sublots are specified (page 130) and registered, and one or more labels are printed (page 131) for sublot identification.

Sublot Split: Filter Sublots

In order to split a sublot you have to pick the sublot concerned from the list of available sublots. Since the number of sublots in the system can be very large, it is often useful to apply a filter and thus reduce the number of sublots listed.



The **Scanner ready** icon in the navigation bar indicates that your scanner is ready for use. After you have successfully scanned a label, the system will move on to the **Select Sublot** step.

1. The system provides the following operators to refine your search criteria. Please note that all search criteria are case-sensitive.

- ▣ **Starts with**

One or more characters or numbers the results will start with.

- ▣ **Containing**

A sequence of one or more characters or numbers the results will contain.

- ▣ **Ends with**

One or more characters or numbers the results will end with.

- ▣ **Equal to**

Only if there is an exact match for the given sequence of characters or numbers will a search result be displayed. If you select this operator and do not type a criterion, no results will be found.

- ▣ **Not equal to**

All search results other than the exact match for the given sequence of characters or numbers will be displayed.

- ▣ **Greater than**

For numbers: the search results will be higher than the given number.

For letters: the search results will consider all letters in the alphabet after the given letter.

For dates: the search results will be later than the given date.

- ▣ **Greater than or equal**

For numbers: the search results will be equal or higher than the given number.

For letters: the search results will consider all letters in the alphabet with and after the given letter.

For dates: the search results will be identical to or later than the given date.

▣ **Less than**

For numbers: the search results will be lower than the given number.

For letters: the search results will consider all letters in the alphabet before the given letter.

For dates: the search results will be earlier than the given date.

▣ **Less than or equal**

For numbers: the search results will be equal or lower than the given number.

For letters: the search results will consider all letters in the alphabet before and with the given letter.

For dates: the search results will be identical to or earlier than the given date.

▣ *** (asterisk)**

Select the asterisk operator without a criterion to see all objects.

2. Choose the criterion by which you wish to sort the list. The system provides the following sorting operators to influence the order in which your search results are listed:

▣ **Ascending**

For numbers: smallest number first and rising.

For letters: from A to Z.

▣ **Descending**

For numbers: highest number first and going down.

For letters: from Z to A.

▣ **Disabled**

No Sorting.



3. Tap the **Next** button to complete the workflow step and move to the next.

Sublot Split: Select a Sublot

The system lists the sublots according to your search and sorting criteria.



The **Scanner ready** icon in the navigation bar indicates that your scanner is ready for use. The system will show the scanned sublot as selected in the list.

1. Select the sublot you are about to split.
2. Decide whether you need to perform an even or an uneven split.



Tap the **Even split** button to split the sublot into even parts.



Tap the **Uneven split** button to split some parts off the base sublot.

Perform an Even or Uneven Split

When you split a sublot you need to decide how many new sublots you will require and how much material the new sublots are to contain. On the basis of these requirements you can either perform an even or an uneven split on your selected sublot.

SUBLOT SPLIT: PERFORM AN EVEN SPLIT

You have decided to split your base sublot into even parts. Proceed as follows:

1. In the **Count** box, type the number of new sublots into which you wish to split your base sublot.
2. Tap the **Split** button. The system now lists the new sublots.



3. Tap the **Next** button to complete the workflow step and move to the next.

SUBLOT SPLIT: PERFORM AN UNEVEN SPLIT

You have decided to split your base subplot into uneven parts. Proceed as follows:

1. In the **Count** box, type the number of new sublots you wish to split from your base subplot.
2. In the **Split quantity** box, type the quantity and unit of measure you wish each of the new sublots to have.
3. Tap the **Split** button. The system lists the new sublots. The **Quantity** box now displays the quantity of material that remains in your base subplot.
4. If you need to create further subplot splits, repeat the first three steps.



5. Tap the **Next** button to complete the workflow step and move to the next.

Sublot Split: Print Labels

You have created one or more new sublots. The last step in the **Sublot Split** process is to print the labels for your sublots.

1. Type the number of labels you need to print per subplot in the respective box.
2. If you print more than one label per subplot, select whether you want to collate them or not.

TIP

Please note that the default system only allows one label per subplot. For this reason, the **Number of labels per subplot box** is preset with 1 and read-only. The **Printing Sequence** options are consequently disabled.

3. Tap the **Print labels** button.
4. If there was a problem during label printing that renders the printed labels unusable, tap the **Reprint labels** button to print again.



5. Tap the **OK** button to complete the step and finish the workflow.

Inventory Correction

The **Inventory Correction** workflow provides all functions and actions you need for correcting the quantities of inventoried sublots.

What Is Inventory Correction?

An inventory correction becomes necessary if, during an inventory check for example, it turns out that a subplot contains more or less material than registered in the system. In these cases, the subplot quantity can be changed to the actual quantity of material contained in the subplot. The subplot whose quantity needs to be correct is selected (page 134) and its actual quantity (page 134) is registered.

Inventory Correction: Filter Sublots

In order to correct the quantity of a subplot you have to pick the subplot from the list of available sublots. Since the number of sublots in the system can be very large, it is often useful to apply a filter and thus reduce the number of sublots listed.



The **Scanner ready** icon in the navigation bar indicates that your scanner is ready for use. After you have successfully scanned a label, the system will move on to the **Change Quantity** step.

1. The system provides the following operators to refine your search criteria. Please note that all search criteria are case-sensitive.
 - ▣ **Starts with**
One or more characters or numbers the results will start with.
 - ▣ **Containing**
A sequence of one or more characters or numbers the results will contain.
 - ▣ **Ends with**
One or more characters or numbers the results will end with.
 - ▣ **Equal to**
Only if there is an exact match for the given sequence of characters or numbers will a search result be displayed. If you select this operator and do not type a criterion, no results will be found.
 - ▣ **Not equal to**
All search results other than the exact match for the given sequence of characters or numbers will be displayed.

■ **Greater than**

For numbers: the search results will be higher than the given number.

For letters: the search results will consider all letters in the alphabet after the given letter.

For dates: the search results will be later than the given date.

■ **Greater than or equal**

For numbers: the search results will be equal or higher than the given number.

For letters: the search results will consider all letters in the alphabet with and after the given letter.

For dates: the search results will be identical to or later than the given date.

■ **Less than**

For numbers: the search results will be lower than the given number.

For letters: the search results will consider all letters in the alphabet before the given letter.

For dates: the search results will be earlier than the given date.

■ **Less than or equal**

For numbers: the search results will be equal or lower than the given number.

For letters: the search results will consider all letters in the alphabet before and with the given letter.

For dates: the search results will be identical to or earlier than the given date.

■ *** (asterisk)**

Select the asterisk operator without a criterion to see all objects.

2. Choose the criterion by which you wish to sort the list. The system provides the following sorting operators to influence the order in which your search results are listed:

▣ **Ascending**

For numbers: smallest number first and rising.
For letters: from A to Z.

▣ **Descending**

For numbers: highest number first and going down.
For letters: from Z to A.

▣ **Disabled**

No Sorting.



3. Tap the **Next** button to complete the workflow step and move to the next.

Inventory Correction: Select a Sublot

The system lists the sublots according to your search and sorting criteria.



The **Scanner ready** icon in the navigation bar indicates that your scanner is ready for use. After you have successfully scanned a label, the system will move on to the **Change Quantity** step.

1. Select the subplot you are about to correct.



2. Tap the **Next** button to complete the workflow step and move to the next.

INVENTORY CORRECTION: CHANGE THE QUANTITY

You have selected a subplot to correct its quantity. Proceed as follows:

1. In the **Quantity** box, type the actual quantity and unit of measure of the subplot.
2. Tap the **Set quantity** button. The system now registers the new quantity for the subplot.



3. Tap the **OK** button to complete the step and finish the workflow.

A

- Abort and reactivate unit procedure (Cockpit action) • 11
- Abort operation (Cockpit action) • 11
- About dialog • 26
- Action (Cockpit) • 67
 - Event-triggered operation • 71
 - Operation • 70
 - Unit procedure • 68
- Actions (Cockpit button) • 11
- Alarm (common phase indicator) • 20
- Append (Cockpit action) • 11
- Appendable unit procedure • 74
- Appendable unit procedures • 74
- Appendable workflows • 75
- Auto-lock mechanism • 26

B

- Back (EBR framework button) • 8
- Back (framework button) • 8
- Basic operation • 49
 - Aborting an operation • 81
 - Aborting and reactivating a unit procedure • 79
 - Aborting phases • 65
 - Appending workflows • 73
 - Cockpit action • 67
 - Detaching a unit procedure • 78
 - Detaching an operation • 80
 - Electronic signature • 81
 - Escalation handling • 77
 - Exception handling • 55
 - Help access • 88
 - Login • 50
 - Logout • 52
 - Password change • 52
 - Pausing unit procedures • 72
 - Repairing phases • 65
 - Signature • 81

Start • 50

- Batch Processing • 91
- Bookmark (common phase button) • 20

C

- Campaign weighing • 93
- Cancel all runs (Cockpit action) • 11
- Cancel current workflow (framework button) • 8
- Cancel order step (common workflow button) • 20
- Cancel run (Cockpit action) • 11
- Cancel workflow (cockpit action) • 12
- Change password • 52
- Change station • 50
- Change user • 26
- Checkmark (common phase icon) • 20
- Cockpit • 26
 - Searching • 26
 - Sorting • 26
- Cockpit (framework button) • 8
- Cockpit button • 11
 - Abort and reactivate unit procedure (action) • 11
 - Abort operation (action) • 11
 - Actions • 11
 - Append (action) • 11
 - Cancel all runs (action) • 11
 - Cancel run (action) • 11
 - Cancel workflow (action) • 12
 - Currently active operation (button) • 12
 - Detach operation • 12
 - Event-triggered operation run before expiry (button) • 12
 - Event-triggered operation run before overdue (button) • 12
 - Event-triggered operation run before reminder (button) • 12
 - Exceptions (action) • 13
 - New run (action) • 13
 - Open in viewer (button) • 13

- Operation running at different station (button) • 13
- Pause unit procedure (action) • 14
- Remove template (action) • 14
- Resume unit procedure (action) • 14
- Running operation (button) • 15
- Running unit procedure of a one-click startable workflow (button) • 16
- Running unit procedure of a one-click startable workflow marked by Cockpit search (button) • 17
- Running unit procedure of a one-click startable workflow with treatment ID (button) • 16
- Running unit procedure of Batch Processing (button) • 15
- Running unit procedure of Batch Processing marked by Cockpit search (button) • 17
- Running unit procedure of Batch Processing with treatment ID (button) • 15
- Running unit procedure of Workflow Processing (button) • 15
- Sort ascending (button) • 17
- Sort by timestamp (button) • 17
- Sort by treatment ID (button) • 17
- Sort descending (button) • 18
- Startable batch orders (button) • 18
- Startable material for processing (button) • 18
- Startable operation (button) • 18
- Startable workflow (button) • 18
- Startable workflows (button) • 18
- Template (button) • 18
- Cockpit icon • 11
 - Run due reminder (icon) • 14
 - Run overdue (icon) • 14
 - Unit procedure paused (icon) • 19
- Commenting an exception • 63
- Consume subplot (specific phase button) • 24
- Control bar (screen area) • 2
- Conventions (typographical) • 1
- Currently active operation (Cockpit button) • 12

D

- Detach operation (Cockpit button) • 12

- Detail Information Window • 44
- Device configuration • 26
- Documentation • 50
- Double signature • 81
- Due (framework graphic) • 9
- Due reminder (framework graphic) • 9

E

- EBR server state • 26
- EBR workflow • 91
- Enable (common phase button) • 20
- Error (common phase icon) • 20
- Escalation handling • 77
- Even split (specific workflow button) • 24
- Event-triggered operation (actions) • 71
- Event-triggered operation run before expiry (Cockpit button) • 12
- Event-triggered operation run before overdue (Cockpit button) • 12
- Event-triggered operation run before reminder (Cockpit button) • 12
- Exception (framework button) • 9
- Exception at value (common phase icon) • 20
- Exception handling • 55
 - Bookmarking an exception • 64
 - Recording a comment • 63
 - Recording a system-triggered exception • 62
 - Recording a user-triggered free-text exception • 60
 - Recording a user-triggered situation-specific exception • 61
 - System-triggered exception • 55
 - User-triggered exception • 57
- Exception Window • 39
 - Operation-specific • 41
 - Overview • 41
 - Phase-specific • 40
- Exceptions (Cockpit action) • 13
- Execution Window • 37

F

- Finish due reminder (framework graphic) • 9
- Finish overdue (framework graphic) • 9

First page of document (common phase button) • 21

Framework button • 8

Back • 8

Back (EBR) • 8

Cancel current workflow • 8

Cockpit • 8

Exception • 9

Help • 9

Mark notification as read • 19

Mark notification as unread • 19

Navigator • 9

Next • 9

OK • 10

Open Cockpit • 10

Page down • 10

Page up • 10

Pin notification • 19

Remove notification • 19

Running at current work center • 10

Switch • 10

Unpin notification • 19

Framework graphic • 8

Due • 9

Due reminder • 9

Finish due reminder • 9

Finish overdue • 9

Overdue • 10

Paused unit procedure • 10

Free-text exception • 57

H

Header bar (screen area) • 2

Help (framework button) • 9

Help access • 88

I

Identify batch manually (common workflow button) • 21

Identify load carrier manually (common workflow button)
• 21

Identify subplot manually (common workflow button) • 21

Input box • 81

Mandatory property • 81

K

Keyboard operation • 45

Cockpit navigation • 47

Execution Window navigation • 47

Framework shortcut • 45

Navigation bar • 46

L

Last page of document (common phase button) • 21

Loading (common phase icon) • 21

Locking the system • 26

Login • 50

Change password • 52

Change station • 50

Password • 50

Logout • 52

M

Manual scale (indicator in Dispense and Weighing) • 24

Mark notification as read (framework button) • 19

Mark notification as unread (framework button) • 19

Material Processing • 93

Move all down (common workflow button) • 21

Move all up (common workflow button) • 22

Move selected down (common workflow button) • 22

Move selected up (common workflow button) • 22

N

Navigator • 34

Navigator (framework button) • 9

New batch (specific workflow button) • 24

New run (Cockpit action) • 13

Next (framework button) • 9

Next page in document (common phase button) • 22

Notification panel • 42

Notification panel area (screen area) • 2

O

OK (framework button) • 10

Open Cockpit (framework button) • 10
 Open in viewer (Cockpit button) • 13
 Operation (actions) • 70
 Operation aborting • 81
 Operation detaching • 80
 Operation running at different station (Cockpit button) • 13
 Overdue (framework graphic) • 10

P

Page down (framework button) • 10
 Page up (framework button) • 10
 Password • 52
 Password change • 52
 Pause unit procedure (Cockpit action) • 14
 Paused unit procedure (framework graphic) • 10
 Phase aborting • 65
 Phase activity ongoing (common phase icon) • 22
 Phase button (common) • 20

- Bookmark • 20
- Enable • 20
- First page of document • 21
- Last page of document • 21
- Next page in document • 22
- Previous page in document • 22
- Scroll down in document • 23
- Scroll up in document • 23
- Toggle unit of measure • 23

 Phase button (specific) • 23

- Consume subplot • 24
- Release scale • 24

 Phase icon (common) • 20

- Alarm • 20
- Checkmark • 20
- Error • 20
- Exception at value • 20
- Loading • 21
- Phase activity ongoing • 22

 Phase icon (specific) • 23

- Manual scale (indicator in Dispense and Weighing) • 24
- Scale loaded (indicator in Dispense or Weighing) • 24

Scan container (scan action indicator in Dispense or Weighing) • 24
 Scan printer (scan action indicator in Dispense or Weighing) • 25
 Scan scale (scan action indicator in Dispense or Weighing) • 25
 Scan scale with value (scan action indicator in Dispense or Weighing) • 25
 Scan subplot (scan action indicator in Dispense or Weighing) • 25
 Scan undefined scale (scan action indicator in Dispense or Weighing) • 25
 Phase repairing • 65
 Phase-specific exception • 57
 Pin notification (framework button) • 19
 Previous page in document (common phase button) • 22
 Process (EBR) • 91

- Batch Processing • 91
- Material Processing • 93

 Production process • 91

- Batch Processing • 91
- Material Processing • 93
- Workflow Processing • 95

R

Reference selector • 85

- Search criteria • 85
- Single-reference • 85
- Sort criteria • 85

 Register at station • 26
 Release scale (specific phase button) • 24
 Remove notification (framework button) • 19
 Remove template (Cockpit action) • 14
 Resume unit procedure (Cockpit action) • 14
 Run due reminder (Cockpit icon) • 14
 Run overdue (Cockpit icon) • 14
 Running at current work center (framework button) • 10
 Running operation (Cockpit button) • 15
 Running unit procedure of a one-click startable workflow (Cockpit button) • 16

Running unit procedure of a one-click startable workflow marked by Cockpit search (Cockpit button) • 17

Running unit procedure of a one-click startable workflow with treatment ID (Cockpit button) • 16

Running unit procedure of Batch Processing (Cockpit button) • 15

Running unit procedure of Batch Processing marked by Cockpit search (Cockpit button) • 17

Running unit procedure of Batch Processing with treatment ID (Cockpit button) • 15

Running unit procedure of Workflow Processing (Cockpit button) • 15

S

Scale loaded (indicator in Dispense or Weighing) • 24

Scan container (scan action indicator in Dispense and Weighing) • 24

Scan printer (scan action indicator in Dispense or Weighing) • 25

Scan scale (scan action indicator in Dispense or Weighing) • 25

Scan scale with value (scan action indicator in Dispense or Weighing) • 25

Scan subplot (scan action indicator in Dispense or Weighing) • 25

Scan undefined scale (scan action indicator in Dispense or Weighing) • 25

Scanner ready - mandatory (common workflow icon) • 22

Scanner ready - optional (common workflow icon) • 23

Screen area • 2

- Control bar • 2
- Header bar • 2
- Navigation bar • 7
- Notification panel area • 2
- Status bar • 7
- Work area (EBR) • 2
- Work area for pre-defined workflows • 7

Scroll down in document (common phase button) • 23

Scroll up in document (common phase button) • 23

Searching the Cockpit • 26

Select batch (workflow button) • 25

Selected (common workflow icon) • 23

Selection Window • 37

Signature • 81

Single signature • 61

Skip current step (common workflow button) • 23

Sort ascending (Cockpit button) • 17

Sort by timestamp (Cockpit button) • 17

Sort by treatment ID (Cockpit button) • 17

Sort descending (Cockpit button) • 18

Sorting the Cockpit • 26

Start PharmaSuite • 50

Startable batch orders (Cockpit button) • 18

Startable material for processing (Cockpit button) • 18

Startable operation (Cockpit button) • 18

Startable workflow (Cockpit button) • 18

Startable workflows (Cockpit button) • 18

Switch (framework button) • 10

System-triggered exception • 55

T

Template (Cockpit button) • 18

Toggle unit of measure (common phase icon) • 23

U

Uneven split (specific workflow button) • 26

Unit procedure (actions) • 68

Unit procedure aborting and reactivating • 79

Unit procedure detaching • 78

Unit procedure paused (Cockpit icon) • 19

Unit procedure pausing • 72

Unpin notification (framework button) • 19

Unselected (common workflow icon) • 23

User-triggered exception • 57

W

Work area (screen area) • 2

Workflow appending • 73

Workflow button (common) • 20

- Cancel order step • 20
- Identify batch manually • 21
- Identify load carrier manually • 21

Identify subplot manually •	21
Move all down •	21
Move all up •	22
Move selected down •	22
Move selected up •	22
Skip current step •	23
Workflow button (specific) •	23
Even split •	24
New batch •	24
Select batch •	25
Uneven split •	26
Workflow icon (common) •	20
Scanner ready - mandatory •	22
Scanner ready - optional •	23
Selected •	23
Unselected •	23
Workflow icon (specific) •	23
Workflow Processing •	95